
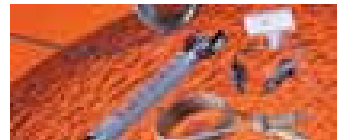
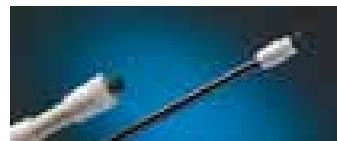
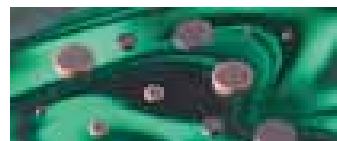
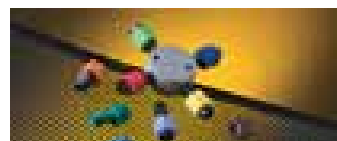
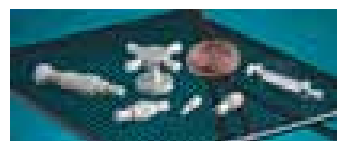
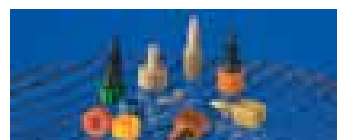
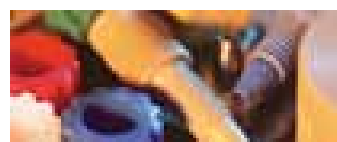


A photograph of a Vici Jour chromatography manifold assembly. The assembly consists of three white cylindrical manifold blocks connected by a network of colored tubing. The tubing includes red, blue, green, and black lines, some of which are thin-bore and others thicker. Various colored fittings (green, orange, blue, red, black, yellow) are used to connect the tubes to the manifold ports. The entire setup is arranged in a circular pattern on a light gray background.

**VICI Jour<sup>®</sup>**

**Chromatography  
Fluid/Gas Transfer**

CATALOG 15 INT

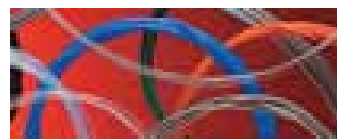

**2** Products for UHPLC

**4** Tubing

**16** Fittings introduction

**18** High pressure polymeric fittings

**30** High pressure stainless fittings

**36** Nanovolume® HPLC fittings

**40** Low pressure fittings

**50** Filters and mobile phase filters

**58** Safety products

**72** Frits

**76** Columns and guard columns

**80** Sample loops

**84** Accessories

**92** General reference

**108** Index by product number

**110** Index by name of product

**112** Order information

**Europe, Middle East, Asia,  
Africa, Australia, New Zealand,  
and other countries**

#### VICI AG International

##### Sales

Phone: Int + 41·41·925-6200

Fax: Int + 41·41·925-6201

Email: [order@vici.ch](mailto:order@vici.ch)

##### Technical/Service

Phone: Int + 41·41·925-6200

Fax: Int + 41·41·925-6201

Email: [support@vici.ch](mailto:support@vici.ch)

Parkstrasse 2

CH-6214 Schenkon / Switzerland

[info@vici.ch](mailto:info@vici.ch)

**North America, South America,  
and Central America  
(except Canada)**

#### VICI Precision Sampling

##### Sales

Toll-free: 800-828-1653

Phone: 225-927-1128

Fax: 225-923-1331

Email: [precsamp@vici.com](mailto:precsamp@vici.com)

##### Technical/Service

Phone: 225-927-1128

Fax: 225-923-1331

Email: [precsamp@vici.com](mailto:precsamp@vici.com)

P.O. Box 15886

Baton Rouge, Louisiana / USA

[precsamp@vici.ch](mailto:precsamp@vici.ch)

#### Canada

#### VICI Valco Instruments Canada

##### Sales

Toll-free: 866-297-2626

Phone: 613-342-2600

Fax: 613-342-0111

Email: [canada@vici.com](mailto:canada@vici.com)

##### Technical/Service

Phone: 613-342-2600

Fax: 613-342-0111

Email: [canada@vici.com](mailto:canada@vici.com)

26 Water Street East

Bockville, Ontario / Canada K6V 1A1

[canada@vici.com](mailto:canada@vici.com)

## VICI AG International

VICI Jour has been a part of VICI AG International in Switzerland since 2004. The complete Jour line of accessories and consumables for fluid transfer in GC and LC is available from VICI AG International, VICI Valco Canada, VICI Precision Sampling, and through our international distributor network.

VICI AG International, in Schenkon, Switzerland, is an independent site for the manufacturing of Valco and Cheminert valves and the handling of all VICI product lines in Europe, Africa, the Middle East, and the Far East.

VICI AG International has obtained the ISO certificate to the standard of 9001:2008, and is the first ISO-certified manufacturer of high pressure and low pressure valves for scientific equipment.

## VICI Jour product line

Jour products include a complete line of tubing and capillaries in a wide variety of diameters and materials, high and low pressure fingertight fittings and connection accessories, column hardware, high and low pressure filters, cutters and other tools, backpressure regulators, and much more.

We offer accessories and consumables for fluid transfer in GC and LC, including:

- Fittings - high pressure polymeric and stainless, Nanovolume® PEEK and stainless, and low pressure
- Tubing and sample loops
- Safety products and mobile phase handling systems
- Filters and frits
- PEEK columns and guard columns



## PRODUCTS FOR UHPLC

VICI Jour now offers high quality, high pressure fittings, filters, and tubing for use in Ultra-High Performance Liquid Chromatography (UHPLC).

### Page 15 **Stainless steel tubing – premium grade** **UHPLC**

- Seamless Type 316 stainless steel
- Variety of sizes for UHPLC, HPLC, and GC applications
- Smooth internal surface
- Soft annealed OD for easy ferrule connections
- Available in 3 m and 10 m lengths

### Page 19 **PEEK one-piece hex-head fittings** **NEW** **UHPLC**

- Made of special PEEK HT™ for better pressure resistance
- One-piece – no ferrule required
- Biocompatible



### Page 22 **PEEK ferrules** **NEW** **UHPLC**

- Biocompatible 1/16" ferrule for ultra-high pressures
- Alternative to titanium ferrules in applications with titanium valves
- Excellent in combination with stainless nuts



### Pages 38-39 **Nanovolume® 360 micron fittings** **UHPLC**

- For direct connection of 360 µm tubing
- Work with metal, fused silica, or PEEK
- Up to 20,000 psi with metal tubing
- Snap-in rotating ferrule for “one-piece” fitting with no tubing twist

Our high pressure fittings permit direct connection of 360 micron OD fused silica, PEEK, stainless, or electroformed nickel tubing without the use of liners.



### Pages 31, 32 **Valco® stainless steel nuts and ferrules** **UHPLC**

- Excellent fit to Valco fitting details
  - Different lengths available
- These fittings are for use with all standard Valco internal fittings and most valves. They may be used with fittings from other manufacturers as well. Valco metal ferrules cut a ring near the end of the tube, preventing tube release at high pressures without significantly deforming and restricting the tube interior.



### Page 32 **Stainless steel plugs and caps** **UHPLC**

- For 1/32", 1/16", and 1/8" OD tubing

Plugs close off an unused port in a valve, manifold, or HPLC column. Caps close off a tube with a ferrule and internal nut already made up.




**Page 34 Stainless steel unions**
**UHPLC**

- Choice of four bores: 0.15 mm, 0.25 mm, 0.50 mm, or 0.75 mm
- For 1/16" OD tubing


**Page 34 Stainless steel reducing unions, imperial**
**UHPLC**

These reducing unions join 1/16" OD to 1/8" OD tubing.


**Page 35 Stainless steel reducing unions, imperial to metric**
**NEW**
**UHPLC**

- For 1 mm, 2 mm, and 3 mm OD tubing
- With these new reducing unions, metric tubing can be connected to imperial fittings, valves, and chromatographic systems.


**Page 33 Stainless steel tees and crosses**
**UHPLC**

- For 1/16" OD tubing
- 0.25 mm, 0.50 mm, and 0.75 mm bores


**Page 56 Stainless steel precolumn filter**
**NEW**
**UHPLC**

- Fewer connections compared to in-line filters
- Works with 1/16" column connections from all manufacturers
- Three porosity options: 0.5  $\mu$ m, 2  $\mu$ m, and 5  $\mu$ m

The design of this precolumn filter reduces the number of connections compared to a standard in-line filter, increasing the stability and reducing the potential of leakage.


**Page 57 Stainless steel in-line filter**
**UHPLC**

- Minimized hold-up volume
- Two bores for analytical and semi-prep flow rates
- Three porosity options: 0.5  $\mu$ m, 2  $\mu$ m, and 5  $\mu$ m

The stainless steel in-line filter has virtually no hold-up volume and can be used in various applications without significant band broadening or loss of efficiency.


**Page 90 Stainless steel starter kit**
**UHPLC**

The kit consists of stainless fittings, ferrules, tubing, and solvent filters, as well as a union, our cutting-wheel type metal tubing cutter, and other useful items like our ValvTool. All boxed together for a very competitive price.

# Tubing

## METAL AND POLYMERIC

We offer chromatography grade tubing in a variety of polymers as well as stainless steel. Outside diameters are 360 µm, 1/32", 1/16", 1/8", and 1/4".

## Tubing index

### PEEK tubing

Natural PEEK.....page 5

Available in ODs from 360 micron to 1/4", in custom lengths (priced per meter) and packages of 3 and 10 meters.

Striped color-coded PEEK ..... 6

A stripe is added to the outside of the tubing, so dye never contacts the fluid stream.

Dual-layer PEEK..... 7

With an inner layer of natural PEEK and an outer layer of color-coded PEEK, there is no possibility of stream contents coming in contact with coloring agents.

Solid color-coded PEEK..... 8

Color coding identifies the ID or your PEEK tubing. Each ID is coded with a different color.

### Other polymeric tubing

ETFE..... 9

FEP

Solid color-coded ..... 10

Natural ..... 11

PFA ..... 9

PTFE

Solid color-coded ..... 12

Natural ..... 12

Thermo-formed polymeric tubing ..... 13

Impermeable tubing ..... 13

FEP coated with PVDF combines the chemical resistance of FEP with the tight, impermeable structure of PVDF.

### Stainless steel tubing

Pre-cut SS tubing

Pre-cut tubing ..... 14

Flexible pre-cut tubing ..... 14

1/32" tubing with 1/16" ends.

Bulk SS tubing

Premium grade, seamless ..... 15

Standard grade, welded seam ..... 15

### DISCLAIMER

The maximum holding pressure for any type of connections involving tubing and a ferrule varies considerably with the tubing material, the ferrule material, the clearance between tubing OD and ferrule ID, and the shape of the fitting detail.

### RELATED PRODUCTS

SolviFlex

coiled tubing ..... 70

### Tubing tools and accessories

Easy-Flange kit ..... 85

Easy-Flange combi kit . 85

PEEK tubing elbows ... 86

Tubing cutters for

polymeric tubing

Clean-cut..... 86

Guillotine..... 86

Tubing cutters for

stainless tubing

Cutting wheel type .. 87

Plier type ..... 87

Tubing clip – the LC

tubing organizer..... 86



**PEEK TUBING TOLERANCES**

Due to tighter tolerances, VICI Jour PEEK tubing is the best choice for chromatography applications in which exact flow rates or tube volume are of importance. The tighter tolerance ensures a more consistent internal volume of transfer lines.

Tubing OD	OD (±) (mm)	ID (±) (mm)
360 µm	0.025	0.013
1/32"	0.025	0.013
1/16" *	0.025	0.013
1/16" **	0.025	0.025
1/8"	0.05	0.05
* For ID up to 0.25 mm		
** For ID over 0.25 mm		

**OPTIONS**

Contact your distributor for more information about:

- Longer bulk tubing lengths
- Customized ODs and IDs (30 m minimum quantity)

**TECH TIP**

- A clean burr-free perpendicular cut can be achieved with our Clean-Cut tubing cutter. . . . . 86
- To bend PEEK tubing at the optimal radius, use our tubing elbows. . . 86

**MORE INFORMATION**

About PEEK® . . . . . 99  
 Chemical resistance  
   PEEK . . . . . 96-98  
   PEEK at higher temperatures . . . . . 96  
 Maximum working pressure . . . . . 94  
 Tubing ID chart:  
   volume per length . . . 93  
 Tubing tolerances . . . . 95

**RELATED PRODUCTS**

Clean-cut tubing cutter . . . . . 86  
 Easy-Flange kit . . . . . 85  
 PEEK fittings  
   High pressure . . . 19-29  
   Low pressure . . . 42, 44-47  
 PEEK tubing elbows . . . 86

VICI offers natural and color-coded PEEK® tubing. PEEK tubing has the strength required to withstand continuous use at HPLC pressure without swelling or bursting, and is not affected by halide salts, high strength buffers, or other aggressive mobile phases that corrode stainless steel. The polymer surface will not leach metal ions into the eluent or extract metal-sensitive components from the sample. Note however that dichloromethane, THF, and DMSO may cause swelling in PEEK, and concentrated nitric and sulphuric acid will attack PEEK.

**PEEK tubing – natural**

Qty/pkg:		3 m	10 m	per meter
OD	ID	Product No.	Product No.	Product No.
360 µm	75 µm	JR-T-6060-M3	JR-T-6060-M10	—
1/32"	0.13 mm	JR-T-5993-M3	JR-T-5993-M10	—
1/16"	0.064 mm	JR-T-5998-M3	JR-T-5998-M10	—
1/8"	0.75 mm	JR-T-6004-M3	JR-T-6004-M10	—
	1.59 mm	JR-T-60041-M3	JR-T-60041-M10	—
	2.00 mm	JR-T-60042-M3	JR-T-60042-M10	—
1/4"	3.17 mm	—	—	JR-T-6006

**SPECS**

- Maximum recommended working temp: 100°C (continuous) for 1/16" OD tubing with ID up to 0.75 mm
- Max working pressures: see page 94
- Tubing tolerances: see page 95



**Natural PEEK tubing**  
 shown with ODs of 1/8", 1/16", 1/32", and 360 µm

**WHICH PEEK TUBING FOR WHICH APPLICATION?**

Application	Tubing
Standard HPLC*	1/16" OD x 0.25 mm ID
High pressure semi-prep LC*	1/8" OD
Agilent 1100 LC systems**	1/32"
Most capillary systems	360 µm

\* Smaller ID for low flow rates/higher ID for high flow rates. Note that a low linear flow rate (cm/min) causes peak broadening.

\*\* For some 1100 LC systems with high pressure flow paths

**TECH TIP: STRAIGHTENING PEEK TUBING**

1. Start with stainless steel tubing:  
ID slightly larger than OD of PEEK  
Length = 2 cm shorter than PEEK
2. Slide PEEK tubing into stainless tubing, extending 1 cm out each end of the stainless
3. Place sleeved PEEK tubing into oven:  
220°C for 30 minutes, or  
180°C for 60 minutes
4. Allow sleeved tubing to cool in oven to ambient temperature
5. Once tubing is cooled, remove PEEK tubing from the SS sleeve
6. Check for straightness
7. Repeat if necessary



## Striped color-coded PEEK tubing

## SPECS

- Maximum recommended working temp: 100°C (continuous) for 1/16" OD tubing with ID up to 0.75 mm
- Max working pressures: see page 94
- Tubing tolerances: see page 95

We recommend striped PEEK tubing for easy color-coded identification of tubing ID, with maximum chemical resistance and biocompatibility. The stripe is added to the outside, so dye never contacts the fluid stream.

All colorants used in the manufacture of this tubing are RoHS-compliant. (Restriction of Hazardous Substances)

Qty/pkg:			1 m	3 m
Tubing OD	Color	ID (μm)	Product No.	Product No.
360 μm	Orange	50	JR-T-3650-M1	JR-T-3650-M3
	Green	75	JR-T-3675-M1	JR-T-3675-M3
	Grey	100	JR-T-36100-M1	JR-T-36100-M3
	Red	125	JR-T-36125-M1	JR-T-36125-M3
	Black	150	JR-T-36150-M1	JR-T-36150-M3



			1 m	3 m
	Color	ID (mm)	Product No.	Product No.
1/32"	Orange	0.050	JR-T-5702-M1	JR-T-5702-M3
	Pink	0.064	JR-T-5764-M1	JR-T-5764-M3
	Green	0.075	JR-T-5703-M1	JR-T-5703-M3
	Black	0.100	JR-T-5704-M1	JR-T-5704-M3
	Red	0.13	JR-T-5705-M1	JR-T-5705-M3
	Yellow	0.18	JR-T-5707-M1	JR-T-5707-M3
	Orange	0.20	JR-T-5608-M1	JR-T-5608-M3
	Blue	0.25	JR-T-5710-M1	JR-T-5710-M3
	Orange	0.50	JR-T-5720-M1	JR-T-5720-M3

			3 m	5m
	Color	ID (mm)	Product No.	Product No.
1/16"	Orange	0.050	JR-T-5802-M3	JR-T-5802-M5
	Pink	0.064	JR-T-5998C-M3	JR-T-5998C-M5
	Green	0.075	JR-T-5803-M3	JR-T-5803-M5
	Black	0.100	JR-T-5804-M3	JR-T-5804-M5
	Red	0.13	JR-T-5999-M3	JR-T-5999-M5
	Yellow	0.18	JR-T-6000-M3	JR-T-6000-M5
	Blue	0.25	JR-T-6001-M3	JR-T-6001-M5
	Orange	0.50	JR-T-6002-M3	JR-T-6002-M5
	Green	0.75	JR-T-6003-M3	JR-T-6003-M5
	Grey	1.00	JR-T-60031-M3	JR-T-60031-M5
	Black	1.40	JR-T-60032-M3	JR-T-60032-M5

			10 m	25m
	Color	ID (mm)	Product No.	Product No.
	Orange	0.050	JR-T-5802-M10	JR-T-5802-M25
	Pink	0.064	JR-T-5998C-M10	JR-T-5998C-M25
	Green	0.075	JR-T-5803-M10	JR-T-5803-M25
	Black	0.100	JR-T-5804-M10	JR-T-5804-M25
	Red	0.13	JR-T-5999-M10	JR-T-5999-M25
	Yellow	0.18	JR-T-6000-M10	JR-T-6000-M25
	Blue	0.25	JR-T-6001-M10	JR-T-6001-M25
	Orange	0.50	JR-T-6002-M10	JR-T-6002-M25
	Green	0.75	JR-T-6003-M10	JR-T-6003-M25
	Grey	1.00	JR-T-60031-M10	JR-T-60031-M25
	Black	1.40	JR-T-60032-M10	JR-T-60032-M25

## TECH TIP

VICI Jour sample loops for analytical applications are made of PEEK striped color-coded tubing.

PEEK loops ..... 81

## OPTIONS

Contact your distributor or VICI about:

- Longer bulk tubing lengths
- Customized ODs and IDs (30 m minimum quantity)

## MORE INFORMATION

About PEEK® ..... 99

Chemical resistance

PEEK ..... 96-98

PEEK at higher temperatures ..... 96

Maximum working pressure ..... 94

Tubing ID chart: volume per length ... 93

Tubing tolerances ..... 95

## RELATED PRODUCTS

Clean-cut

tubing cutter ..... 86

Easy-Flange kit ..... 85

PEEK fittings

High pressure .... 19-29

Low pressure . 42, 44-47

PEEK tubing elbows ... 86



### Dual layer color-coded PEEK tubing

- Inner layer of natural PEEK, outer layer of color-coded PEEK
- No possibility of stream coming in contact with coloring agents
- Pressure performance equivalent to single layer tubing
- 1/16" OD, IDs from 0.10 to 0.75 mm

Our dual layer PEEK tubing eliminates any concern that a critical sample stream could be contaminated by pigments used to color code the tubing. It looks like any other color-coded tubing at first glance, but a closer look reveals that the pigmented layer\* surrounds a separate but integrally-bonded inner layer of natural PEEK. Only natural PEEK is in contact with the fluid stream.

The ID tolerance of  $\pm 0.025$  mm (.001") is equivalent to or better than that of tubing typically sold as premium grade.

\*All colorants used in the manufacture of this tubing are RoHS-compliant. (Restriction of Hazardous Substances)

#### SPECS

- Maximum recommended working temp: 100°C (continuous)
- Max working pressures: see page 94
- Tubing tolerances: see page 95

Qty/pkg:					3 m	10 m	25 m
Tubing OD	Color	ID (mm)	bar	psi	Product No.	Product No.	Product No.
1/16"	Black	0.10	460	6700	JR-TD-5804-M3	JR-TD-5804-M10	JR-TD-5804-M25
	Red	0.13	420	6100	JR-TD-6007-M3	JR-TD-6007-M10	JR-TD-6007-M25
	Yellow	0.18	400	5800	JR-TD-6008-M3	JR-TD-6008-M10	JR-TD-6008-M25
	Blue	0.25	386	5600	JR-TD-6009-M3	JR-TD-6009-M10	JR-TD-6009-M25
	Orange	0.50	350	4500	JR-TD-6010-M3	JR-TD-6010-M10	JR-TD-6010-M25
	Green	0.75	240	3500	JR-TD-6011-M3	JR-TD-6011-M10	JR-TD-6011-M25

COLORED PEEK  
OUTER LAYER

VIRGIN NATURAL PEEK  
INNER LAYER

#### OPTIONS AND MORE INFORMATION

See facing page.

#### RELATED PRODUCTS

Solid color-coded  
PEEK tubing ..... 8

## Solid color-coded PEEK tubing


## SPECS

- Maximum recommended working temp: 100°C (continuous) for 1/16" OD tubing with ID up to 0.75 mm
- Max working pressures: see page 94
- Tubing tolerances: see page 95

Color-coding identifies the ID of your PEEK tubing – each ID is coded with a different color. Solid colored tubing is made of pre-dyed raw material.

All colorants used in the manufacture of this tubing are RoHS-compliant. (Restriction of Hazardous Substances)

Qty/pkg:			1 m	3 m	10 m
Tubing OD	Color	ID (μm)	Product No.	Product No.	Product No.
360 μm	Orange	50	JR-T-3750-M1	JR-T-3750-M3	—
	Natural	75	—	JR-T-6060-M3	JR-T-6060-M10
	Green	75	JR-T-3775-M1	JR-T-3775-M3	—
	Grey	100	JR-T-37100-M1	JR-T-37100-M3	—
	Red	125	JR-T-37125-M1	JR-T-37125-M3	—
	Black	150	JR-T-37150-M1	JR-T-37150-M3	—
			1 m	3 m	10 m
	Color	ID (mm)	Product No.	Product No.	Product No.
1/32"	Orange	0.050	JR-T-5602-M1	JR-T-5602-M3	—
	Pink	0.064	JR-T-5664-M1	JR-T-5664-M3	—
	Green	0.075	JR-T-5603-M1	JR-T-5603-M3	—
	Black	0.100	JR-T-5604-M1	JR-T-5604-M3	—
	Natural	0.13	—	JR-T-5993-M3	JR-T-5993-M10
	Red	0.13	JR-T-5605-M1	JR-T-5605-M3	—
	Purple	0.15	JR-T-5606-M1	JR-T-5606-M3	—
	Natural	0.18	—	JR-T-5994-M3	JR-T-5994-M10
	Yellow	0.18	JR-T-5607-M1	JR-T-5607-M3	—
	Orange	0.20	JR-T-5608-M1	JR-T-5608-M3	—
	Natural	0.25	—	JR-T-5995-M3	JR-T-5995-M10
	Blue	0.25	JR-T-5610-M1	JR-T-5610-M3	—
	Natural	0.38	—	JR-T-5715-M3	JR-T-5715-M10
	Orange	0.50	JR-T-5620-M1	JR-T-5620-M3	—
			3 m	10 m	25 m
	Color	ID (mm)	Product No.	Product No.	Product No.
1/16"	Natural	0.064	JR-T-5998-M3	JR-T-5998-M10	JR-T-5998-M25
	Black	0.100	JR-T-6020-M3	JR-T-6020-M10	JR-T-6020-M25
	Red	0.13	JR-T-6007-M3	JR-T-6007-M10	JR-T-6007-M25
	Purple	0.15	JR-T-6030-M3	JR-T-6030-M10	JR-T-6030-M25
	Yellow	0.18	JR-T-6008-M3	JR-T-6008-M10	JR-T-6008-M25
	Blue	0.25	JR-T-6009-M3	JR-T-6009-M10	JR-T-6009-M25
	Grey	0.38	JR-T-6040-M3	JR-T-6040-M10	JR-T-6040-M25
	Orange	0.50	JR-T-6010-M3	JR-T-6010-M10	JR-T-6010-M25
	Green	0.75	JR-T-6011-M3	JR-T-6011-M10	JR-T-6011-M25

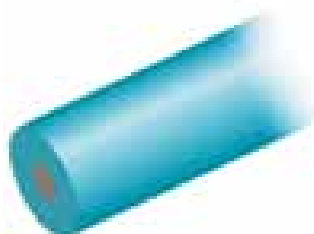

## TECH TIP

See our dual layer PEEK tubing for solid colored PEEK tubing with no possibility of stream coming into contact with coloring agents. . . . . 7

## OPTIONS

Contact your distributor or VICI about:


- Longer bulk tubing lengths
- Customized ODs and IDs (30 m minimum quantity)

## MORE INFORMATION

About PEEK® . . . . . 99  
 Chemical resistance  
   PEEK . . . . . 96-98  
   PEEK at higher temperatures . . . . . 96  
 Maximum working pressure . . . . . 94  
 Tubing ID chart:  
   volume per length . . . 93  
 Tubing tolerances . . . . 95

## RELATED PRODUCTS

Clean-cut  
   tubing cutter . . . . . 86  
 Easy-Flange kit . . . . . 85  
 PEEK fittings  
   High pressure . . . 19-29  
   Low pressure . . . 42, 44-47  
 PEEK tubing elbows . . . 86

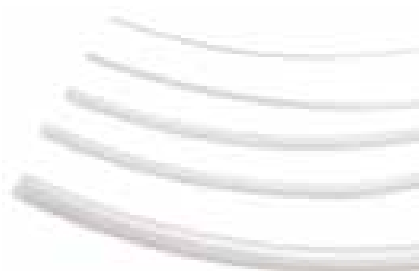

## ETFE tubing

ETFE (ethyltrifluoroethylene) tubing is particularly suitable for HPLC, since it has a higher burst pressure as compared to PTFE. 1/16" OD x 0.25 mm ID tubing will withstand up to 186 bar (2700 psi) with water as mobile phase. However, it is important to note that when used with organic solvents, ETFE tubing may swell, reducing that maximum pressure to 70 bar (1000 psi).

This is the ideal polymer for applications with chloride salts, which typically corrode stainless steel.

### SPECS

- Max working pressures: see page 94
- Tubing tolerances: see page 95

Qty/pkg:		3 m	10 m	25 m
Tubing OD	ID (mm)	Product No.	Product No.	Product No.
1/32"	0.25	JR-T-084-M3	JR-T-084-M10	JR-T-084-M25
1/16"	0.18	JR-T-078-M3	JR-T-078-M10	JR-T-078-M25
	0.25	JR-T-080-M3	JR-T-080-M10	JR-T-080-M25
	0.50	JR-T-082-M3	JR-T-082-M10	JR-T-082-M25
	0.75	JR-T-083-M3	JR-T-083-M10	JR-T-083-M25
1/8"	1.00	JR-T-085-M3	JR-T-085-M10	JR-T-085-M25
	1.59	JR-T-086-M3	JR-T-086-M10	JR-T-086-M25
		1 m	3 m	
	ID (mm)	Product No.	Product No.	
3.0 mm	1.0	JR-T-087-M1	JR-T-087-M3	
6.0 mm	4.0	JR-T-088-M1	JR-T-088-M3	

### TECH TIP

- When selecting tubing for low pressure applications, give careful consideration to gas permeability characteristics.
- Use PFA instead of PTFE if gas permeability is an issue.

### OPTIONS

Contact your distributor or VICI about:

- Longer tubing lengths at competitive bulk prices
- Customized ODs and IDs (30 m minimum quantity)

### MORE INFORMATION

About ETFE. . . . . 99  
 Chemical resistance 96-98  
 Maximum working pressure . . . . . 94  
 Tubing ID chart: volume per length . . . 93  
 Tubing tolerances . . . . 95

### RELATED PRODUCTS

Clean-cut tubing cutter . . . . . 86  
 Easy-Flange kit . . . . . 85  
 Low pressure fittings . . . . . 40-47  
 SolviFlex coiled PFA tubing . . . . . 70

## PFA tubing


PFA (perfluoroalkoxy) tubing from VICI Jour has excellent chemical stability and mechanical strength. 1/16" OD high purity tubing has a maximum pressure rating of up to 100 bar (1450 psi), while 1/8" OD tubing can be used at up to 72 bar (1050 psi).

Always use PFA instead of PTFE if gas permeability is an issue.

### SPECS

- Max working pressures: see page 94
- Tubing tolerances: see page 95

Qty/pkg:		1 m	3 m	
Tubing OD	ID (mm)	Product No.	Product No.	
360 µm	50 µm	JR-T-43650-M1	JR-T-43650-M3	
	75 µm	JR-T-43675-M1	JR-T-43675-M3	
	100 µm	JR-T-436100-M1	JR-T-436100-M3	
		3 m	10 m	25 m
	ID (mm)	Product No.	Product No.	Product No.
1/16"	0.50 mm	JR-T-4001-M3	JR-T-4001-M10	JR-T-4001-M25
	0.75 mm	JR-T-4002-M3	JR-T-4002-M10	JR-T-4002-M25
	1.00 mm	JR-T-4007-M3	JR-T-4007-M10	JR-T-4007-M25
1/8"	1.59 mm	JR-T-4003-M3	JR-T-4003-M10	JR-T-4003-M25
		per meter		
	ID (mm)	Product No.		
1/8"	1.70 mm	JR-T-4005		
3/16"	3.20 mm	JR-T-4009		
1/4"	4.35 mm	JR-T-4004		

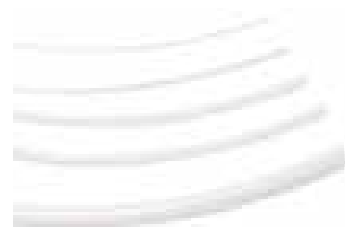

## FEP tubing, solid color-coded

**SPECS**

- Maximum recommended working temp: 80°C (continuous)
- Max working pressures: see page 94
- Tubing tolerances: see page 95

FEP (FluorinatedEthylenePropylene) tubing is chemically inert to most HPLC solvents, and is recommended for low pressure applications. Typically 1/16" OD tubing holds up to 117 bar (1700 psi), and 1/8" OD holds up to 70 bar (1000 psi).

Solid colors identify the tubing ID.

We recommend FEP for ion chromatography applications.



Qty/pkg:		3 m		10 m		25 m	
Tubing OD	Color	ID (mm)	Product No.	Product No.	Product No.	Product No.	Product No.
1/32"	Black	0.10	JR-T-6604C-M3	JR-T-6604C-M10	JR-T-6604C-M25	JR-T-6604C-M25	JR-T-6604C-M25
	Red	0.13	JR-T-6605C-M3	JR-T-6605C-M10	JR-T-6605C-M25	JR-T-6605C-M25	JR-T-6605C-M25
	Purple	0.15	JR-T-6606C-M3	JR-T-6606C-M10	JR-T-6606C-M25	JR-T-6606C-M25	JR-T-6606C-M25
	Yellow	0.18	JR-T-6607C-M3	JR-T-6607C-M10	JR-T-6607C-M25	JR-T-6607C-M25	JR-T-6607C-M25
	Blue	0.25	JR-T-6610C-M3	JR-T-6610C-M10	JR-T-6610C-M25	JR-T-6610C-M25	JR-T-6610C-M25
	Orange	0.50	JR-T-6620C-M3	JR-T-6620C-M10	JR-T-6620C-M25	JR-T-6620C-M25	JR-T-6620C-M25
		3 m		10 m		25 m	
	Color	ID (mm)	Product No.	Product No.	Product No.	Product No.	Product No.
1/16"	Black	0.10	JR-T-6704C-M3	JR-T-6704C-M10	JR-T-6704C-M25	JR-T-6704C-M25	JR-T-6704C-M25
	Red	0.13	JR-T-6705C-M3	JR-T-6705C-M10	JR-T-6705C-M25	JR-T-6705C-M25	JR-T-6705C-M25
	Purple	0.15	JR-T-6706C-M3	JR-T-6706C-M10	JR-T-6706C-M25	JR-T-6706C-M25	JR-T-6706C-M25
	Yellow	0.18	JR-T-6707C-M3	JR-T-6707C-M10	JR-T-6707C-M25	JR-T-6707C-M25	JR-T-6707C-M25
	Blue	0.25	JR-T-6812C-M3	JR-T-6812C-M10	JR-T-6812C-M25	JR-T-6812C-M25	JR-T-6812C-M25
	Orange	0.50	JR-T-6802C-M3	JR-T-6802C-M10	JR-T-6802C-M25	JR-T-6802C-M25	JR-T-6802C-M25
	Green	0.75	JR-T-6803C-M3	JR-T-6803C-M10	JR-T-6803C-M25	JR-T-6803C-M25	JR-T-6803C-M25
	Black	1.00	JR-T-6806C-M3	JR-T-6806C-M10	JR-T-6806C-M25	JR-T-6806C-M25	JR-T-6806C-M25

**TECH TIP**

When selecting tubing for low pressure applications, give careful consideration to gas permeability characteristics.

**OPTIONS**

Contact your distributor or VICI about:

- Longer tubing lengths at competitive bulk prices
- Customized ODs and IDs (30 m minimum quantity)

**MORE INFORMATION**

About FEP . . . . . 99  
 Chemical resistance 96-98  
 Maximum working pressure . . . . . 94  
 Tubing ID chart: volume per length . . . 93  
 Tubing tolerances . . . . . 95

**RELATED PRODUCTS**

Clean-cut tubing cutter . . . . . 86  
 Easy-Flange kit . . . . . 85  
 Low pressure fittings . . . . . 40-47

**FEP tubing, natural color**

FEP (FluorinatedEthylenePropylene) tubing is chemically inert to most HPLC solvents, and is recommended for low pressure applications. Typically 1/16" OD tubing holds up to 117 bar (1700 psi), and 1/8" OD holds up to 70 bar (1000 psi).

We recommend FEP for ion chromatography applications.

**SPECS**

- Maximum recommended working temp: 80°C (continuous)
- Max working pressures: see page 94
- Tubing tolerances: see page 95

Qty/pkg:		1 m	3 m	10 m
Tubing OD	ID (mm)	Product No.	Product No.	Product No.
1/32"	0.10	JR-T-6794-M1	JR-T-6794-M3	JR-T-6794-M10
	0.13	JR-T-6605-M1	JR-T-6605-M3	JR-T-6605-M10
	0.15	JR-T-6606-M1	JR-T-6606-M3	JR-T-6606-M10
	0.18	JR-T-6607-M1	JR-T-6607-M3	JR-T-6607-M10
	0.25	JR-T-6610-M1	JR-T-6610-M3	JR-T-6610-M10
	0.38	JR-T-6615-M1	JR-T-6615-M3	JR-T-6615-M10
	0.50	JR-T-6620-M1	JR-T-6620-M3	JR-T-6620-M10
		3 m	10 m	25 m
	ID (mm)	Product No.	Product No.	Product No.
1/16"	0.10	JR-T-6704-M3	JR-T-6704-M10	JR-T-6704-M25
	0.13	JR-T-6705-M3	JR-T-6705-M10	JR-T-6705-M25
	0.15	JR-T-6706-M3	JR-T-6706-M10	JR-T-6706-M25
	0.18	JR-T-6707-M3	JR-T-6707-M10	JR-T-6707-M25
	0.20	JR-T-6708-M3	JR-T-6708-M10	JR-T-6708-M25
	0.25	JR-T-6812-M3	JR-T-6812-M10	JR-T-6812-M25
	0.30	JR-T-6702-M3	JR-T-6702-M10	JR-T-6702-M25
	0.50	JR-T-6802-M3	JR-T-6802-M10	JR-T-6802-M25
1/8"	0.75	JR-T-6803-M3	JR-T-6803-M10	JR-T-6803-M25
	1.00	JR-T-6806-M3	JR-T-6806-M10	JR-T-6806-M25
1/8"	1.59	JR-T-6804-M3	JR-T-6804-M10	JR-T-6804-M25
	2.10	JR-T-6813-M3	JR-T-6813-M10	JR-T-6813-M25
3/16"	3.18	JR-T-6835-M3	JR-T-6835-M10	JR-T-6835-M25
1 mm	0.50	JR-T-6105-M3	JR-T-6105-M10	JR-T-6105-M25
2 mm	1.00	JR-T-6210-M3	JR-T-6210-M10	JR-T-6210-M25
3 mm	1.00	JR-T-6310-M3	JR-T-6310-M10	JR-T-6310-M25
	2.00	JR-T-6320-M3	JR-T-6320-M10	JR-T-6320-M25
Qty/pkg:		per meter		
Tubing OD	ID (mm)	Product No.		
1/4"	4.35	JR-T-6809		

**OPTIONS AND MORE INFORMATION**

See facing page.

## TUBING

## PTFE tubing, solid color-coded

## SPECS

- Maximum recommended working temp: 100°C (continuous)
- Maximum working pressures: see page 94
- Tubing tolerances: see page 95

PTFE (polytetrafluoroethylene) tubing is chemically inert and suitable for lower pressure applications. Typically 1/16" OD tubing holds up to 62 bar (900 psi) and 1/8" up to 35 bar (500 psi). Maximum continuous operating temperature is 100°C.

Black PTFE, available in two IDs, is recommended for light sensitive applications such as enzymatic reactions, luminescent reactions, connecting to a flow cell, and radioactive HPLC.

Solid colors identify the tubing ID.

Qty/pkg:		3 m		10 m	25 m
OD	Color	ID (mm)	Product No.	Product No.	Product No.
1/16"	Blue	0.25	JR-T-4011C-M3	JR-T-4011C-M10	JR-T-4011C-M25
	Black	0.25	JR-T-6808-M3	JR-T-6808-M10	JR-T-6808-M25
		0.50	JR-T-6811-M3	JR-T-6811-M10	JR-T-6811-M25
	Orange	0.50	JR-T-4183C-M3	JR-T-4183C-M10	JR-T-4183C-M25
	Green	0.75	JR-T-4136C-M3	JR-T-4136C-M10	JR-T-4136C-M25



## PTFE tubing, natural color

## SPECS

- Maximum recommended working temp: 100°C (continuous)
- Max working pressures: see page 94
- Tubing tolerances: see page 95

PTFE tubing is chemically inert and suitable for lower pressure applications. Typically 1/16" OD tubing holds up to 62 bar (900 psi) and 1/8" up to 35 bar (500 psi). Maximum continuous operating temperature is 100°C.

Qty/pkg:		3 m		5 m	10 m	25 m
Tubing OD	ID (mm)	Product No.	Product No.	Product No.	Product No.	Product No.
1/16"	0.18	JR-T-6805-M3	JR-T-6805-M5	JR-T-6805-M10	JR-T-6805-M25	
	0.25	JR-T-4011-M3	JR-T-4011-M5	JR-T-4011-M10	JR-T-4011-M25	
	0.50	JR-T-4183-M3	JR-T-4183-M5	JR-T-4183-M10	JR-T-4183-M25	
	0.75	JR-T-4036-M3	JR-T-4036-M5	JR-T-4036-M10	JR-T-4036-M25	
	1.00	JR-T-6807-M3	JR-T-6807-M5	JR-T-6807-M10	JR-T-6807-M25	
1/8"	1.59	JR-T-6800-M3	JR-T-6800-M5	JR-T-6800-M10	JR-T-6800-M25	
	2.40	JR-T-4037-M3	JR-T-4037-M5	JR-T-4037-M10	JR-T-4037-M25	
		3 m	10 m	25 m		
	ID (mm)	Product No.	Product No.	Product No.		
2 mm	1.70	JR-T-6801-M3	JR-T-6801-M10	JR-T-6801-M25		
4 mm	3.00	JR-T-4039-M3	JR-T-4039-M10	—		
Qty/pkg:		per meter				
Tubing OD	ID (mm)	Product No.				
in.	mm					
3/16"	4.76	3.76	JR-T-4041			
1/4"	6.35	4.35	JR-T-6810			



## TECH TIP

Use PFA instead of PTFE if gas permeability is an issue.

## TECH TIP

Use black PTFE for light sensitive applications such as enzymatic reactions, luminescent reactions, connection to a flow cell, and radioactive HPLC.

## OPTIONS

Contact your distributor or VICI about:


- Longer tubing lengths at competitive bulk prices
- Customized ODs and IDs (30 m minimum quantity)

## MORE INFORMATION

About PTFE . . . . . 99  
 Chemical resistance 96-98  
 Maximum working pressure . . . . . 94  
 Tubing ID chart:  
   volume per length . . . 91  
 Tubing tolerances . . . . 95

## RELATED PRODUCTS

Clean-cut  
 tubing cutter . . . . . 86  
 Easy-Flange kit . . . . . 85  
 Flanged fittings . . . . 40, 43  
 Flangeless fittings . . 41-42

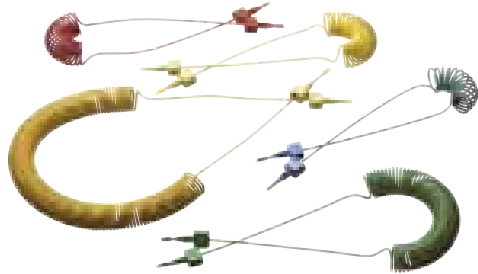

### Thermo-formed polymeric tubing


- Available in a variety of polymers
- Provides axial flexibility
- Telephone cables or customer-specified shapes

VICI Jour offers thermo-formed tubing in various geometric shapes and polymers. The form will maintain its shape so long as the temperature does not exceed 100°C.

Straight tubing, spirals, or even more complex forms – just define your form and material on a drawing for us.

#### FOR PRICING

Contact VICI or your local distributor for more information or a quotation.

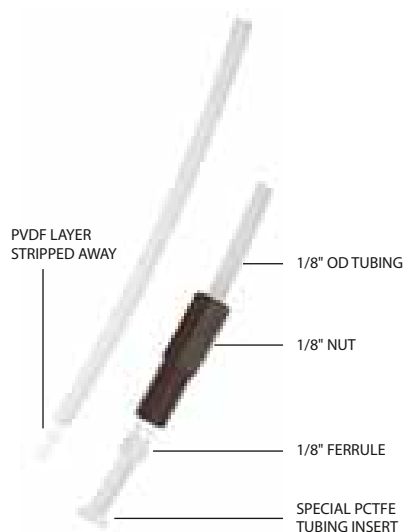

### Impermeable polymeric tubing

- FEP coated with PVDF
- Prevents “re-gassing”
- Excellent chemical resistance

PTFE and FEP have very good chemical resistance, but are too permeable to be suitable for mobile phase transfer lines. PVDF, on the other hand, lacks chemical resistance but is a very tight polymer with no observable micropores. Our impermeable tubing is produced by co-extruding PVDF over FEP to combine the best features of both polymers, ensuring that degassed solvent will stay degassed until it reaches the HPLC pump.

The PVDF layer can easily be stripped away with our Clean-Cut tubing cutter or a razor blade, leaving a length of exposed FEP to extend into the solvent reservoir.

When making connections, a special PCTFE tubing insert fits into the ID of the tubing and seals against the bottom of the port, preventing the solvent from contacting the PVDF layer at the end of the tube.

#### SPECS

- Low pressure only
- 1/8" OD, 1.65 mm ID
- OD and ID tolerances:  $\pm 0.05$  mm (.002")
- Materials  
Tubing: FEP coated with PVDF  
Insert: PCTFE  
Ferrule: ETFE  
Fitting: PPS

#### TECH TIP

Use FEP or PFA tubing inside the solvent reservoir, since PVDF has limited chemical resistance.

#### TECH TIP

Always use the specially-designed tubing insert to prevent solvent contact with the PVDF layer.

#### OPTIONS

Contact us for more info about longer tubing lengths at competitive bulk prices.


#### MORE INFORMATION

About FEP ..... 99  
About PVDF ..... 100  
Chemical resistance 96-98  
Tubing ID chart:  
volume per length... 93  
Tubing tolerances .... 95

#### RELATED PRODUCTS

Clean-cut  
tubing cutter ..... 86  
Color-sleeve  
fingertight adapters . 17

Qty/pkg:		
	3 m	10 m
	Product No.	Product No.
Impermeable tubing, 1/8" OD x 1.65 mm ID	JR-T-6130-M3	JR-T-6130-M10
Additional spare parts	Qty/pkg	Product No.
PCTFE tubing inserts	10	JR-6141-10
1/8" PPS nuts	10	JR-55071-10
1/8" ETFE ferrules	10	JR-051-10
Fitting kit		JR-6140
contains:		
PCTFE tubing inserts	2	
1/8" PPS nuts	2	
1/8" ETFE ferrules	2	

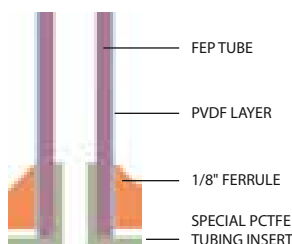



## TUBING

## Pre-cut stainless steel tubing

UHPLC

## SPECS

- Pressure rating varies with ID; minimum is 365 bar (5300 psi) (see page 93)
- Tolerances
 

OD	±0.05mm (.002")
ID	±0.025mm (.001")

- For zero dead volume connections
- Steam-cleaned and ready to use
- Soft annealed OD for easy ferrule connections
- Color bands for quick identification of tubing ID

VICI Jour offers pre-cut 1/16" stainless steel tubing in lengths most commonly required in UHPLC and HPLC systems. The tubes are cut to square, burr-free ends for zero dead volume connections. Colored bands allow quick ID identification.

VICI pre-cut tubing is specially cleaned with steam from micro-filtered deionized water, which removes both organic and inorganic contaminants. Longer tubes are cleaned with liquid hexane. Softener-free polyethylene caps are used to prevent any contamination during shipment and storage.

Length:		50 mm		100 mm	200 mm	300 mm
OD	Color	ID (mm)	Product No.	Product No.	Product No.	Product No.
1/16"	Red	0.13	JR-T-97005	JR-T-97010	JR-T-97015	JR-T-97020
	Yellow	0.18	JR-T-97025	JR-T-97030	JR-T-97035	JR-T-97040
	Blue	0.25	JR-T-97045	JR-T-97055	JR-T-97060	JR-T-97065
	Orange	0.50	JR-T-97075	JR-T-97080	JR-T-97085	JR-T-97090
	Green	0.75	JR-T-97095	JR-T-97100	JR-T-97105	JR-T-97110



## Pre-cut flexible stainless steel tubing

UHPLC

## SPECS

- Pressure rating varies with ID; minimum is 365 bar (5300 psi) (see page 93)
- Tolerances
 

OD	±0.05mm (.002")
ID	±0.025mm (.001")

- Seamless 1/32" tubing with 1/16" laser-welded ends
- For zero dead volume connections
- Steam-cleaned and ready to use
- Soft annealed OD for easy ferrule connections
- Color bands for quick identification

Our flexible pre-cut stainless steel tubing is made of 1/32" premium grade 316 stainless steel tubing with ends laser-welded into 1/16" tubing of the same premium grade. Available in lengths most commonly required in UHPLC and HPLC systems.




The tubes are cut to square, burr-free ends for zero dead volume connections. Colored bands allow quick ID identification.

VICI pre-cut tubing is specially cleaned with steam from micro-filtered deionized water, which removes both organic and inorganic contaminants. Longer tubes that cannot be steam-cleaned are cleaned with liquid hexane. Softener-free polyethylene caps are used to prevent any contamination during shipment and storage.

Length:		105 mm		150 mm	200 mm	280 mm	400 mm
OD	Color	ID (mm)	Product No.	Product No.	Product No.	Product No.	Product No.
1/16"	Red	0.13	JR-T-96005	JR-T-96010	JR-T-96015	JR-T-96020	JR-T-96025
	Yellow	0.18	JR-T-96105	JR-T-96110	JR-T-96115	JR-T-96120	JR-T-96125
	Blue	0.25	JR-T-96205	JR-T-96210	JR-T-96215	JR-T-96220	JR-T-96225
	Orange	0.50	JR-T-96305	JR-T-96310	JR-T-96315	JR-T-96320	JR-T-96325
		600 mm		700 mm	800 mm	900 mm	
		ID (mm)	Product No.	Product No.	Product No.	Product No.	
1/16"	Red	0.13	JR-T-96030	JR-T-96035	JR-T-96040	—	
	Yellow	0.18	JR-T-96130	JR-T-96135	JR-T-96140	JR-T-96145	
	Blue	0.25	JR-T-96230	JR-T-96235	JR-T-96240	—	
	Orange	0.50	JR-T-96330	JR-T-96335	JR-T-96340	—	

## TECH TIP

Contact your local distributor or VICI for custom lengths.

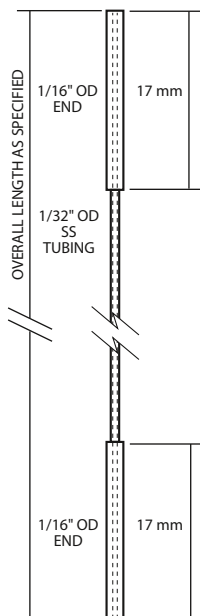

## Stainless steel tubing – premium grade

UHPLC

- Seamless Type 316 stainless steel
- Variety of sizes for UHPLC, HPLC, and GC applications
- Smooth internal surface
- Soft annealed OD for easy ferrule connections

Our premium grade seamless stainless steel tubing is the base for all high pressure and ultra-high pressure lines in your equipment.

Note that this tubing is roughly cut with tubing pliers so the ends must be properly cut before use.

Qty/pkg:		3 m	10 m
Tubing OD	ID (mm)	Product No.	Product No.
1/32"	0.13	JR-TSS.505-M3	JR-TSS.505-M10
	0.18	JR-TSS.507-M3	JR-TSS.507-M10
	0.25	JR-TSS.510-M3	JR-TSS.510-M10
	0.50	JR-TSS.520-M3	JR-TSS.520-M10
1/16"	0.13	JR-TSS105-M3	JR-TSS105-M10
	0.25	JR-TSS110-M3	JR-TSS110-M10
	0.38	JR-TSS115-M3	JR-TSS115-M10
	0.50	JR-TSS120-M3	JR-TSS120-M10
	0.75	JR-TSS130-M3	JR-TSS130-M10
	1.00	JR-TSS140-M3	JR-TSS140-M10
1/8"	0.75	JR-TSS230-M3	JR-TSS230-M10
	1.00	JR-TSS240-M3	JR-TSS240-M10
	1.52	JR-TSS260-M3	JR-TSS260-M10
	1.70	JR-TSS267-M3	JR-TSS267-M10
	2.16	JR-TSS285-M3	JR-TSS285-M10

## SPECS

- Pressure rating (values shown are minimum; rating varies with ID; (see page 93))
 

OD	bar	psi
1/32"	400	5800
1/16"	365	5300
1/8"	315	4600
- Tolerances
 

OD	±0.05mm (.002")
ID	±0.025mm (.001")

## TECH TIP

For zero dead volume connections, we strongly recommend pre-cut stainless steel tubing . . . 14

## NOTE

ID dimensions in mm are converted from nominal inch dimensions.

## OPTIONS

Contact us for more info about longer tubing lengths at competitive bulk prices.

## MORE INFORMATION

About Type 316 stainless steel . . . . . 100  
 Maximum working pressure . . . . . 93  
 Tubing ID chart: volume per length . . . 93

## RELATED PRODUCTS

Stainless steel tubing cutter . . . . . 87  
 Stainless fittings . . . . . 30-35

## Stainless steel tubing – standard grade

UHPLC

- Type 316 stainless steel, welded under inert gas
- Variety of sizes for HPLC and GC applications
- Smooth internal surface
- Soft annealed OD for easy ferrule connections

Qty/pkg:		3 m	10 m
Tubing OD	ID (mm)	Product No.	Product No.
1/16"	0.13	JR-T-625-04-M3	JR-T-625-04-M10
	0.18	JR-T-625-05-M3	JR-T-625-05-M10
	0.25	JR-T-625-10-M3	JR-T-625-10-M10
	0.50	JR-T-625-20-M3	JR-T-625-20-M10
	0.75	JR-T-625-30-M3	JR-T-625-30-M10
	1.00	JR-T-625-40-M3	JR-T-625-40-M10
1/8"	2.1	JR-T-626-00-M3	JR-T-626-00-M10
Qty/pkg:		per meter	
Tubing OD	ID (mm)	Product No.	
1/4"	4.65	JR-T-628-00	

## SPECS

- Pressure rating (values shown are minimum; rating varies with ID; (see page 93))
 

OD	bar	psi
1/16"	315	4600
1/8"	315	4600
- Tolerances
 

OD	±0.05mm (.002")
ID	±0.05mm (.002")


# Fittings

## HIGH AND LOW PRESSURE




### High pressure polymeric fittings... page 18

High pressure polymeric fittings are useful in applications involving proteins, peptides, and other samples of biological origin where metal systems may interact with samples.

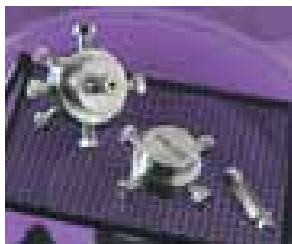


### High pressure stainless fittings.... page 30

Our stainless steel fittings are typically compression fittings, in which a ferrule is compressed onto the tube as a nut is tightened. They offer the best reliability in high pressure situations.

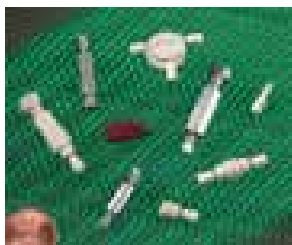

### Nanovolume® HPLC fittings ..... page 36

With bore sizes of 100 and 150 microns, (.004" -.006"), the minimal transfer volume contributed by Nanovolume® components makes them especially beneficial in applications with flow rates in the µl/min range, when the transfer volume can be critical.

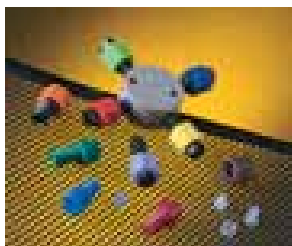

### Low pressure fittings..... page 40

Low pressure fittings are ideally suited for flow injection analysis, low pressure liquid chromatography, and stream sampling devices. They may be safely used at pressures up to 500 psi and temperatures to 50°C.

#### MORE INFORMATION

##### Tools for fittings

Ferrule removal kit.....	89
Tweezer.....	89
ValvTool.....	88
Wrench.....	89






Color-It fingertight adapters

- For high and low pressure
- Fit all hex-head PEEK and SS fittings for 1/16" tubing
- Available in six colors for easy solvent line identification

**SPECS**  
■ For 1/4" hex-head nuts  
■ Material: Polypropylene

Color-It snap-on fingertight adapters are used with for our 1/4" PEEK and stainless hex-head nuts, converting them into color-coded fingertight fittings in an instant. Available in six colors.

Color	Qty/pkg	Product No.
Blue	5	JR-55010-5
Yellow	5	JR-55011-5
Green	5	JR-55012-5
Black	5	JR-55013-5
White	5	JR-55014-5
Red	5	JR-55015-5
Assorted	2 each	JR-55016-12

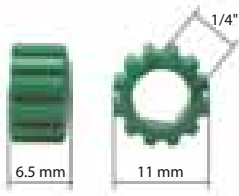


Color-Sleeve fingertight adapters

- For high and low pressure
- Designed for our PPS fittings for 1/16" tubing
- Available in twelve colors for easy solvent line identification

**SPECS**  
■ For 1/4" hex-head nuts  
■ Material: Polypropylene

This adapter is useful for color coding different solvent lines, and turns VICI Jour PPS nuts into fingertight fittings in an instant.

Color	Qty/pkg	Product No.	Color	Qty/pkg	Product No.
 Natural	5 10	JR-55080-5 JR-55080-10	 Brown	5 10	JR-55088-5 JR-55088-10
 White	5 10	JR-55081-5 JR-55081-10	 Green	5 10	JR-55089-5 JR-55089-10
 Dark grey	5 10	JR-55082-5 JR-55082-10	 Blue	5 10	JR-55090-5 JR-55090-10
 Black	5 10	JR-55083-5 JR-55083-10	 Purple	5 10	JR-55091-5 JR-55091-10
 Lavender	5 10	JR-55084-5 JR-55084-10	Assorted		
 Red	5 10	JR-55085-5 JR-55085-10	All 12 colors	2 each	JR-5508X-24
 Yellow	5 10	JR-55086-5 JR-55086-10	Blue, red, green, yellow, black, white	2 each	JR-5508X-12
 Orange	5 10	JR-55087-5 JR-55087-10			

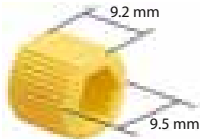

**MORE INFORMATION**  
Hex head nuts  
PEEK ..... 19-23  
Stainless..... 31  
PPS  
High pressure... 19, 25  
Low pressure ..... 41

# High pressure fittings

POLYMERIC

Our high pressure fittings are rated up to 350 bar (5000 psi) with fingertight nuts, well beyond the burst strength of most PEEK tubing. These fittings are manufactured from high quality inert polymers to exacting tolerances.

## High pressure polymeric fittings index

### One-piece fittings

PEEK hex-head	
Over 5,000 psi	<b>NEW</b> ... page 19
5,000 psi	19
PEEK fingertight	
Color-coded	20
Natural	20
PPS hex-head	19
PEEK no-twist one piece fittings	21

### Nuts and ferrules

PEEK ferrules, over 5,000 psi	<b>NEW</b> ... 22
Single ferrule connections	22
Double ferrule connections	23
1/4-28 connections	25

### Plugs and caps

PEEK plugs	24
PEEK caps	24

### PEEK unions

Unions	26
Bulkhead unions	26
Reducing unions	27

### PEEK tees and crosses

Standard tees and crosses	28
Biocompatible mixing tee	28

### PEEK adapters

PEEK port and Luer adapters	29
-----------------------------	----

### DISCLAIMER

The maximum holding pressure for any type of connections involving tubing and a ferrule varies considerably with the tubing material, the ferrule material, the clearance between tubing OD and ferrule ID, and the shape of the fitting detail.

### RELATED PRODUCTS

ChromBox kit	91
HPLC parts kit with PEEK and stainless fittings and tubing	
Fingertight adapters	
Color-It	17
Color-Sleeve	17
Nanovolume® high pressure fittings	36-39
PEEK starter kits	90
Tools for fittings	88-89



**PEEK one-piece hex-head fittings****NEW****UHPLC**

- Made of special PEEK™ HT for better pressure resistance
- One-piece – no ferrule required
- Biocompatible
- Handy tightening tool included

These fittings are machined from special PEEK material to reach a higher pressure rating than standard PEEK fittings.

Sold in packages of 10.

Includes a tightening tool (*page 89*) for easy installation.

Tubing ID	Thread	Product No.
1/32"	6-40	JR-5525-10
1/16"	10-32	JR-5520-10

**SPECS**

- Material:  
PEEK, special black
- Threads:  
1/32" 6-40  
1/16" 10-32
- Pressure ratings, fingertight, on SS tubing:  
1/32" 625 bar (9,000 psi)  
1/16" 700 bar (10,000 psi)



**PEEK one-piece hex-head fittings****10-32**

- Fingertight up to 350 bar (5000 psi)
- Excellent chemical and mechanical properties
- One-piece – no ferrule required
- Biocompatible
- For any 10-32 high pressure fitting

PEEK one-piece hex-head fittings are easy to use with the ValvTool or with the Color-It color-coded fingertight adapters. The 10-32 thread will fit virtually all standard high pressure female connections, and will withstand pressures of up to 350 bar (5000 psi). Available in short and long versions.

Sold in packages of 5.

Description	Product No.
Short	JR-55100-5
Long	JR-55110-5

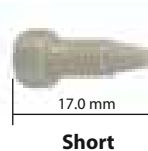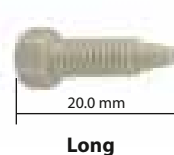**Short****Long****SPECS**

- For 1/16" OD tubing
- Material:  
PEEK
- 10-32 threads
- Tolerance:  
±0.05 mm (.002")
- Pressure rating:  
<350 bar (5000 psi)  
(Varies with tubing material and ID)

**TECH TIP**

Coloration pigments are carefully chosen so as to contain no heavy metals or hazardous chemicals.

**MORE INFORMATION**

About PEEK ..... 99  
About PPS..... 99  
Chemical resistance 96-98

**RELATED PRODUCTS**

1/16" OD tubing for high pressure applications ..... 5-8  
PEEK starter kit ..... 90  
Color-Sleeve adapters . 17



1/32" and 1/16" tightening tools for PEEK fittings. .... 89

**PPS one-piece hex-head fittings****10-32**

- Economically-priced alternative to PEEK
- Excellent chemical and mechanical properties
- One-piece – no ferrule required
- Slim design – finger-tight with Color-Sleeve adapter


The VICI Jour PPS one-piece hex-head (fingertight with adapter) fitting offers excellent chemical and pressure resistance.

Sold in packages of 5.

Product No.
JR-58000-5



Shown with Color-Sleeve adapter, sold separately

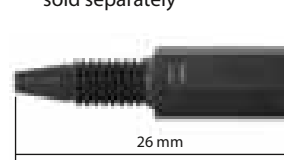**SPECS**

- For 1/16" OD tubing
- Material:  
Black PPS (Polyphenylene sulfide)
- 10-32 threads
- Tolerance:  
±0.05 mm (.002")
- Pressure rating:  
<350 bar (5000 psi)  
(Varies with tubing material and ID)

## PEEK one-piece fingertight fittings – color-coded

10-32

**SPECS**


- For 1/16" OD tubing
- Material: PEEK
- 10-32 threads
- Tolerance:  $\pm 0.05$  mm (.002")
- Pressure rating: <350 bar (5000 psi)  
(Varies with tubing material and ID)

- Fingertight up to 350 bar (5000 psi)
- Excellent chemical and mechanical properties
- One-piece – no ferrule required
- Biocompatible

Our PEEK fingertight fittings are available in six colors, for easy identification of tubing lines. They are rated to 350 bar (5000 psi) and can be used in virtually any type of 10-32 thread HPLC fitting detail on the market.

Sold in packages of 5.

Color	Product No.
Natural	JR-55020-5
Black	JR-55021-5
Red	JR-55022-5
Yellow	JR-55023-5
Blue	JR-55024-5
Green	JR-55025-5

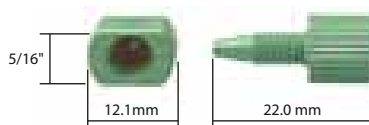

## PEEK one-piece fingertight fittings

10-32


**SPECS**

- For 1/16" OD tubing
- Material: PEEK
- 10-32 threads
- Tolerance:  $\pm 0.05$  mm (.002")
- Pressure rating: <350 bar (5000 psi)  
(Varies with tubing material and ID)

- Natural PEEK versions fingertight up to 350 bar (5000 psi)
- One-piece – no ferrule required
- Biocompatible
- Four styles, machined for extra precision


Sold in packages of 5.

Description	Product No.
Standard	JR-5502-5
Short	JR-5504-5
Long	JR-5507-5
Narrow hex-head/fingertight	JR-5508-5

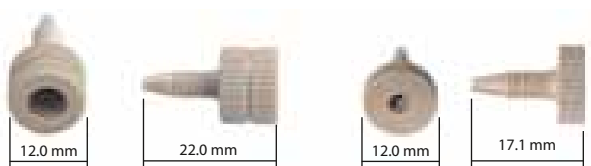

Standard

Short

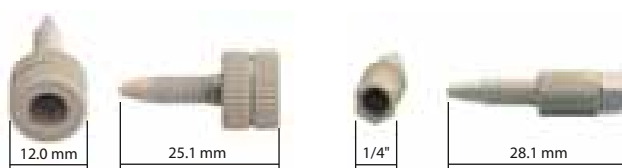

Long

Narrow hex-head/fingertight

**TECH TIP**

Coloration pigments are carefully chosen so as to contain no heavy metals or hazardous chemicals.

**MORE INFORMATION**

About PEEK ..... 99  
About CTFE ..... 99  
Chemical resistance 96-98



**PEEK no-twist "one-piece" fittings**

10-32

- Snap-in ferrule rotates freely
- Choice of ferrule materials
- Choice of fitting lengths

No-twist fittings offer the convenience of a one-piece fitting while solving a problem inherent to such designs. In other one-piece designs, the ferrule rotates against the fitting detail, creating particulates. The no twist design has a separate ferrule that snaps into the nut, so it's attached but still free to avoid rotation during tightening.

**SPECS**

- For 1/16" OD tubing
- Materials:  
Nut: PEEK  
Ferrule: PEEK, glass-filled PEEK, or CTFE
- 10-32 threads
- Pressure rating:  
<350 bar (5000 psi)  
(Varies with tubing material and ID)





**OPTIONS**

- Optional ferrule materials available: FEP, PFA, PTFE, and glass-filled PTFE. Contact us for availability.
- Some 1/32" versions are available. Call for details.


Since the ferrule is not machined onto the nut, it can be made from a different material. We offer several PEEK nuts, with ferrules of PEEK, CTFE, or glass-filled PEEK.

Sold in packages of 5.

<b>Ferrule material:</b>	<b>PEEK</b>	<b>Glass-filled PEEK</b>	<b>CTFE</b>
<i>Description</i>	<i>Product No.</i>	<i>Product No.</i>	<i>Product No.</i>
Fingertight	JR-ZNF1FPK-5	JR-ZNF1FPKG-5	JR-ZNF1FKF-5
Hex, short	JR-ZNF1PK-5	JR-ZNF1PKG-5	JR-ZNF1KF-5
Hex, medium	JR-MZNF1PK-5	JR-MZNF1PKG-5	JR-MZNF1KF-5
Hex, long	JR-LZNF1PK-5	JR-LZNF1PKG-5	JR-LZNF1KF-5
Replacement ferrules	JR-ZNF1PK-F-5	JR-ZNF1PKG-F-5	JR-ZNF1KF-F-5

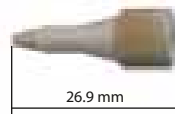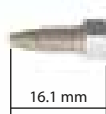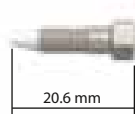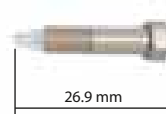**Fingertight****Hex, short****Hex, medium****Hex, long****RELATED PRODUCTS**

1/16" OD tubing  
for high pressure  
applications . . . . . 5-8  
PEEK starter kit . . . . . 90  
ValvTool for tightening  
hex-head fittings . . . . 88  
Color-It fingertight  
adapters. . . . . 17

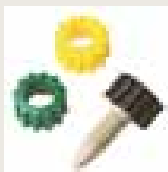

PEEK ferrules **NEW** **UHPLC**

SPECS

- For 1/16" OD tubing
- Material: PEEK, special black
- Pressure rating on SS tubing and with SS nut: (preliminary) 1170 bar (17,000 psi)

- Biocompatible 1/16" ferrule for ultra-high pressures
- Alternative to titanium ferrules in applications with titanium valves
- Excellent in combination with stainless nuts

Sold in packages of 5.

Product No.  
JR-5010-5



PEEK nuts with single ferrules 10-32

SPECS



- Material: PEEK
- Pressure rating: <280 bar (4000 psi) (Varies with tubing material and ID)


- Reuseable nut, replaceable ferrule
- For 1/32", 1/16", and 1/8" tubing
- Suitable for pressures to 280 bar (4000 psi)

Fingertight nuts have a knurled surface designed to allow sufficient sealing force to be applied to the ferrule without wrenches. Hex style nuts allow wrench tightening; however, since they are polymeric, they can break, and are recommended for use only when space is limited and fingers will not fit.


Sold in packages of 5.



	Tubing OD	Description	Thread	Product No.
Nuts	1/32"	Fingertight	6-40	JR-ZN.5FPK-5
	1/16"	Fingertight	10-32	JR-5570-5
	1/8"	Hex-head	5/16-24	JR-ZN2PK-5
Single ferrules	1/32"			JR-ZF.5PK-5
	1/16"			JR-5003-5
	1/8"			JR-ZF2PK-5

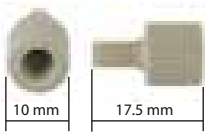

1/16" fingertight

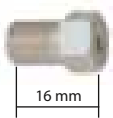

1/8" hex-head

MORE INFORMATION

About PEEK ..... 99  
Chemical resistance 96-98

RELATED PRODUCTS

1/16" OD tubing for high pressure applications ..... 5-8  
PEEK starter kit ..... 90  
PEEK tubing ..... 5-8  
Stainless steel nuts .... 31  
ValvTool for tightening hex-head fittings .... 88  
Color-It fingertight adapters..... 17



**PEEK nuts with double ferrules**

10-32

- High pressure rating - up to 420 bar (6000 psi)
- Ideal for 1/16" PEEK tubing
- Fit any 10-32 fitting

These nuts are for use with VICI Jour double ferrules, which grip the tubing twice for added stability. The fingertight version may be hand-tightened for operation at up to 420 bar (6000 psi). The hex-head version can be used for connections that are difficult to reach, or are closely spaced. The fittings work in virtually any female 1/16" fitting.





Sold in packages of 5.

**SPECS**

- For 1/16" OD tubing
- Material: PEEK
- 10-32 threads
- Pressure rating: <420 bar (6000 psi)  
(Varies with tubing material and ID)

Description	Product No.
1/16" nuts    Hex-head, short	JR-5511-5
Hex-head, long	JR-5510-5
Fingertight	JR-5580-5
Universal fingertight*	JR-5503-5*
1/16" ferrule    Double	JR-5004-5

\*Limited to existing stock.

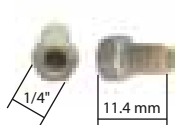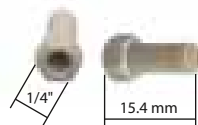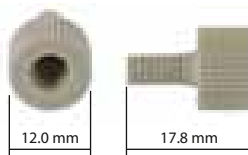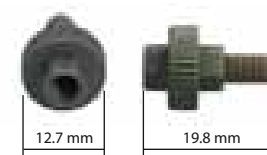**Hex-head, short****Hex-head, long****Fingertight****Universal fingertight**

### PEEK plugs and caps


#### SPECS

- Material: PEEK
- Pressure rating: <350 bar (5000 psi)  
(Varies with tubing material and ID)


- For plugging ports and capping solvent lines
- Available in all common thread sizes

Sold individually.


	Threads	Product No.
Plugs	6-40	JR-ZP5FPK
	10-32	JR-408
	fingertight	JR-ZP1FPK
	hex	JR-MZP1PK
	long hex	JR-LZP1PK
	1/4-28	JR-409
Caps	M6	JR-410
	6-40	JR-ZC.5FPK
	10-32	JR-ZC1PK
	1/4-28	JR-CCPK-1

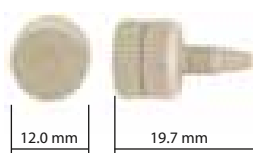


**10-32 plug**  
(JR-408)

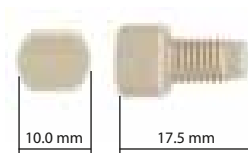


**1/4-28 plug**

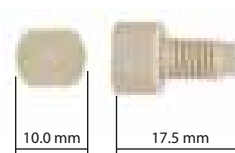

**M6 plug**

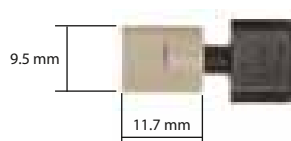

**10-32 cap**

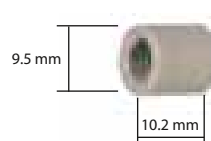

**1/4-28 cap**



### Nut, ferrule, and bushing

1/4-28

- High pressure rating – up to 480 bar (7000 psi) with 1/16" tubing
- For 1/4-28 connections
- Available for 1/32" and 1/16" tubing


This combination of PPS nut, PEEK ferrule, and stainless steel bushing allows a high pressure connection in a 1/4-28 flat-bottomed fitting detail.

Sold in packages of 10.

Order nuts and bushing/ferrules separately.

Description	Tubing OD	Product No.
PPS hex-head nut	1/32" or 1/16"	JR-55050-10



Bushing and ferrule	1/32"		JR-FPKBSS.5-10
---------------------	-------	--	----------------



	1/16"		JR-FPKBSS1-10
--	-------	--	---------------

### SPECS

- For 1/32" and 1/16" OD tubing
- Materials  
Nut: PPS  
Ferrule: PEEK  
Bushing: SS
- 1/4-28 connection
- Pressure rating  


OD	bar	psi
1/32"	175	2500
1/16"	480	7000

 (Varies with tubing material and ID)

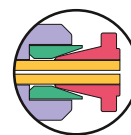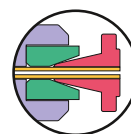

### LOW PRESSURE 1/4-28 NUTS AND FERRULES

1/4-28 PPS hex-head nuts are available with three different sealing designs for flangeless and flanged fittings ..... 41

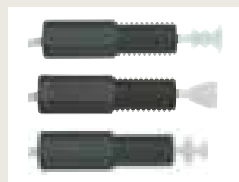

### TECH TIP

To make up high pressure connection with nut, ferrule, and bushing, first slide the nut on the tubing, then the bushing, and then the PEEK ferrule as shown in the product photo.

## PEEK unions

10-32

**SPECS**

- For 1/16" OD tubing
- Material: PEEK
- 10-32 threads
- Pressure rating: <350 bar (5000 psi)  
(Varies with tubing material and ID)

- Easy connection of transfer/solvent lines
- Low dead volume
- Completely inert and biocompatible


VICI Jour PEEK unions are usable with most common organic and inorganic HPLC solvents without risking corrosion, metal contamination, or protein degradation that might occur with stainless steel.

Small diameter bores are used to ensure low dead volume connections. These unions have the standardized Valco fitting details with smooth transition of the taper into the threaded part. Connections to these PEEK unions can be made with any 10-32 fitting that uses polymeric ferrules.

PEEK unions are available in body-only versions as well as complete with our PEEK one-piece fingertight fittings ready to connect to any 1/16" OD tubing.

	Bore (mm)	Qty/pkg	Product No.
Union	0.25	1	JR-1061
with fittings	0.50	1	JR-1066
	0.75	1	JR-1067
Body only	0.25	1	JR-061
	0.50	1	JR-066
	0.75	1	JR-067

Replacement fittings	Qty/pkg	Product No.
PEEK 1/16" one-piece fingertight fittings, 10-32, black	5	JR-55021-5

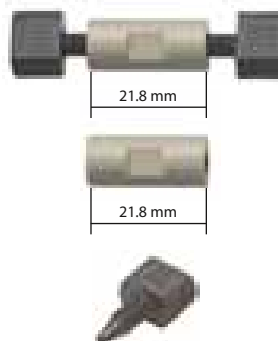

## PEEK bulkhead unions

10-32

**SPECS**

- For 1/16" OD tubing
- Material: PEEK
- 10-32 threads
- Pressure rating: <350 bar (5000 psi)  
(Varies with tubing material and ID)

- Designed for panel mounting
- Easy connection of transfer/solvent lines
- Low dead volume
- Completely inert and biocompatible


Bulkhead unions are designed to be mounted through an instrument panel or on a bracket. Typically the mounting nut goes inside the instrument, so that the long threaded portion will be out of sight.

Bulkhead unions require a minimum 10 mm panel hole, 9 mm maximum thickness.

Includes fittings.

Bore (mm)	Qty/pkg	Product No.
0.25	1	JR-ZBU1CFPK
0.50	1	JR-ZBU1MFPK
0.75	1	JR-ZBU1FPK

Replacement fittings	Qty/pkg	Product No.
PEEK 1/16" one-piece fingertight fittings, 10-32, black	5	JR-55021-5
JAM-Nut replacement mounting nut for bulkhead fittings	1	JR-CBN

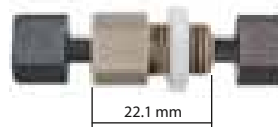**RELATED PRODUCTS**

Low pressure PEEK  
1/4-28 unions. . . . . 44  
PEEK one-piece  
fingertight fittings 19-21

**PEEK reducing unions**

- Connects 1/16" to 1/32" tubing, or 1/8" to 1/16" tubing
- Easy connection of transfer/solvent lines
- Completely inert and biocompatible

VICI Jour PEEK reducing unions are usable with most common organic and inorganic HPLC solvents without risking corrosion, metal contamination, or protein degradation that might occur with stainless steel. These reducing unions have the standardized Valco fitting details with smooth transition of the taper into the threaded part.

Includes fittings.

<i>Tubing OD</i>	<i>Bore (mm)</i>	<i>Qty/pkg</i>	<i>Product No.</i>
1/16" to 1/32"	0.25	1	JR-ZRU1.5FPK
1/8" to 1/16"	0.75	1	JR-ZRU21FPK

<b>PEEK replacement fittings</b>	<i>Qty/pkg</i>	<i>Product No.</i>
----------------------------------	----------------	--------------------

For 1/16" to 1/32" reducing unions


1/16"	One-piece fitting, black	5	JR-55021-5
1/32"	Fingertight nut	1	JR-ZN.5FPK
1/32"	Ferrule	1	JR-ZF.5PK

For 1/8" to 1/16" reducing unions

1/8"	Hex-head PEEK nut	5	JR-ZN2PK-5
1/8"	PEEK ferrule	5	JR-ZF2PK-5
1/16"	One-piece fitting, black	5	JR-55021-5

**SPECS**

- Material: PEEK
- Threads: 1/32": 6-40  
1/16": 10-32  
1/8": 5/16-24
- Pressure rating: <350 bar (5000 psi)  
(Varies with tubing material and ID)

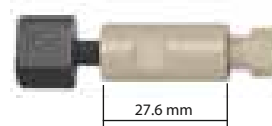**RELATED PRODUCTS**

Low pressure PEEK  
1/4-28 reducing  
unions ..... 45



## PEEK tees and crosses

10-32

### SPECS

- For 1/16" OD tubing
- Material: PEEK
- 10-32 threads
- Body size: 25.0 mm OD 12.0 mm thick
- Pressure rating: <350 bar (5000 psi) (Varies with tubing material and ID)

- Easy connection of transfer/solvent lines
- Low dead volume
- Completely inert and biocompatible

The 100% PEEK flowpath makes VICI Jour PEEK tees and crosses usable with most common organic and inorganic HPLC solvents without risking corrosion, metal contamination, or protein degradation that might occur with stainless steel.

Small diameter through-holes are used to ensure low dead volume connections. These unions have the standardized Valco fitting details with smooth transition of the taper into the threaded part. Connections to these PEEK tees and crosses can be made with any 10-32 fitting that uses polymeric ferrules.

PEEK tees and crosses are available in body-only versions as well as complete with our PEEK one-piece fingertight fittings ready to connect to any 1/16" OD tubing.

	Bore (mm)	Qty/pkg	With fittings Product No.	Body only Product No.
Tee	0.25	1	JR-1030	JR-2030
	0.50	1	JR-1032	JR-2032
	0.75	1	JR-1033	JR-2033
Cross	0.25	1	JR-1040	JR-2040
	0.50	1	JR-1042	JR-2042
	0.75	1	JR-1043	JR-2043
<b>Replacement fittings</b>		<b>Qty/pkg</b>	<b>Product No.</b>	
PEEK 1/16" one-piece fingertight fittings, 10-32, black		5	JR-55021-5	



## Biocompatible mixing tee

10-32

### SPECS

- For 1/16" OD tubing
- Material: PEEK
- 10-32 threads
- Pressure rating: <350 bar (5000 psi) (Varies with tubing material and ID)

- High pressure gradient mixing
- Post column reactions
- Completely biocompatible
- Built-in protein-friendly titanium filter, 10 µm


A unique, more thorough mixing results from the turbulence generated by the introduction angle of the solvents as well as by the flow through our titanium filter element. The biocompatible tee can be used for high pressure gradient mixing as well as for post column reactions.



Mixing tee with fittings

Product No.

JR-9000-0670

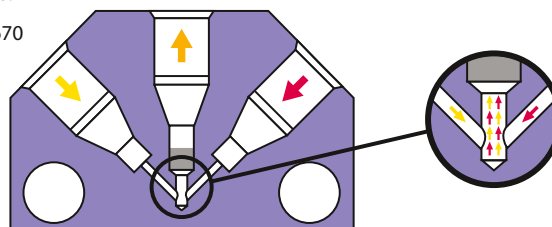

The two solvent streams are well-mixed in the vortex mixing chamber

### TECH TIP

Use our color-coded one-piece fingertight fittings for easy solvent line identification. . . . . 20



### MORE INFORMATION

About PEEK . . . . . 99  
Chemical resistance 96-98

### RELATED PRODUCTS

Low pressure 1/4-28 tees and crosses . . . . 46  
Low pressure 1/4-28 mixing tees . . . . . 47

**PEEK adapters**

- Connect various combinations of thread types, male and female
- Includes male and female Luer adapters
- Low dead volume
- For any 10-32 high pressure fitting

We offer male and female Luer adapters as well as adapters to connect standard LC fittings to standard LC ports.

The 100% PEEK flowpath makes VICI Jour PEEK adapters usable with most common organic and inorganic HPLC solvents without risking corrosion, metal contamination, or protein degradation that might occur with stainless steel.


Sold individually.

**SPECS**


- Material: PEEK
- Bore: 0.40 mm
- Pressure ratings:  
M6, 1/4-28 <70 bar (1000 psi)  
Luer <35 bar (502 psi)  
(Varies with tubing material and ID)

Description		Product No.
1/4-28 female to	10-32 male	JR-0612
	10-32 female	JR-062
	M6 male	JR-0622
	Luer female	JR-0625
	Luer male	JR-0626
10-32 female to	1/4-28 male	JR-0617
	M6 female (Pharmacia*)	JR-5595
	Luer female	JR-0623
Waters female to	10-32 male	JR-0611
M6 female to	10-32 male	JR-5505
Female Luer to	1/4-28	JR-CFLAPK
	10-32	JR-ZUFLPK
Male Luer to	1/4-28	JR-CMLAPK


**1/4-28 female to:**

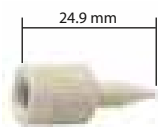


**10-32 male**

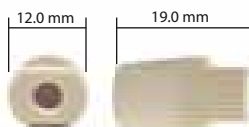


**10-32 female**

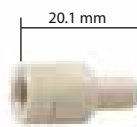

**M6 male**

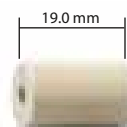


**Luer female**

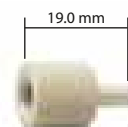


**Luer male**


**10-32 female to:**

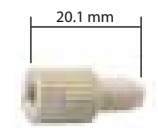

**1/4-28 male**

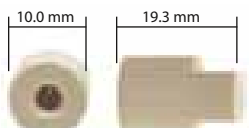


**M6 female  
(Pharmacia)**

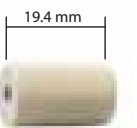


**Luer female**

**Waters female to:**

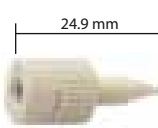

**10-32 male**

**M6 female to:**

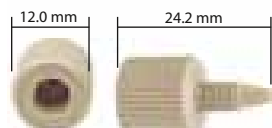

**10-32 male**

\* Pharmacia is a division of GE Healthcare.

# High pressure fittings

STAINLESS STEEL

Our stainless steel fittings are typically compression fittings, in which a ferrule is compressed onto the tube as a nut is tightened. They offer the best reliability in high pressure situations.

## High pressure stainless fittings index

### Nuts

Valco nuts ..... page 31

Valco nuts are available for 1/32",  
1/16", and 1/8" tubing, in several  
lengths.

Rheodyne® compatible nuts .....31

These nuts are sized for 1/16"  
tubing, and come in three lengths.

### Ferrules

Valco .....32

Rheodyne compatible.....32

Plugs and caps .....32

### Tees and crosses

Tees .....33

Crosses .....33

### Unions and reducing unions

Unions.....34

Reducing unions

Imperial .....34

Imperial to metric **NEW** .....35

### DISCLAIMER

The maximum holding pressure for any type of connections involving tubing and a ferrule varies considerably with the tubing material, the ferrule material, the clearance between tubing OD and ferrule ID, and the shape of the fitting detail.

### MORE INFORMATION

ChromBox kit..... 91

HPLC parts kit with  
PEEK and stainless  
fittings and tubing

Fingertight adapters

Color-It ..... 17

Color-Sleeve..... 17

Nanovolume® high  
pressure fittings...36-39

Stainless starter kits ... 90

Tools for fittings ....88-89



**Valco® stainless nuts****UHPLC**

- Excellent fit to Valco fitting details
- Different lengths available
- For 1/32", 1/16", and 1/8" OD tubing




These nuts are for use with all standard Valco internal fittings and most valves. They may be used with fittings from other manufacturers as well. The medium and long types are ideal for metal Valco and Cheminert valves with more than 8 ports, or where the fitting head may be otherwise inaccessible.

Sold in packages of 5.

Tubing OD	Length	Product No.
1/32"	Standard	JR-ZN.5-5
1/16"	Standard	JR-ZN1-5
	Medium	JR-MZN1-5
	Long	JR-LZN1-5
1/8"	Standard	JR-ZN2-5

**SPECS**

- For 1/32", 1/16", and 1/8" OD tubing
- Material: Type 303 SS
- Threads: 1/32": 6-40  
1/16": 10-32  
1/8": 5/16-24
- Pressure rating with premium SS tubing and ZF1S6 ferrule: 1375 bar (20,000 psi)

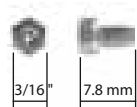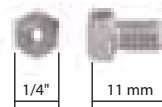
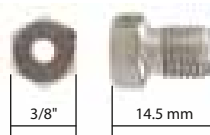
**1/32" standard****1/16" standard****1/16" medium****1/16" long****1/8" standard****Rheodyne® compatible stainless nuts**

- Excellent fit to Rheodyne fitting details
- Different lengths available
- For 1/16" OD tubing

These nuts are for use with metal Rheodyne valves. Three different length nuts are available, in order to accommodate any possible connecting situation. The standard nut is suitable for easily accessible connections (couplers, detectors, columns, etc.) For very tight and difficult to reach areas, we recommend the long and extra long nut.

Sold in packages of 5.

Tubing OD	Length	Product No.
1/16"	Short	JR-56-5
	Long head	JR-57-5
	Extra long	JR-58-5

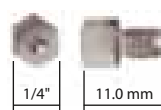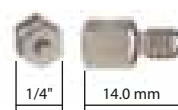**Short****Long head****Extra long****SPECS**

- For 1/16" OD tubing
- Material: Type 316 SS
- Threads: 10-32
- Pressure rating: <500 bar (7,250 psi) (Varies with tubing material and ID)

**TECH TIP**

Using nuts with different lengths will create more space, for better wrench access.

**MORE INFORMATION**

About stainless steel... 100

**RELATED PRODUCTS**

Ferrules

Rheodyne ..... 32

Valco ..... 32

Stainless steel

starter kit..... 90


Stainless tubing .... 14-15

ValvTool for tightening

hex-head fittings .... 88

Color-It fingertight

adapters..... 17

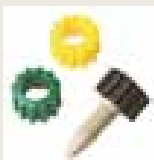

## Valco® stainless ferrules

UHPLC




### SPECS

- For 1/32", 1/16", and 1/8" OD tubing
- Material: Type 316 SS
- Pressure rating: 1375 bar (20,000 psi) (Varies with tubing material and ID)

- Excellent fit to Valco fitting details
- For 1/32", 1/16", and 1/8" OD tubing

These ferrules are for use with all standard Valco internal fittings and most valves, and with those of most other manufacturers as well. Valco metal ferrules cut a ring near the end of the tube, preventing tube release at high pressures without significantly deforming and restricting the tube interior. (However, if the hardness of the tubing is equal to or greater than that of the ferrule, deformation of the tube rather than a cut ring is likely.)

Sold in packages of 5.

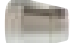Size		Product No.
1/32"		JR-ZF5-5
1/16"		JR-ZF1-5
1/8"		JR-ZF2-5

## Rheodyne® compatible stainless ferrules



### SPECS

- For 1/16" OD tubing
- Material: Type 316 SS
- Pressure rating: <500 bar (7,250 psi) (Varies with tubing material and ID)

- Excellent fit to Rheodyne fitting details
- For 1/16" OD tubing

These ferrules are for use with metal Rheodyne valves.

Sold in packages of 5.

Description		Product No.
Standard ferrule		JR-59-5
2-step Rheodyne-type ferrule		JR-60-5

## Stainless steel plugs

UHPLC


### SPECS

- For 1/16", and 1/8" OD tubing
- Material: 1/16": SS 316, 1/8": SS 313
- Threads: 1/16": 10-32, 1/8": 5/16-24

Use a plug to close off an unused port in a valve, manifold, or HPLC column. The standard stainless plug consists of a zero volume nut with a ferrule made up on a piece of solid rod.

Sold individually.

Tubing OD	Product No.
1/32"	JR-ZP5
1/16"	JR-ZP1
1/8"	JR-ZP2

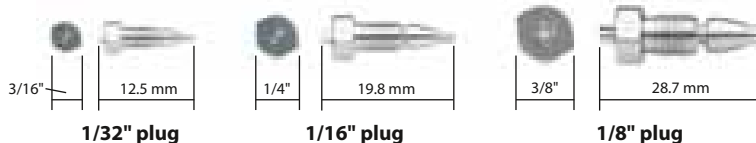

## Stainless steel caps

UHPLC


### SPECS

- For 1/16", and 1/8" OD tubing
- Material: Type 303 SS
- Threads: 1/16": 10-32, 1/8": 5/16-24

Use a cap to close off a tube with a ferrule and internal nut already made up. A cap is essentially a piece of hex stock with a zero volume fitting detail machined into it, but with no through-hole.

Sold individually.

Tubing OD	Product No.
1/32"	JR-ZC5
1/16"	JR-ZC1
1/8"	JR-ZC2

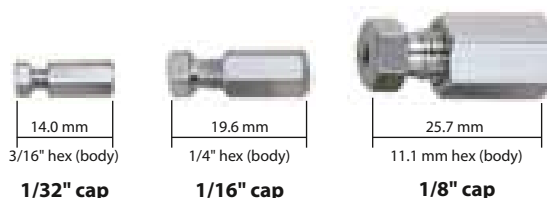

**RELATED PRODUCTS**  
UHPLC PEEK ferrules . . . 22

## HIGH PRESSURE FITTINGS ■ STAINLESS STEEL

- Easy connection of solvent lines
- Low dead volume
- For 1/16" OD tubing

VICI stainless steel tees and crosses are for joining solvent lines. They are precision machined with a focus on the concentricity of the bore as well as on the finish of the surfaces to give an absolute leak-free connection.

These tees and crosses have internal threads, which we recommend for most applications. Internal threads make stronger connections and internal fitting details provide the lowest dead volume for high performance instrumentation.



## Stainless tees

UHPLC

VICI stainless steel tees join three solvent lines.

Three nuts and ferrules ship with each tee.

Description	Bore (mm)	Qty/pkg	Product No.
1/16" tee	0.25	1	JR-ZT1C
	0.50	1	JR-ZT1M
	0.75	1	JR-ZT1
Replacement nut		5	JR-ZN1-5
Replacement ferrule		5	JR-ZF1S6-5

## SPECS

- For 1/16" OD tubing
- Material:  
Body Type 316 SS  
Nuts Type 303 SS  
Ferrules 316 SS
- Pressure rating:  
1375 bar  
(20,000 psi)  
(Varies with  
tubing material  
and ID)



## Stainless crosses

UHPLC

VICI stainless steel crosses join four solvent lines.

Four nuts and ferrules ship with each tee.

Description	Bore (mm)	Qty/pkg	Product No.
1/16" cross	0.25	1	JR-ZX1C
	0.50	1	JR-ZX1M
	0.75	1	JR-ZX1
Replacement nut		5	JR-ZN1-5
Replacement ferrule		5	JR-ZF1S6-5

## SPECS

- For 1/16" OD tubing
- Material:  
Body Type 316 SS  
Nuts Type 303 SS  
Ferrules 316 SS
- Pressure rating:  
1375 bar  
(20,000 psi)  
(Varies with  
tubing material  
and ID)

## MORE INFORMATION

About stainless steel... 100

## RELATED PRODUCTS

Open end wrench ..... 89

Stainless steel  
starter kit..... 90

ValvTool for tightening  
hex-head fittings .... 88

- Easy connection of solvent lines
- Low dead volume

VICI stainless steel unions and reducing unions are precision machined with a focus on the concentricity of the bore as well as on the finish of the surfaces to give an absolute leak-free connection.

These unions have internal threads, which we recommend for most applications. Internal threads make stronger connections and internal fitting details provide the lowest dead volume for high performance instrumentation.

## Stainless unions

UHPLC

### SPECS

- For 1/16" OD tubing
- Material:  
Body Type 316 SS  
Nuts Type 303 SS  
Ferrules 316 SS
- Pressure rating:  
1375 bar  
(20,000 psi)  
(Varies with  
tubing material  
and ID)

- Choice of four bores
- For 1/16" OD tubing

VICI stainless steel unions join two pieces of 1/16" OD metal tubing.

Two nuts and ferrules ship with each union.

Description	Bore (mm)	Qty/pkg	Product No.
1/16" union	0.15	1	JR-ZU1XC
	0.25	1	JR-ZU1CS6
	0.50	1	JR-ZU1MS6
	0.75	1	JR-ZU1S6
Replacement nut		5	JR-ZN1-5
Replacement ferrule		5	JR-ZF1S6-5



## Stainless reducing unions

Imperial

UHPLC


### SPECS

- For joining 1/8" to 1/16" OD tubing
- Material:  
Body Type 316 SS  
Nuts Type 303 SS  
Ferrules 316 SS
- Pressure rating:  
1375 bar  
(20,000 psi)  
(Varies with  
tubing material  
and ID)

These reducing unions join 1/16" OD to 1/8" OD tubing.

Two nuts and ferrules ship with each union.

Description	Size	Qty/pkg	Product No.
Reducing union	1/8" to 1/16"	1	JR-ZRU21
Replacement nut	1/8"	5	JR-ZN2-5
	1/16"	5	JR-ZN1-5
Replacement ferrule	1/8"	5	JR-ZF2S6-5
	1/16"	5	JR-ZF1S6-5

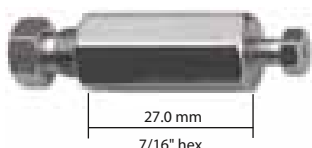

### RELATED PRODUCTS

Nanovolume® high pressure fittings . . . 36-39  
 Stainless starter kit . . . 90  
 HPLC parts kit with stainless fittings and tubing  
 ChromBox kit. . . . . 91  
 HPLC parts kit with PEEK and stainless fittings and tubing





### Stainless reducing unions **NEW** Imperial to metric **UHPLC**

- Ultra-high pressure
- For 1 mm, 2 mm, and 3 mm OD tubing

With these new reducing unions, metric tubing can be connected to imperial fittings, valves, and chromatographic systems. The unions are manufactured from Type 303 stainless steel.

For easy identification, the metric end is marked by a groove.




Two nuts and ferrules ship with each reducing union.

#### SPECS

- Material:
  - Body Type 303 SS
  - Nuts Type 303 SS
  - Ferrules 316 SS
- Tube reducing OD
  - 1/32" to 1.0 mm
  - 1/16" to 2.0 mm
  - 1/8" to 3.0 mm
- Pressure rating:
  - 1375 bar
  - (20,000 psi)
  - (Varies with tubing material and ID)

Description	Bore (mm)	Threads	Qty/pkg	Product No.
Reducing union				
1/32" to 1.0 mm	0.25	6-40 to 10-32	1	JR-096
1/16" to 2.0 mm	0.75	10-32	1	JR-082
1/8" to 3.0 mm	0.75	5/16-24	1	JR-086

Description	Size	Threads	Qty/pkg	Product No.
Replacement nut	1/32"	6-40	5	JR-ZN1-5
	1/16"	10-32	5	JR-ZN1-5
	1/8"	5/16-24	5	JR-ZN2-5
	for 1.0 mm	10-32	5	JR-ZN1-5
	for 2.0 mm	10-32	5	JR-80-5
	for 3.0 mm	5/16-24	5	JR-ZN2-5
Replacement ferrule	1/32"		5	JR-ZF5S6-5
	1/16"		5	JR-ZF1S6-5
	1/8"		5	JR-ZF2S6-5
	1.0 mm	<b>NEW</b>	5	JR-097-5
	2.0 mm	<b>NEW</b>	5	JR-083-5
	3.0 mm	<b>NEW</b>	5	JR-087-5

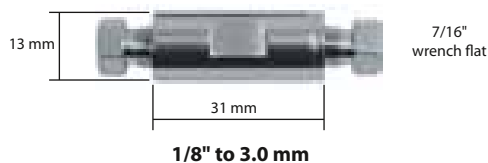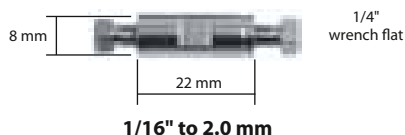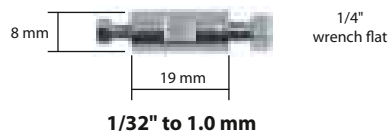

**Metric ferrules**

#### WRENCH SIZE FOR BODY FLAT

Wrench size	Reducing union	P/N
1/4"	1/32" to 1.0 mm	JR-096
1/4"	1/16" to 2.0 mm	JR-082
7/16"	1/8" to 3.0 mm	JR-086

# High pressure fittings

NANOVOLUME® FITTINGS ■ FOR 360 MICRON AND 1/32" TUBING

VICI Nanovolume® fittings generally have bore sizes of 100-150 µm (.004" - .006"), with some as small as 50 µm (.002"). The minimal transfer volume contributed by Nanovolume® components makes them especially beneficial in applications with flow rates in the µl/min range, when the transfer volume can be critical.

## High pressure Nanovolume® fittings index

### 1/32" Nanovolume® Fittings

Nuts, ferrules, unions, reducing unions . . . . . 37

1/32" fittings, with 100 µm or 150 µm bore, are ideal for high resolution capillary chromatography. Rated at 5,000 psi with fingertight nuts, they will remain leak-tight well beyond the burst strength of most PEEK tubing. These fittings are machined from high quality inert polymers to the same exacting tolerances as our popular Valco zero dead volume fittings, and the taper angle and detail design conform to the industry standard established by the Valco line.

### 360 Micron Nanovolume® Fittings

For PEEK or FS tubing, <10K psi . . . . . 38

For fused silica tubing, 10K psi . . . . . 39

For metal tubing, <20K psi . . . . . 39

Our newest high pressure fittings permit direct connection of 360 micron OD fused silica, PEEK, stainless, or electroformed nickel tubing without the use of liners. The ferrule snaps into the nut so that the fitting is "one-piece," but the ferrule remains free to rotate as the nut is tightened so that the tube doesn't twist. Because of the compact size and fine 2-56 threads, a leak-free connection that seals at pressures in excess of 20,000 psi can be easily formed with the available manual tool.

#### DISCLAIMER

The maximum holding pressure for any type of connections involving tubing and a ferrule varies considerably with the tubing material, the ferrule material, the clearance between tubing OD and ferrule ID, and the shape of the fitting detail.





- 100 and 150 µm bores available
- Nuts, ferrules, and unions
- For PEEK tubing

VICI Nanovolume® fittings are designed for capillary LC, used in HTS techniques, proteomics, genomics, and other newer developments in industrial research requiring small sample size for reduction of analysis time.

VICI Nanovolume connectors include our one-piece fingertight fittings for 1/32" OD tubing, featuring our

patented\* collapsible ferrule. PEEK ferrules for use with our external nuts also feature the collapsible, grooved design. To avoid potential confusion, all fittings utilizing the Cheminert collapsible ferrule are made of black PEEK; fittings with a standard Valco ZDV fitting detail are natural PEEK.

VICI Nanovolume fittings are as easy to handle as our standard analytical fittings and unions. The only difference is their size.

\*U.S. patent no. 6,575,501.

## Nanovolume® fittings


1/32"

The internal nut includes a collapsible ferrule as an integral part of the fitting; the external nut must be used with the separate ferrule listed below.

Description	Tubing OD	Product No.
One-piece nut/ferrule	1/32"	JR-C-NNFFPK
External nut	1/32"	JR-C-EN.5FPKB
	1/16"	JR-C-EN1FPKB
Ferrule, collapsible	1/32"	JR-ZGF.5PK
	1/16"	JR-ZGF1PK
<b>Bore size:</b>		
	Tubing OD	100 micron Product No. 150 micron Product No.
Union (external to internal)	1/32"	JR-C-NEU.5XFPK JR-C-NEU.5FPK
Reducing union	1/16" to 1/32"	— JR-C-NERU1FPK

### SPECS

- For 1/32" OD PEEK
- Material: PEEK
- Threads:  
Nut/ferrule 6-40  
External nut 1/32" 8-36  
1/16" 10-32
- Pressure rating:  
<350 bar (5,000 psi)  
(Varies with tubing material and ID)
- Tolerances:  
±0.05mm (.002")

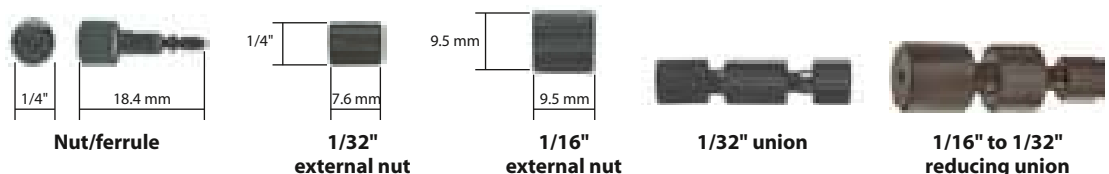

## Liners for 1/32" connectors

for use with fused silica tubing

Use these natural PEEK liners to adapt 1/32" connectors to the most common sizes of fused silica tubing. The 27 mm liners are for internal nuts with collapsible ferrules. 10 mm liners are for use with external nuts. Sold in packages of 5.

### MORE INFORMATION

360 µm Nanovolume® fittings .....38-39

### CONVERSIONS

50 µm	= .002"
100 µm	= .004"
150 µm	= .006"
0.25 mm	= .010"
0.50 mm	= .020"
0.75 mm	= .030"

### 27 mm liners

Use with internal nuts JR-C-NNFFPK or JR-C-NNFLFPK

For tubing OD	Prod No
125 - 175 µm	JR-C-NL.15L-5
175 - 225 µm	JR-C-NL.20L-5
225 - 275 µm	JR-C-NL.25L-5
275 - 325 µm	JR-C-NL.30L-5
325 - 375 µm	JR-C-NL.35L-5

### 10 mm liners

Use with external nut JR-C-EN.5FPKB

For tubing OD	Prod No
125 - 175 µm	JR-C-NL.15S-5
175 - 225 µm	JR-C-NL.20S-5
225 - 275 µm	JR-C-NL.25S-5
275 - 325 µm	JR-C-NL.30S-5
325 - 375 µm	JR-C-NL.35S-5

## Nanovolume® fittings ■ for 360 µm tubing

### HIGH PRESSURE FITTINGS ■ NANO VOLUME®

- For direct connection of 360 µm tubing
- Work with metal, fused silica, or PEEK
- Up to 20,000 psi with metal tubing
- Snap-in rotating ferrule for “one-piece” fitting with no tubing twist
- Eliminate use of troublesome liners

Our high pressure fittings permit direct connection of 360 micron OD fused silica, PEEK, stainless, or electroformed nickel tubing without the use of liners. The ferrule snaps into the nut so that the fitting is “one-piece”, but the ferrule remains free to rotate as the nut is tightened so that the tube doesn’t twist. Because of the compact size and fine 2-56 threads, a leak-free connection that seals at pressures in excess of 20,000 psi can be easily formed with the available manual tool.

VICI Nanovolume® 360 µm fittings are dedicated for use with either fused silica, metal, or PEEK tubing; components cannot be mixed or used with a different tubing material.

Tees and crosses offer a choice of three bore sizes, and feature a “quick-mount” base with adhesive backing to make sure that the fitting is stable and fragile tubing doesn’t get broken. There is also a quick-mount PEEK union.



### For PEEK or fused silica tubing — up to 10,000 psi

360 µm

#### SPECS

- For PEEK or FS tubing, 360 µm
- Material: PEEK
- Threads: 2-56
- Pressure rating: <689.5 bar (10,000 psi) (Varies with tubing material and ID)

These fittings are constructed from premium grade natural PEEK material. They are intended for use with PEEK or fused silica tubing at pressures up to 10,000 psi, or the maximum pressure for which the tubing is rated, whichever is lower. Quick-mount versions have integral base with double stick tape to secure fittings to a surface.

	Product No.
Nut/ferrule	JR-C360NFPK
Cap	JR-C360CPK
Plug	JR-C360PPK
Tightening tool	JR-C360ET

Bore size:	50 micron Product No.	100 micron Product No.	150 micron Product No.
Union	JR-C360UPK2	JR-C360UPK4	JR-C360UPK6
Reducing union, 1/16" to 360 µm	—	—	JR-C360RU1PK6
Union, quick mount	JR-C360QUPK2	JR-C360QUPK4	JR-C360QUPK6
Tee, quick mount	JR-C360QTPK2	JR-C360QTPK4	JR-C360QTPK6
Cross, quick mount	JR-C360QXPK2	JR-C360QXPK4	JR-C360QXPK6



360 µm union



1/16" to 360 µm  
reducing union



Union,  
quick mount



Tee,  
quick mount



Cross,  
quick mount

#### TECH TIP

360 micron fittings are dedicated for use with either fused silica, metal, or PEEK tubing. Components cannot be mixed or used with a different tubing material.

#### TECH TIP

Grooves on the fitting body indicate the bore size:





Grooves indicate bore size  
1 groove: 50 µm  
2 grooves: 100 µm  
no grooves: 150 µm

## For fused silica tubing — 10,000 psi and above


360 µm

These fittings are constructed from HPLC grade stainless steel, with stainless steel nut and a special ferrule which is precision machined from electroformed nickel. For optimal sealing characteristics, the ferrule is gold plated.

	Product No.		
Nut/ferrule 	JR-C360NFFS		
Cap 	JR-C360CFS		
<b>Bore size:</b>	<b>50 micron</b>	<b>100 micron</b>	<b>150 micron</b>
	Product No.	Product No.	Product No.
Union	JR-C360UFS2	JR-C360UFS4	JR-C360UFS6
Reducing union, 1/32" to 360 µm	JR-C360RU.5FS2	JR-C360RU.5FS4	JR-C360RU.5FS6
Reducing union, 1/16" to 360 µm	—	—	JR-C360RU1FS6

### SPECS

- For FS tubing, 360 µm
- Materials  
Body and nut:  
HPLC grade SS  
Ferrule:  
Gold-plated  
nickel
- Threads: 2-56
- Pressure rating:  
≥ 689.5 bar  
(10,000 psi)  
(Varies with  
tubing material  
and ID)

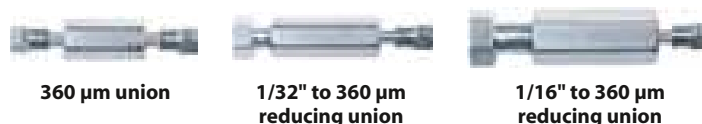

## For metal tubing — up to 20,000 psi


360 µm

Our highest pressure Nanovolume® fittings are constructed of HPLC grade stainless steel, including stainless steel nut and ferrule. These fittings are optimized for use with stainless or electroformed nickel tubing.

	Product No.		
Nut/ferrule 	JR-C360NFS6		
Cap 	JR-C360C		
<b>Bore size:</b>	<b>50 micron</b>	<b>100 micron</b>	<b>150 micron</b>
	Product No.	Product No.	Product No.
Union	JR-C360US62	JR-C360US64	JR-C360US66
Reducing union, 1/32" to 360 µm	JR-C360RU.5S62	JR-C360RU.5S64	JR-C360RU.5S66
Reducing union, 1/16" to 360 µm	—	—	JR-C360RU1S66

### SPECS

- For metal tubing, 360 µm
- Materials  
HPLC grade SS
- Threads: 2-56
- Pressure rating:  
<1375 bar  
(20,000 psi)  
(Varies with  
tubing material  
and ID)

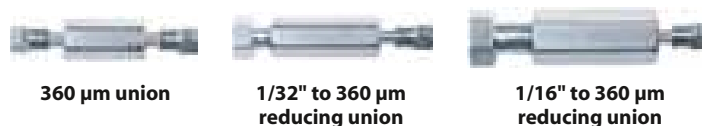

### MORE INFORMATION

1/32" Nanovolume®  
fittings ..... 37

### CONVERSIONS

50 µm	= .002"
100 µm	= .004"
150 µm	= .006"
0.25 mm	= .010"
0.50 mm	= .020"
0.75 mm	= .030"

# Low pressure fittings

POLYMERIC

Our low pressure fittings are ideally suited for flow injection analysis, low pressure liquid chromatography, and stream sampling devices. They may be safely used at pressures up to 500 psi and temperatures to 50°C.

## Low pressure polymeric fittings index

### Nuts, ferrules, and washers

#### for 1/4-28 flat-bottom fitting details

PPS hex head nuts	
Flangeless .....	page 41
Flanged .....	41
Flangeless ferrules	
Collapsible .....	41
Inverted .....	41
Flangeless nuts	
PPS hex-head .....	41
PEEK fingertight nuts .....	42
Color-coded polyacetal fingertight nuts ...	42
Flanged color-coded polypropylene	
fingertight nuts .....	43

### Unions and reducing unions

Union body, polyamide .....	43
PEEK unions .....	44
PEEK bulkhead unions .....	44
PEEK reducing unions .....	45
PEEK reducing unions, imperial to metric <b>NEW</b>	45

### Tees, crosses, and manifolds

Tees .....	46
Biocompatible mixing tee .....	47
Crosses .....	46
Manifolds with 5 and 9 ports .....	46

### Plugs

10-32 and 1/4-28 polypropylene plugs .....	48
--	----

### Adapters, polypropylene .....

48  
A wide offering of adapters, including male and female Luer to barbed, 10-32, and 1/4-28 male and female, plus many other combinations.

### Micro valves for LC and GC .....

49  
Simplify your liquid or gas handling application. The unique design of the fitting detail allows a leak-free seal with no potential for rotor damage from overtightening.

### Peristaltic tubing adapter set .....

49  
Connect soft tubing to PTFE and PTFE-like tubing with this biocompatible adapter set, which works with peristaltic tubing up to 4.50 mm OD.

### FITTING DESIGNS

VICI Jour offers two types of low pressure fittings for 1/4-28 flat-bottomed fitting details:


- **Flangeless** connections utilize several ferrule designs. The ferrule grips the tubing as the fitting is tightened, without significantly reducing the tubing ID.
- **Flanged** tube end fittings are retained on polymeric tubing by a flange formed around a polypropylene washer. The washer never comes in contact with the mobile phase.

### DISCLAIMER

The maximum holding pressure for any type of connections involving tubing and a ferrule varies considerably with the tubing material the ferrule material, the clearance between tubing OD and ferrule ID, and the shape of the fitting detail.



- Economically priced alternative to PEEK
- Excellent chemical resistance
- For 1/16" or 1/8" OD tubing
- Fingertight when used with Color-Sleeve adapters

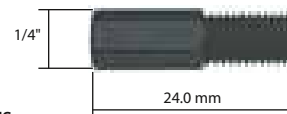


Our PPS 1/4-28 fitting series offers excellent chemical and pressure resistance. Three different sealing options are available.

### PPS hex-head nuts and collapsible ferrules (flangeless)

1/4-28

Patented\* collapsible ferrules are not included. Order separately.

\* Patent No. 6,575,501

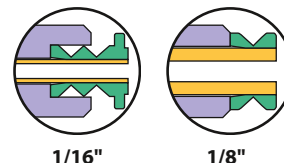

Description	Material	Tubing OD	Qty/pkg	Product No.
Nut for use with collapsible ferrule	PPS	1/16"	10	JR-55050-10
		1/8"	10	JR-55051-10
Collapsible ferrules	CTFE	1/16"	10	JR-CFL-CB1KF
		1/8"	10	JR-CFL-CB2KF
	PEEK	1/16"	10	JR-CFL-CB1PK

#### SPECS

- For 1/16" and 1/8" OD tubing
- Threads: 1/4-28
- Pressure rating:  

OD	Bar	Psi
1/16"	105	1500
1/8"	35	500

 (Varies with tubing material and ID)
- \* Patent No. 6,575,501

### HIGH PRESSURE 1/4-28 NUTS AND FERRULES

A PPS nut, PEEK ferrule, and stainless steel bushing allow a high pressure connection in a 1/4-28 flat-bottomed fitting detail... 25



### TECH TIP

Color-Sleeve adapters turn PPS nuts into fingertight fittings and color-code different solvent lines. .... 17



### CHEMICAL RESISTANCE OF PPS

A = suitable  
B = marginal

C = not recommended

Chemical class	Resistance
Acids, organic	A
Acids, inorganic	B/C*
Aldehydes	A
Alcohols	A
Bases	B
Esters	A
Halogenated org.	A/B
Hydrocarbons	A
Ketones	A

\* Concentrated halogen acids  
See detailed chemical resistance chart .... 96-98


### MORE INFORMATION

About polymers ... 99-100

### PPS hex-head nuts and inverted ferrules (flangeless)

1/4-28

Ferrules are not included. Order separately.

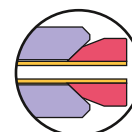

Description	Material	Tubing OD	Qty/pkg	Product No.
Nut for use with inverted ferrule	PPS	1/32"	10	JR-55073-10
		1/16"	10	JR-55070-10
		1/8"	10	JR-55071-10
		1.0 mm	10	JR-55073-10
		2.0 mm	10	JR-55072-10
		3.0 mm	10	JR-55074-10
Inverted ferrules	ETFE	1/32"	NEW 10	JR-032-10
		1/16"	10	JR-041-10
		1/8"	10	JR-051-10
		1.0 mm	NEW 10	JR-031-10
		2.0 mm	NEW 10	JR-042-10
		3.0 mm	NEW 10	JR-049-10

#### SPECS

- For 1/32"-1/8" OD tubing
- Threads: 1/4-28
- Pressure rating:  

OD	Bar	Psi
1/16"	105	1500
1/8"	35	500

 (Varies with tubing material and ID)

### PPS hex-head nuts and washers (flanged connections)

1/4-28

Washers are not included. Order separately.



Description	Material	Tubing OD	Qty/pkg	Product No.
Nut for flange/washer connection	PPS	1/16"	10	JR-55060-10
		1/8"	10	JR-55061-10
Washer	PP	1/16"	10	JR-201926-10
		1/8"	10	JR-201928-10

#### SPECS

- For 1/16" and 1/8" OD tubing
- Threads: 1/4-28
- Pressure rating:  

OD	Bar	Psi
1/16"	105	1500
1/8"	35	500

 (Varies with tubing material and ID)

Easy-Flange kit (see page 85)

JR-201540



## PEEK fingertight nuts (flangeless)

1/4-28

## SPECS

- For 1/16" and 1/8" OD tubing
- Materials  
Nut: PPS  
Ferrules: ETFE
- 1/4-28 threads
- Pressure rating:  

OD	Bar	Psi
1/16"	105	1500
1/8"	35	500

 (Varies with tubing material and ID)




- Compatible with all 1/4-28 flat-bottom ports
- Excellent chemical resistance
- For 1/16" or 1/8" OD tubing
- True fingertight – no tools required

PEEK flangeless nuts in combination with our inverted ferrule do not require a flanging tool. A hand-tightening procedure makes leak-free seals.

With this fitting system, the inert ETFE ferrule is the only part besides the tubing which comes into contact with the fluid stream – the nut is never in contact.


Ferrules are not included. Order separately.

Tubing OD:		1/16"		1/8"
Description	Material	Qty/pkg	Product No.	Product No.
Nut for use with inverted ferrule	PEEK	10	JR-20124-10	JR-20116-10
Inverted ferrules	ETFE	10	JR-041-10	JR-051-10

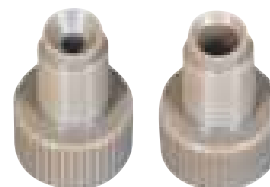

1/16"

1/8"

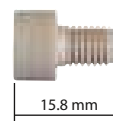

## Color-coded fingertight nuts (flangeless)

1/4-28

## SPECS

- For 1/16" and 1/8" OD tubing
- Materials  
Nut: Polyacetal  
Ferrules: ETFE
- 1/4-28 threads
- Pressure rating:  

OD	Bar	Psi
1/16"	70	1000
1/8"	35	500


 (Varies with tubing material and ID)



- Polyacetal – economical alternative to PEEK
- Compatible with all 1/4-28 flat-bottom ports
- For 1/16" or 1/8" OD tubing
- True fingertight – no tools required


Polyacetal flangeless nuts in combination with the inverted ferrule do not require a flanging tool. A hand-tightening procedure makes leak-free seals. With this fitting system, the inert ETFE ferrule is the only part besides the tubing which comes into contact with the fluid stream – the nut is never in contact. Six colors are available.

Tubing OD:			1/16"	1/8"	
Description	Color	Qty/pkg	Product No.	Product No.	
Nuts		Green	10	JR-040GR-10	JR-050GR-10
		Blue	10	JR-040BL-10	JR-050BL-10
		Yellow	10	JR-040YL-10	JR-050YL-10
		Red	10	JR-040RE-10	JR-050RE-10
		White	10	JR-040WH-10	JR-050WH-10
		Black	10	JR-040BK-10	JR-050BK-10
1" long nuts	White	10	JR-040WL-10	JR-050WL-10	
Inverted ETFE ferrules		10	JR-041-10	JR-051-10	

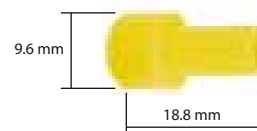


## TECH TIP

The inert ETFE ferrule used with these nuts is the only part besides the tubing which is in contact with the fluid stream.

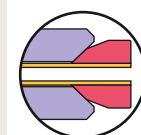

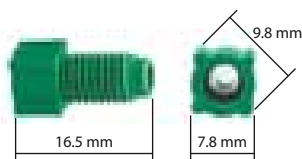

**Color-coded fingertight nuts (flanged connections)**

1/4-28

- For flanged connections
- Compatible with all 1/4-28 flat-bottom ports
- For 1/16" or 1/8" OD tubing
- Fiberglass-reinforced PP with PP washers

For those systems requiring or preferring flanged connections, use these fiberglass-reinforced fittings with our Easy-Flange flanging tool and polypropylene washers to make a clean connection. The flange on the tubing ensures that the washer never comes in contact with the mobile phase.


Six colors are available.

**SPECS**

- For 1/16" and 1/8" OD tubing
- Material: Fiberglass-reinforced polypropylene
- 1/4-28 threads
- Pressure rating: <30 bar (435 psi) (Varies with tubing material and ID)

**Tubing OD:**

Description	Color	Qty/pkg	1/16" Product No.	1/8" Product No.
Nuts	 Green	10	JR-202123-10	JR-202143-10
	 Blue	10	JR-202122-10	JR-202142-10
	 Yellow	10	JR-202124-10	JR-202144-10
	 Red	10	JR-202125-10	JR-202145-10
	 White	10	JR-201094WH-10	JR-201580WH-10
	 Black	10	JR-201094BK-10	JR-201580BK-10
1" long nuts	White	10	JR-202120-10	JR-202119-10
	Black	10	JR-202118-10	JR-202121-10
Polypropylene washer		10	JR-201926-10	JR-201928-10

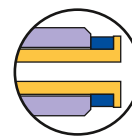**Union body for flanged connections**

1/4-28

- Through bore low pressure union
- For flanged 1/4-28 flat-bottom connections
- For 1/16" or 1/8" OD tubing
- Completely inert and biocompatible


Flanged fittings are not included.

Order PPS hex-head flanged nuts or color-coded fingertight flanged nuts separately.



Description	Qty/pkg	Product No.
Union body, polyamide	5	JR-060-5

**SPECS**

- For 1/16" and 1/8" OD tubing
- Material: Polyamide
- 1/4-28 threads

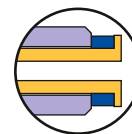**TECH TIP**

Use our Easy-Flange kit for forming flanges on tubing ..... 85

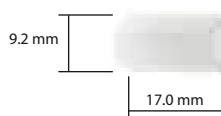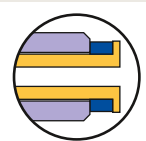

## PEEK unions

1/4-28

## SPECS

- For 1/16" and 1/8" OD tubing
- Materials  
Body: PEEK  
Nut: PPS  
Ferrules: ETFE
- 1/4-28 threads
- Pressure rating:  

OD	Bar	Psi
1/16"	105	1500
1/8"	35	500

 (Varies with tubing material and ID)

- For easy connection of low pressure solvent lines
- Low dead volume
- For 1/16" or 1/8" OD tubing
- Flat-bottomed 1/4-28 fitting details
- Completely inert and biocompatible

These unions facilitate the quick and reliable connection of two pieces of 1/16" or 1/8" OD tubing. Use with color-coded adapters for easy solvent line identification.

These are available as the body only, or complete with fittings.


A union with fittings includes:

- PEEK low pressure union
- 2 PPS flangeless hex-head nuts
- ETFE inverted ferrules
- Color-Sleeve adapters



	Bore (mm)	Tubing OD	Qty/pkg	Product No.
Union, with fittings	0.50	1/16"	1	JR-1064
	0.75	1/8"	1	JR-1068
	1.30	1/8"	1	JR-1065
Body only	0.50	1/16"	1	JR-064
	0.75	1/8"	1	JR-068
	1.30	1/8"	1	JR-065

Spare parts	Tubing OD	Qty/pkg	Product No.
PPS flangeless nuts	1/16"	10	JR-55070-10
	1/8"	10	JR-55071-10
Inverted ETFE ferrules	1/16"	10	JR-041-10
	1/8"	10	JR-051-10

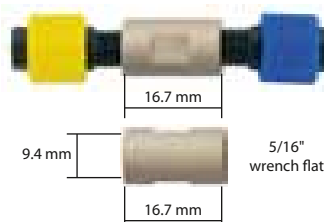


## TECH TIP


While our PEEK unions are typically used with inverted ferrules, collapsible ferrules or flanged connections can also be used.



**Standard:  
Inverted ferrule**

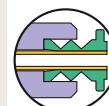

**Collapsible  
ferrule**

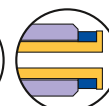

**Flanged  
ferrule**

## PEEK bulkhead unions

1/4-28

## SPECS

- For 1/16" and 1/8" OD tubing
- Materials  
Body: PEEK  
Nut: PPS  
Ferrules: ETFE
- 1/4-28 threads
- Pressure rating:  

OD	Bar	Psi
1/16"	105	1500
1/8"	35	500


 (Varies with tubing material and ID)

- Panel mounting (10 mm hole required, max 9 mm thick)
- For easy connection of low pressure solvent lines
- Low dead volume
- For 1/16" or 1/8" OD tubing
- Completely inert and biocompatible

Bulkhead unions are designed to be mounted through an instrument panel or on a bracket.

They are sold complete with fittings, and include:

- 2 PPS flangeless hex-head nuts
- 2 ETFE inverted ferrules
- 2 Color-Sleeve adapters

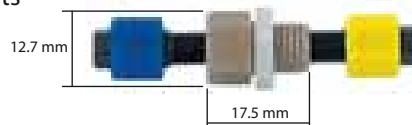

	Bore (mm)	Tubing OD	Qty/pkg	Product No.
Bulkhead union, with fittings	0.75	1/16"	1	JR-CBUMPK
	1.50	1/8"	1	JR-CBULPK

Spare parts	Tubing OD	Qty/pkg	Product No.
PPS flangeless nuts	1/16"	10	JR-55070-10
	1/8"	10	JR-55071-10
Inverted ETFE ferrules	1/16"	10	JR-041-10
	1/8"	10	JR-051-10



## TECH TIP

Use our impermeable tubing to avoid gas diffusion in the solvent lines ..... 13

## RELATED PRODUCTS

Color-coded fingertight flanged nuts ..... 43  
 Color-Sleeve adapters, in 12 colors ..... 17  
 High pressure PEEK unions ..... 26  
 PPS hex-head flanged nuts ..... 41

**PEEK reducing unions****NEW***Imperial to metric, 1/4-28*

- Connects 1/16", 1/8", 1 mm, 2 mm, and 3 mm OD polymeric tubing
- Completely inert and biocompatible

With these new reducing unions, metric tubing can be connected to imperial fittings on valves and LC systems.

A reducing union with fittings includes:

- PEEK low pressure reducing union
- 1 PPS flangeless hex-head nut, imperial size noted
- 1 PPS flangeless hex-head nut, metric size noted
- ETFE inverted ferrules
- Color-Sleeve adapters

Tubing OD	Bore (mm)	Qty/pkg	Product No.
Reducing union, with fittings			
1/32" to 1 mm	0.50	1	JR-0640
1/16" to 2 mm	0.75	1	JR-0680
1/8" to 3 mm	1.30	1	JR-0650

Spare parts	Tubing OD	Qty/pkg	Product No.
PPS flangeless nuts	1/16"	10	JR-55070-10
	1/8"	10	JR-55071-10
	1.0 mm/ 1/32"	10	JR-55073-10
	2.0 mm	10	JR-55072-10
Inverted ETFE ferrules	1/32"	10	JR-032-10
	1/16"	10	JR-041-10
	1/8"	10	JR-051-10
	1.0 mm	10	JR-031-10
	2.0 mm	10	JR-042-10
	3.0 mm	10	JR-049-10



3.0 mm




2.0 mm



1.0 mm

**Metric ferrules****SPECS**

- For 1/32"–1/8" OD tubing, flangeless connections
- Materials  
Body: PEEK  
Nuts: PPS  
Ferrules: ETFE
- 1/4-28 threads
- Pressure rating: ≤48 bar (700 psi)  
(Varies with tubing material and ID)

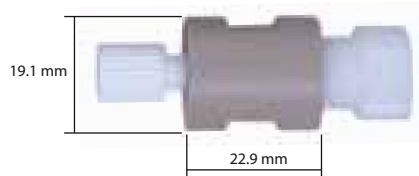**TECH TIP**

While our PEEK bulkhead unions are typically used with inverted ferrules, collapsible ferrules or flanged connections can also be used.

**RELATED PRODUCTS**

Color-Sleeve adapters, in 12 colors . . . . . 17

**PEEK reducing unions***Imperial, 1/2-20 to 1/4-28*

- Connects 1/4" and 1/8" OD tubing
- Low dead volume
- Completely inert and biocompatible

PEEK reducing unions are sold complete with fittings, and include:

- One 1/8" flangeless nut, white polypropylene
- One 1/8" ferrule, CTFE
- One 1/4" flangeless nut, polyacetal
- One 1/4" ferrule, CTFE

	Bore (mm)	Tubing OD	Qty/pkg	Product No.
Reducing union with fittings	1.5	1/4" to 1/8"	1	JR-CU4LPK

**Spare parts**

Tubing OD	Description	Qty/pkg	Product No.
1/8"	Flangeless nut, white propylene	1	JR-CFL-2W-S
	Ferrule, CTFE	1	JR-CFL-CB2KF-S
1/4"	Flangeless nut, polyacetal	1	JR-CFL-4D
	Ferrule, CTFE	1	JR-CFL-CB4KF-S

**SPECS**

- For 1/4" to 1/8" OD tubing, flangeless connections
- Materials  
Body: PEEK  
Nuts: 1/8" PP  
1/4" Polyacetal  
Ferrules: CTFE
- Threads: 1/2-20 to 1/4-28
- Pressure rating: <35 bar (500 psi)  
(Varies with tubing material and ID)

PEEK tees, crosses, and manifolds

1/4-28

**SPECS**

- For 1/16" and 1/8" OD tubing
- Materials
  - Body: PEEK
  - Nut: PPS
  - Ferrules: ETFE
- 1/4-28 threads
- Body sizes: see chart below
- Pressure rating:
 

OD	Bar	Psi
1/16"	105	1500
1/8"	35	500

 (Varies with tubing material and ID)

- For 1/16" or 1/8" OD tubing
- Flat-bottomed 1/4-28 fitting details
- Completely inert and biocompatible

These VICI Jour connectors for low pressure applications are machined from PEEK. The bodies have 1/4-28 female ports which accept any flanged or flangeless fittings for 1/8" or 1/16" OD tubing.

Choose smaller through-holes for applications where low flowrate or low dead volume is required. Fittings with 1.50 mm through-holes are preferred for high flowrate applications.

These are available as the body only, or complete with fittings.


The appropriate number of the following are included when fittings are ordered:

- PPS flangeless hex-head nuts
- ETFE inverted ferrules
- Color-Sleeve adapters

Sold individually.



				With fittings	Body only
		Tubing OD	Bore (mm)	Product No.	Product No.
Tees		1/16"	0.50	JR-6029	JR-6034
			0.80	JR-6030	JR-6032
		1/8"	1.50	JR-6031	JR-6033
Crosses		1/16"	0.50	JR-6039	JR-6044
			0.80	JR-6040	JR-6042
		1/8"	1.50	JR-6041	JR-6043
Manifolds	5 ports	1/16"	0.75	JR-C5M1PK	JR-B-C5M1PK
		1/8"	1.50	JR-C5M2PK	JR-B-C5M2PK
	9 ports	1/16"	0.75	JR-C9M1PK	JR-B-C9M1PK
		1/8"	1.50	JR-C9M2PK	JR-B-C9M2PK
<b>Spare parts</b>		Tubing OD	Qty/pkg	Product No..	
PPS flangeless nuts	1/16"	10	JR-55070-10		
	1/8"	10	JR-55071-10		
Inverted ETFE ferrules	1/16"	10	JR-041-10		
	1/8"	10	JR-051-10		

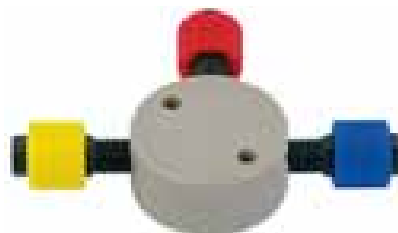


**TECH TIP**


While our PEEK tees and crosses are typically used with inverted ferrules, collapsible ferrules or flanged connections can also be used.



**Standard:  
Inverted ferrule**

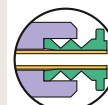

**Collapsible  
ferrule**

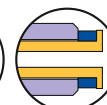

**Flanged  
ferrule**

**TECH TIP**

Use our impermeable tubing to avoid gas diffusion in the solvent lines ..... 13

**RELATED PRODUCTS**

Color-Sleeve adapters, in 12 colors ..... 17  
High pressure PEEK tees ..... 28  
PPS hex-head flanged nuts ..... 41

**BODY SIZES**

	Diameter	Thickness
Tees	25.0 mm	17.1 mm
Crosses	25.0 mm	17.1 mm
Manifolds		
5 ports	25.0 mm	17.1 mm
9 ports	25.0 mm	28.6 mm

**PEEK mixing tees**

1/4-28

- For 1/16" or 1/8" OD tubing
- Mountable
- Completely inert and biocompatible

A unique mixing results from the the turbulence created by the introduction angle of the solvents. The mixing is much more thorough than with a conventional tee.

Typically our inverted ferrules are used in these fittings, but the 1/4-28 flat-bottomed details will also accept collapsible ferrule or flanged fittings.

These are available as the body only, or complete with fittings. A mixing tee with fittings includes:

- 3 PPS flangeless hex-head nuts
- 3 ETFE inverted ferrules
- 3 Color-Sleeve adapters

Sold individually.


**SPECS**

- For 1/16" and 1/8" OD tubing
- Materials  
Body: PEEK  
Nut: PPS  
Ferrules: ETFE
- 1/4-28 threads
- Pressure rating:  

OD	Bar	Psi
1/16"	105	1500
1/8"	35	500

 (Varies with tubing material and ID)
- Mixing volume  

Bore	Volume
0.75 mm	25 µl
1.50 mm	100 µl

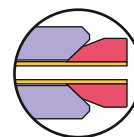

	Tubing OD	Bore (mm)	Product No..
Mixing tee, with fittings	1/16"	0.75	JR-CM1XPK
	1/8"	1.50	JR-CM2XPK
Body only	1/16"	0.75	JR-B-CM1XPK
	1/8"	1.50	JR-B-CM2XPK
<b>Spare parts</b>	<b>Tubing OD</b>	<b>Qty/pkg</b>	<b>Product No..</b>
PPS flangeless nuts	1/16"	10	JR-55070-10
	1/8"	10	JR-55071-10
Inverted ETFE ferrules	1/16"	10	JR-041-10
	1/8"	10	JR-051-10

**RELATED PRODUCTS**

Color-Sleeve adapters,  
in 12 colors . . . . . 17




**Polypropylene plugs**

1/4-28 or 10-32

**SPECS**

- Material: Polypropylene
- 10-32 or 1/4-28 threads
















Sold in packages of 10.

Description		Product No.
10-32 male		JR-072450-10
1/4-28 male		JR-072350-10
1/4-28 male, flat bottom		JR-201999-10


**Polypropylene adapters****SPECS**

- Material: Polypropylene

Sold in packages of 10.

Inlet		Outlet			Product No.
Male	Luer	Barbed	1/16"		JR-070211-10
			3/32"		JR-070212-10
			1/8"		JR-070213-10
	1/4-28	Male	10-32		JR-070224-10
			1/4-28		JR-070225-10
Male	1/4-28	Barbed	1/16"		JR-071123-10
			1/8"		JR-071323-10
Male	10-32	Barbed	1/16"		JR-071124-10
			1/8"		JR-071324-10
Female	Luer	Barbed	1/16"		JR-070111-10
			3/32"		JR-070112-10
			1/8"		JR-070113-10
	10-32	Female Luer			JR-070101-10
					JR-070124-10
					JR-070123-10

Sold individually.

Flat bottom 1/4-28 Barbed 1/8" ID  JR-075013**RELATED PRODUCTS**

High pressure  
PEEK adapters . . . . . 29  
High pressure plugs . . . 24

## Micro valves for GC and LC

1/4-28

- Compact 1" design
- Convenient panel mount, hole diameter 8 mm (0.3125"), 6 mm maximum thickness
- Biocompatible and inert
- Variety of configurations

Simplify your liquid or gas handling application with a VICI Jour sampling micro valve. The unique design of the fitting detail allows a leak-free seal with no potential for rotor damage from overtightening. Internal parts are PEEK and PTFE.

Includes flangeless PPS fittings and CTFE collapsible ferrules for 1/8" OD tubing. Sold individually.



## SPECS

- For 1/16" and 1/8" OD tubing
- Materials
  - Wetted parts: PEEK,
  - Body: POM
  - Nuts: PPS
  - Ferrules: CTFE
- 1/4-28 threads
- 1.50 mm bore
- Pressure rating: see notes at left

MICRO VALVE  
PRESSURE RATINGS

- 200 psi gas  
Tested with helium using 1/16" OD x 0.84 mm ID PTFE tubing with PEEK collapsible ferrules
- 500 psi liquid  
Tested with isopropanol using 1/8" OD x 0.75 mm ID PTFE tubing with CTFE collapsible ferrules

## TECH TIP

While micro valves are typically used with collapsible ferrules, inverted ferrules or flanged connections can also be used.

Standard:  
Collapsible ferruleInverted  
ferruleFlanged  
ferrule

## TECH TIP

Use our impermeable tubing to avoid gas diffusion in the solvent lines . . . . . 13

## RELATED PRODUCTS

Color-Sleeve adapters, in 12 colors . . . . . 17



High pressure PEEK tees . . . . . 28  
PPS hex-head flanged nuts . . . . . 41

		Product No.
"T" flow path	3 ports	JR-660100
	4 ports	JR-660110
180° flow path	2 ports	JR-660200
	4 ports	JR-660210
90° flow path	2 ports	JR-660300
	3 ports	JR-660310
	4 ports	JR-660320

Spare parts	Material	Tubing OD	Qty/pkg	Product No.
Nuts	PPS	1/16"	10	JR-55050-10
		1/8"	10	JR-55051-10
Collapsible ferrules	PEEK	1/16"	10	JR-CFL-CB1PK
	CTFE	1/8"	10	JR-CFL-CB2KF

## "T" flow path



## 180° flow path



## 90° flow path



## Peristaltic tubing adapter set

1/4-28

- For connecting soft tubing to PTFE and PTFE-like tubing
- Works with peristaltic tubing up to 4.50 mm OD
- For 1/16" and 1/8" PTFE and PTFE-like tubing
- Biocompatible
- Low dead volume

This "one-fits-all" set consists of everything needed to connect 1/16" or 1/8" OD PTFE and PTFE-like tubing to soft tubing (such as silicone), up to 4.5 mm in OD, minimum 0.5 mm ID.


The adapter includes:

- 3 different bushings (for the different soft tubing ODs)
- 2 PPS nuts (one each for 1/16" and 1/8" OD tubing)
- 2 ETFE inverted ferrules (one each for 1/16" and 1/8" OD tubing)
- 2 Color-Sleeve adapters

Sold individually.

## SPECS

- For 1/16" and 1/8" OD tubing
- Materials
  - Body: PEEK
  - Bushings: PPS
  - Nut: PPS
  - Ferrules: ETFE
- 1/4-28 threads
- 1.50 mm bore
- Body size: 28 mm long 12 mm diameter.

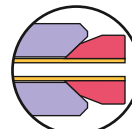

Product No.

PEEK adapter for peristaltic to PTFE tubing connection

JR-PTAS



# Filters

## AND MOBILE PHASE FILTERS

Filters trap foreign particulates from sample and mobile phase streams before they can damage valuable columns and instruments.

### Filters index

#### Low pressure filters

In-line filters .....	page 50
Mobile phase filters	
Last Drop .....	51
Economy Last Drop .....	51
Last Drop biocompatible glass filter .....	52
No-Met biocompatible .....	52
Filter/spargers	
Last Drop .....	53
Stainless steel .....	53

#### High pressure filters

PEEK high pressure filters	
In-line filter kit .....	54
Precolumn filter kit .....	55
Sure-Guard .....	54
Stainless UHPLC filters .....	56
Stainless inline filter .....	57

#### In-line filter cartridges

*Low pressure*

##### SPECS

- Materials
- Cartridge: PTFE/CTFE
- Filter screen: 316 stainless
- 1/4-28 detail
- Dimensions
- OD 5.20 mm
- Thickness 2.03 mm
- Bore 0.8 mm
- Filter surface diameter 2 mm
- Tolerances
- ±0.05 mm (.002")

- Easy to replace
- Compact design

- Fit all 1/4-28 fitting details
- Three porosity options

These convenient-to-use filters can be simply dropped into any 1/4-28 fitting detail, such as in a union. The filters are constructed of PTFE and CTFE, with a Type 316 stainless low-pressure-drop screen. The inner design of the cartridge ensures the equal distribution of the solvent to the screen.

Pore size	Maximum flow rate *	Qty/pkg	Product No.
2 µm	30 ml/min	5	JR-CFE-S2-5
10 µm	30 ml/min	5	JR-CFE-S10-5
75 µm	30 ml/min	5	JR-CFE-S75-5



#### NOTES

\* Flowrate measured with methanol/water (1:1), ultrasonic degassed, helium sparging to prevent regassing.

#### TECH TIP

We recommend using these in-line filter cartridges in combination with our low pressure fittings:

PEEK unions .....	44
Tees and crosses .....	46
Manifolds .....	46
Micro valves .....	49

#### DISCLAIMER

The maximum holding pressure for any type of connections involving tubing and a ferrule varies considerably with the tubing material, the ferrule material, the clearance between tubing OD and ferrule ID, and the shape of the fitting detail.





**Last Drop mobile phase filters**

Low pressure

- No loss of mobile phase
- Biocompatible PTFE frits or stainless steel frits
- Three porosity options
- Two connector types: stepped tubing and 1/4-28 fitting



The Last Drop mobile phase filter allows more analyses per batch of mobile phase and helps reduce hazardous waste. The flat filter element sits parallel to the bottom of the reservoir, allowing the Last Drop to filter all but the last 2% of the mobile phase from the reservoir without drawing air into the system. Compare this with conventional cylindrical filters that can begin to draw air into the system when nearly 10% of the solvent remains in the reservoir.

The Last Drop mobile phase filter consists of a 316 stainless steel, PTFE 2.5 µm, or a hydrophobic PE filter element pressed into an inert PTFE housing. The top of the housing has a 1/4-28 nut and ferrule or a stepped PEEK fitting connector which slips into 1.5, 2.2, or 3.5 mm ID pump inlet lines.

**SPECS**

- Materials
  - Body: PTFE
  - Frits: PTFE/ 316 SS/PE
- Stepped tubing version is for 1.5, 2.2, and 3.5 mm ID tubing
- Fitting version (for 1/8" OD tubing) includes PEEK 1/4-28 nut with ETFE ferrule

Connector type:			Stepped tubing connector	Fitting connector
Filter material	Pore size	Maximum flow rate *	Product No.	Product No.
PTFE	2.5 µm	1.2 ml/min	JR-9000-0520	JR-9000-0520F
	5 µm	2.6 ml/min	JR-9000-0521	JR-9000-0521F
	10 µm	3.5 ml/min	JR-9000-0522	JR-9000-0522F
Polyethylene	10 µm	11 ml/min	JR-9000-0522H	JR-9000-0522HF
Stainless	2 µm	28 ml/min	JR-9000-0530	JR-9000-0530F
	5 µm	30 ml/min	JR-9000-0531	JR-9000-0531F
	10 µm	30 ml/min	JR-9000-0532	JR-9000-0532F

**NOTES**

\* Flowrate measured with methanol/water (1:1), ultrasonic degassed, helium sparging to prevent regassing.

**TECH TIPS**

- Economy Last Drop mobile phase filters are easy to replace - easy to clean. Use ultrasound or replace with a new filter.
- Connect the Economy Last Drop filter with 1/8" OD/ 1.59 mm OD tubing:
  - PFA tubing ..... 9
  - PTFE tubing ..... 12
- We recommend metal-free PTFE or glass filters for sensitive biochromatography applications where metal surfaces may corrode or interact with samples.

**RELATED PRODUCTS**

Last Drop mobile phase filter/spargers ..... 53

**Economy Last Drop mobile phase filters**


Low pressure

- Very competitively priced
- Biocompatible PTFE frits or stainless steel frits
- Three porosity options
- Two diameters



Economy mobile phase filters come in two body diameters, for various dimensions of bottle necks. The filter body is made of PTFE, which is resistant to virtually all common mobile phases.

Frit material	Body OD	Pore size	Maximum flow rate *	Product No.
PTFE	19 mm	2.5 µm	1.2 ml/min	JR-4676-2.5TF
		5 µm	2.6 ml/min	JR-4676-5TF
		10 µm	3.5 ml/min	JR-4676-10TF
	22 mm	2.5 µm	1.2 ml/min	JR-4677-2.5TF
		5 µm	2.6 ml/min	JR-4677-5TF
		10 µm	3.5 ml/min	JR-4677-10TF
Stainless steel	19 mm	2 µm	30 ml/min	JR-4676-2
		10 µm	30 ml/min	JR-4676-10
	22 mm	2 µm	30 ml/min	JR-4677-2
		10 µm	30 ml/min	JR-4677-10

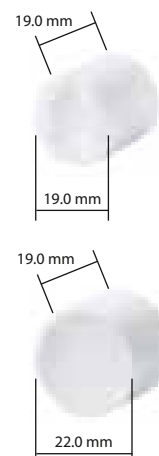

## FILTERS

Last Drop biocompatible glass mobile phase filters *Low pressure*

## SPECS

- Materials
  - Body: Glass
  - Frit: Glass
- For use with 1/8" tubing

- The solution for biocompatibility and high flowrates
- Easy-push tubing connection

Because it has a glass frit, the flowrate can be much higher with this filter than with PTFE mobile phase filters. The glass foot allows withdrawal of nearly all of the mobile phase – less than 2% is left in the bottle.

Includes 1/8" tubing connector.

Frit material	Pore size	Maximum flow rate *	Product No.
Glass	1.0 - 1.6 µm	30 ml/min	JR-9000-0520B
	10-16 µm	50 ml/min	JR-9000-0526G
	40-100 µm	200 ml/min	JR-9000-0528G

## Replacement connectors

For tubing	Size	Qty/pkg	Product No.
1/8" OD	6.8 mm	1	JR-9000-0525GC

No-Met biocompatible mobile phase filters *Low pressure*

## SPECS

- Materials
  - Body: PTFE/ Polyethylene
  - Adapter: PEEK
  - Fitting: PEEK
  - Ferrule: ETFE
- For use with 1/8" tubing

- Very competitively priced
- Inert and biocompatible
- Replacement filters
- Also recommended for IC

Stainless steel in the flowpath is not acceptable in a growing number of applications involving the separation of biomolecules. High salt buffers can corrode stainless steel, and the metal ions released from metallic filters may contaminate or otherwise react with the biomolecules of interest.

The No-Met polyethylene filter is designed for these applications, with inert polymeric fittings for 1/8" tubing and a 20 µm filter effectively eliminating metal contamination from the fluid path. Use them for IC and biochromatography applications.


No-Met filters can be used at flow rates up to 500 ml/min, measured with methanol/water (1:1), ultrasonically degassed. Flow rates can vary with solvent and tubing ID.

The economy version can easily be slipped over 1/8" OD tubing, with no fitting required.

Connection	Filter material	Pore size	Maximum flow rate *	Product No.
1/4-28 nut and ferrule	PTFE	5 µm	2.2 ml/min	JR-32171
	Polyethylene	<20 µm	500 ml/min	JR-32178
Economy slip-on connection	Polyethylene	5 µm	300 ml/min	JR-32174

## Replacement element

Filter material	Pore size	Maximum flowrate *	Product No.
PTFE	5 µm	2.2 ml/min	JR-32172
Polyethylene	<20 µm	500 ml/min	JR-32179

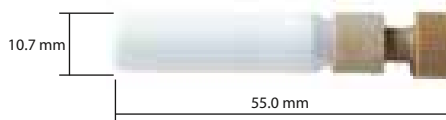

## NOTES

\* Flowrate measured with methanol/water (1:1), ultrasonic degassed, helium sparging to prevent regassing.

## TECH TIPS

- Because they are hydrophobic, No-Met filters may initially require some priming with methanol or acetonitrile.
- Connect the No-Met filter with 1/8" OD/ 1.59 mm OD tubing:
  - PFA tubing ..... 9
  - PTFE tubing ..... 12



### Last Drop filter/sparger


Low pressure

The Last Drop filter/sparger combines filtration and sparging in a single unit. The PTFE housing contains a mobile phase filter with a stainless steel, PTFE, or hydrophobic PE filter element.

Spargers have a porosity of 10 microns.

The filter/sparger includes:

- PEEK tripod connector which slips into 1.5, 2.2, or 3.5 mm ID solvent lines
- 1/4-28 nut and ferrule for the sparging line

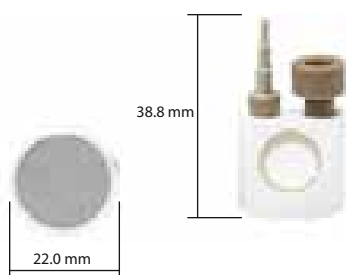

Filter material	Pore size	Maximum flow rate *	Product No.
PTFE	2.5 µm	1.2 ml/min	JR-9000-0602
	5 µm	2.6 ml/min	JR-9000-0603
	10 µm	3.5 ml/min	JR-9000-0604
Polyethylene	10 µm	11 ml/min	JR-9000-0604H
Stainless	2 µm	28 ml/min	JR-9000-0640
	5 µm	30 ml/min	JR-9000-0641
	10 µm	30 ml/min	JR-9000-0642

#### SPECS

- Materials
  - Body: PTFE
  - Frits: PTFE/ 316 SS/PE
- Stepped tubing connector: PEEK
- Nut: PEEK
- Ferrule: ETFE
- Stepped tubing connector for 1.5, 2.2, and 3.5 mm ID tubing
- PEEK 1/4-28 nut with ETFE ferrule for 1/8" OD tubing

### Stainless steel mobile phase filter/helium sparger

Low pressure

- Ideal for helium sparging
- Four porosities
- Fitting and pipe connectors

VICI Jour mobile phase filters protect HPLC systems from small particles in the mobile phase. These filters are made from SS316 with PEEK or PTFE connectors and are suitable for most solvents. The complete line has versions for both analytical and preparative applications.


VICI Jour helium spargers (2 µm versions) offer an inexpensive way to prepare and maintain mobile phases free of dissolved gases. Connected to a regulated supply of helium gas (0-400 mL/min.) they effectively remove dissolved oxygen, nitrogen, and other atmospheric gases from the mobile phase.

#### SPECS

- Materials
  - Body: SS316
  - Pipe: SS316
- Pipe adapter: PTFE
- Fitting adapter: PEEK
- Fittings: PEEK
- Ferrules: ETFE



Fitting type	Tubing/ pipe size	Pore size	Maximum flow rate *	Product No.
1/4-28 fitting	1/16"	2 µm	35 ml/min	JR-367016-2
		10 µm	35 ml/min	JR-367016-10
		20 µm	35 ml/min	JR-367016-20
	1/8"	2 µm	35 ml/min	JR-367008-2
		10 µm	100 ml/min	JR-367008-10
		20 µm	120 ml/min	JR-367008-20
Pipe	1/8"	2 µm	50 ml/min	JR-3675-2
Pipe	1/8"	2 µm	95 ml/min	JR-3678-2
		25 µm	100 ml/min	JR-3678-25

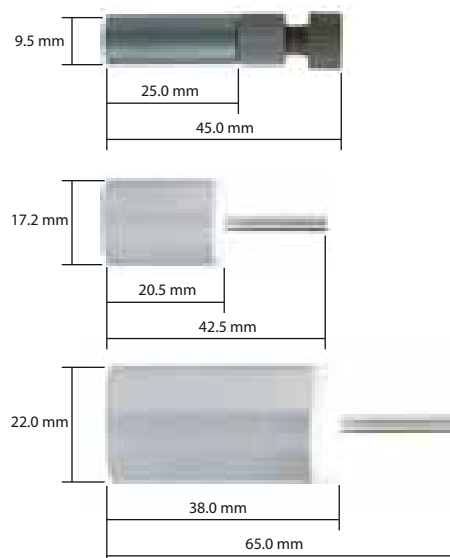

#### TECH TIPS

- We recommend metal free PTFE filters for sensitive biochromatography applications where metal surfaces may corrode or interact with samples.
- We recommend our impermeable tubing to prevent "regassing" of helium-degassed solvents.

#### RELATED PRODUCTS

Last Drop mobile phase filters ..... 51

FILTERS

PEEK in-line filter kit

HPLC

**SPECS**

- For 1/16" OD tubing
- Materials  
Body: PEEK  
Frit: Titanium, Polyethylene
- 0.25/0.4 mm bore
- 10-32 threads
- Pressure rating: <350 bar (< 5000 psi)

- Replaceable metal-free polyethylene frits and biocompatible titanium frits
- 100% biocompatible
- Minimal hold-up volume
- Different porosities

The PEEK in-line filter traps fines and other particles from samples and mobile phases before they damage valuable instruments and columns.

This female-to-female design is made entirely of PEEK, with PEEK-encased titanium or polyethylene filter elements, for biocompatibility and chemical resistance. The design has virtually no hold-up volume and can be used in analytical applications with virtually no band broadening or loss of efficiency.

Sold individually.

Frit material	Pore size	Internal volume (µL) *	Maximum flow rate **	Product No.
Titanium	0.5 µm	9.50	25 ml/min	JR-68247
	2 µm	11.85	30 ml/min	JR-68253
	5 µm	12.7	30 ml/min	JR-68256
	10 µm	13.2	30 ml/min	JR-68257
Polyethylene	10 µm	13.9	30 ml/min	JR-68257PE



Sure-Guard in-line filter kit

HPLC

**SPECS**


- Materials  
Body: PEEK  
Frit: Stainless, Titanium
- 0.4 mm bore
- 10-32 female to 10-32 male threads
- Pressure rating: <350 bar (< 5000 psi)

- Disposable filter guard
- Minimal hold-up volume
- Titanium (biocompatible) and stainless steel frits, permanently pressed in place
- Two porosity options


The VICI Jour Sure-Guard disposable in-line filter offers an easy and inexpensive way to protect valuable columns against fines and particles. It is easily connected directly to any column, with an inlet for 1/16" OD tubes and 10-32 threads. The VICI Jour Sure-Guard can be changed in seconds without tools.

Sold in packages of three. Waters type filter kits are sold individually.

Frit material	Pore size	Internal volume *	Maximum flow rate **	Product No.	Product No. Waters type
Stainless steel	0.5 µm	0.45 µL	30 ml/min	JR-0611-SS05-3	JR-0641-SS05
	2 µm	0.61 µL	30 ml/min	JR-0611-SS2-3	JR-0641-SS2
Titanium	0.5 µm	0.45 µL	30 ml/min	JR-0611-TI05-3	JR-0641-TI05
	2 µm	0.61 µL	30 ml/min	JR-0611-TI2-3	JR-0641-TI2

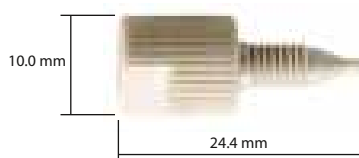

**NOTES / TECH TIPS**  
See facing page.

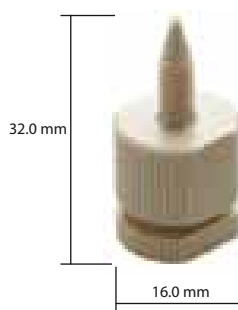**PEEK precolumn filter kit**

HPLC

- 100% biocompatible
- Minimal hold-up volume
- Replaceable metal-free polyethylene frits and biocompatible titanium frits encased in PEEK
- Different porosities

The VICI Jour precolumn filter kit has a standard 10-32 fitting and fits direct to most column types without introducing additional dead volume.

The filter kit protects expensive columns against fine particles, which may otherwise accumulate on the column frit and lead to split peaks and high backpressure.

Sold individually.

**SPECS**

- Materials
  - Body: PEEK
  - Frit: Titanium, Polyethylene
- 0.25/0.4 mm bore
- 10-32 threads
- Pressure rating: <350 bar (< 5000 psi)

Frit material	Pore size	Internal volume (μL) *	Maximum flow rate **	Product No.
Titanium	0.5 μm	11.5	30 ml/min	JR-68258
	2 μm	13.85	30 ml/min	JR-68262
	5 μm	14.7	30 ml/min	JR-68263
	10 μm	15.2	30 ml/min	JR-68264
Polyethylene	10 μm	15.9	30 ml/min	JR-68264P

**Replacement frits**

for PEEK in-line and precolumn filter kits

Replacement frits for PEEK precolumn filter kit and PEEK in-line filter kit (see facing page).

Frit material	Pore size	Internal volume (μL) *	Qty/pkg	Product No.
Titanium	0.5 μm	6.87	5	JR-1125-05P-5
	2 μm	9.24	5	JR-1125-2P-5
	5 μm	10.04	5	JR-1125-5P-5
	10 μm	10.56	5	JR-1125-10P-5
Polyethylene	10 μm	11.31	5	JR-1151-10P-5

**NOTES**

\* The internal volume is the hold-up volume between the end of the connecting tube and the filter 10-32 connector outlet, including the frit volume.

\*\* Flowrate measured with methanol/water (1:1), ultrasonic degassed, helium sparging to prevent regassing.

**TECH TIP**

We recommend 2 μm frits for columns with 5 μm or larger particles, and 0.5 μm frits for smaller particles.

**NOTE**

Polyethylene frits are for single use only.

**MORE INFORMATION**

Calculation of  
frit volume ..... 95  
PEEK-encased  
stainless frits..... 74

## Stainless steel precolumn filter

NEW

UHPLC

## SPECS

- Materials
  - Body: 316 SS
  - Frit: PEEK-encased SS
- 10-32 threads
  - Female for inlet
  - Male for outlet
- Pressure rating with ZN1 nut and ZF1S6 ferrule: 1375 bar (20,000 psi)

- UHPLC
- Minimized hold-up volume
- Three porosity options
- Fewer connections compared to in-line filters
- Works with 1/16" column connections from all manufacturers




VICI Jour's new stainless steel precolumn filter is based on our well-known in-line filter. The design of this precolumn filter reduces the number of connections compared to a standard in-line filter, increasing the stability and reducing the potential of leakage.

The small bore sizes (0.15 mm in the nut/body and 0.13 mm in the tubing) guarantee a very small (virtually zero) hold-up volume.

With minimal dead-volume and available with three porosities, these filters are the best column protectors.

Bore size	Pore size	Internal volume	Product No.
0.13/0.15 mm	0.5 µm	0.61 µl	JR-69230-05
	2 µm	0.66 µ	JR-69230-2
	5 µm	0.68 µ	JR-69230-5

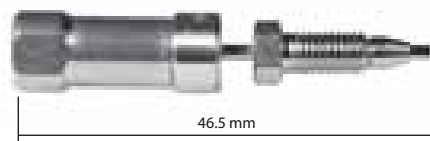

## Inlet nuts and ferrules

Includes one 1/16" stainless nut and one 1/16" stainless ferrule.

Pore size	Product No.
0.5 µm	JR-69330-05
2 µm	JR-69330-2
5 µm	JR-69330-5

## Replacement PEEK-encased SS frits

Pore size	Frit volume	Qty/pkg	Product No.
0.5 µm	0.13 uL	5	JR-1110-05P-5
2 µm	0.18 uL	5	JR-1110-2P-5
5 µm	0.20 ul	5	JR-1110-5P-5



**Stainless steel in-line filter****UHPLC**

- UHPLC (up to 1375 bar / 20,000 psi)
- Minimized hold-up volume
- Two bores for analytical and semi-prep flow rates
- Three porosity options




The stainless steel in-line filter has virtually no hold-up volume and can be used in various applications without significant band broadening or loss of efficiency.

The 0.25 mm bore versions are intended for analytical use between autosampler/injection valve and column. The 0.75 mm bore versions are suitable for higher flow rates in semi-prep systems or for solvent supply from pump to autosampler.

With minimal dead volume and available frit porosities of 0.5, 2 or 5 µm, these in-line filters are the ideal column protectors

Sold individually.

Bore size	Pore size	Internal volume *	Product No.
0.25 mm	0.5 µm	0.2 µL	JR-68230-05
	2 µm	0.3 µL	JR-68230-2
	5 µm	0.3 µL	JR-68230-5
0.75 mm	0.5 µm	1.7 µL	JR-68231-05
	2 µm	1.9 µL	JR-68231-2
	5 µm	1.9 µL	JR-68231-5

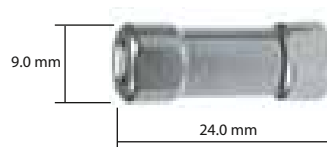**SPECS**

- Materials
  - Body: 316 SS
  - Frit: PEEK-encased 316 stainless
- 0.25 or 0.75 mm bore
- 10-32 female to 10-32 male threads
- Pressure rating with ZN1 nut and ZF1S6 ferrule: 1375 bar (20,000 psi)

**Nuts and ferrules**

Includes one 1/16" stainless nut and one 1/16" stainless ferrule.

Bore size	Pore size	Product No.
0.25 mm	0.5 µm	JR-68330-05
	2 µm	JR-68330-2
	5 µm	JR-68330-5
0.75 mm	0.5 µm	JR-68331-05
	2 µm	JR-68331-2
	5 µm	JR-68331-5

**Replacement PEEK-encased SS frits***for filter above*

Bore size	Pore size	Frit volume	Qty/pkg	Product No.
0.25 mm	0.5 µm	0.13 µL	5	JR-1110-05P-5
	2 µm	0.18 µL	5	JR-1110-2P-5
	5 µm	0.20 µL	5	JR-1110-5P-5
0.75 mm	0.5 µm	0.55 µL	5	JR-1111-05P-5
	2 µm	0.74 µL	5	JR-1111-2P-5
	5 µm	0.80 µL	5	JR-1111-5P-5
	10 µm	0.87 µL	5	JR-1111-10P-5

**NOTES**

\*The internal volume is the hold-up volume between the end of the connecting tube and the filter 10-32 connector outlet, including the frit volume.

# Safety products

## AND MOBILE PHASE HANDLING SYSTEMS

- Save cost on solvents
- Avoid hazardous air in the laboratories
- Increase reproducibility of chromatographic analyses

VICI strives to provide customers with the best solutions for their problems, making the work better and/or easier with respect to environment and health.

## Safety products index

### VICI safety filters

Air inlet valves .....	page 59
Air outlet valve.....	59
Exhaust filters.....	60-61
Carbon dioxide filter .. <b>NEW</b> .....	71
Humidity filter .....	<b>NEW</b> 71

### Caps

Safety cap kits .....	62
Safety caps .....	
DIN50 threads .....	63
GL45 threads .....	63
VICI cap kits.....	64
VICI caps .....	
GL38 threads .....	65
DIN50 threads .....	65
GL45 threads .....	64
S40 threads..... <b>NEW</b> .....	65

### Waste caps

GL45 threads .....	<b>NEW</b> 66
B53 threads .....	68
B83 threads .....	68
GLS80 threads.....	64
DIN60 threads .....	67
S51 threads.....	67
S61 threads.....	67
Eco cap .....	69
Eco cap kit.....	69

### Safety-coated glass bottles .....

### Solviflex tubing .....

Flexible coiled PFA tubing expands and contracts, giving flexibility to move mobile phase bottles on the bench.







### Safety air inlet valves

- No evaporation of volatile compounds into lab air
- Better composition stability of the mobile phases
- Less shift of retention times
- No bacterial contamination
- Fit any VICI cap or VICI safety cap

The air inlet valve has a check valve which blocks vapors of the solution in a bottle, preventing the vapors from escaping into the lab. The mobile phase composition remains stable, resulting in constant retention times.

The check valve regulates the pressure when the negative pressure reaches 1.5 psi or 0.1 bar. This small negative pressure prevents evaporation of the liquid (lower vapor pressure) and still guarantees the stability of the flow rate. A negative pressure which is too high can influence the aspiration of a constant flow.

Several tests have indicated that the 0.2 µm cellulose filter will keep mobile phase liquid free of bacteria and dust particles over time. If you use a filter and still have bacteria growth in your liquid, it is likely due to the presence of bacteria in the deionized water used.

Air inlet valves fit on any VICI safety cap or VICI cap.

Description	Product No.
With 4 mm filter	JR-S-20003
With 15 mm filter	JR-S-20009
Without filter	JR-S-20006

#### Spare parts

Replacement cellulose filter	
0.2 µm x 4 mm	JR-S-20007
0.2 µm x 15 mm	JR-S-20008



#### SPECS

- Connections  
Valve to cap:  
1/4-28  
Filter:  
Female Luer
- Materials  
Housing:  
Polypropylene  
Filter: Cellulose  
O-ring: NBR  
(O-ring is not a wetted part)



### Safety air outlet filters

The VICI safety air outlet filter regulates the pressure during helium sparging.

Air outlet filters fit on any VICI safety cap or VICI cap.

Description	Product No.
Safety air outlet filter	JR-S-20004



#### SPECS

- Connections  
Valve to cap:  
1/4-28  
Filter:  
Female Luer
- Materials  
Housing:  
Polypropylene  
Filter: Cellulose  
O-ring: NBR  
(O-ring is not a wetted part)

#### TECH TIPS

- The air inlet valve is not suitable when helium sparging is used.
- We recommend that safety air inlet valves and safety air outlet valves be replaced every 6-7 months to guarantee correct operation of the check valve. The lifetime of the check valve and the filter depend on the intensity of use and the quality of the ambient air.

#### RELATED PRODUCTS

Safety caps ..... 62-63  
VICI caps ..... 64-65


### SAFETY PRODUCTS

- Prevents pollution of ambient air with volatile compounds
- Large capacity and long lifetime
- Break-through detector available, alerting users to saturation of the filter absorbent
- Universal 1/4-28 connection – fit all VICI safety caps and VICI caps
- Adapter to competitor products available

Our improved exhaust filter is the center of the VICI Jour safety program, ridding ambient air of solvent vapors which could adversely affect the health of lab workers. When placed on the waste bottle of an HPLC system, the specially-treated absorbent filters out nearly all organic solvent vapor evaporating out of the waste solution.

Our new absorbent has an increased capacity over our previous filter and definitely over a competing brand. Excellent capacity leads to longer lifetime of the exhaust filter: up to 10 months\* for methanol, and up to 8 months\* for acetonitrile. The new design also has improved stability at the connecting threads.

The optional break-through detector recognizes trace amounts of acetone, acetonitrile, ethanol, ethyl acetate, methanol, methylene chloride and THF. When the absorbent is saturated, break-through is indicated by a color change in the detector from orange/pink to dark purple. A small charcoal filter absorbs the vapors after they pass the detector.

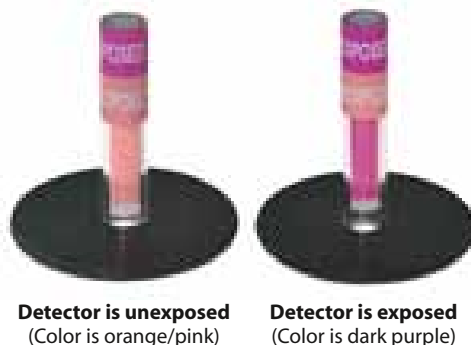

\*Based on results of our tests (water/organics = 50/50; flow rate 1.0 ml/min during the day, 0.5 ml/min overnight, 0.1 ml/min during weekends).

### Safety exhaust filters

1/4-28

#### SPECS

- Connections: 1/4-28
- Materials
  - Filter housing: POM
  - Detector housing: PC
  - Frits: Polyethylene
  - O-ring: NBR (O-ring is not a wetted part)

Description	Product No.
Exhaust filter	
With detector	JR-S-20005
Without detector	JR-S-20002
Adapter	JR-S-40004
Connects	
VICI safety exhaust filter to SCAT® waste caps and A.I.T. smart waste caps	



#### NOTE

Disposal of the saturated filters must be done according to the local regulations for chemical waste.


#### TECH TIPS

- Larger capacity = longer lifetime = lower costs
- As a general rule, halving the flow rate doubles the lifetime at constant composition and temperature. Doubling the flow rate halves it.

#### RELATED PRODUCTS


Carbon dioxide filter... 71  
Humidity filter... 71  
Safety caps... 62-63  
VICI caps... 64-65

See additional information about VICI safety exhaust filters on the facing page.

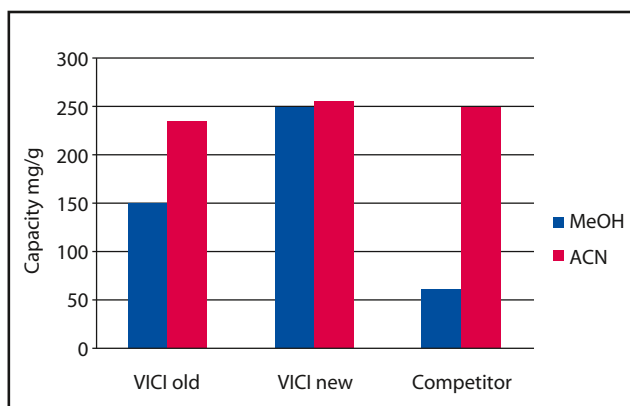


### Optimized absorbent of safety exhaust filters

In our continuous work to improve our products, VICI has optimized the absorbent for increased capacity.

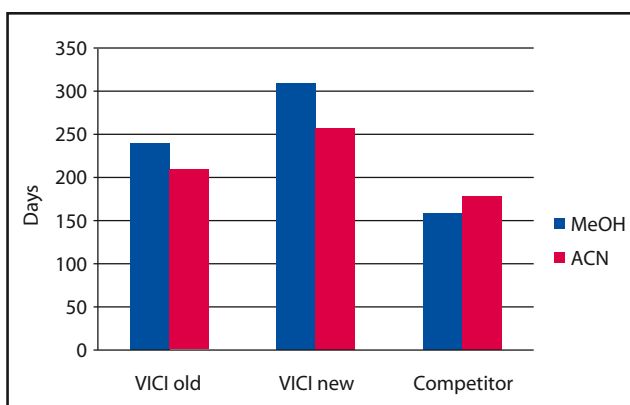

The excellent capacity of our new absorbent leads to longer lifetimes – up to 10 months for methanol and 8 months for acetonitrile.\*

\*Based on results of our tests (water/organics = 50/50; flow rate 1.0 ml/min during the day, 0.5 ml/min overnight, 0.1 ml/min during weekends).

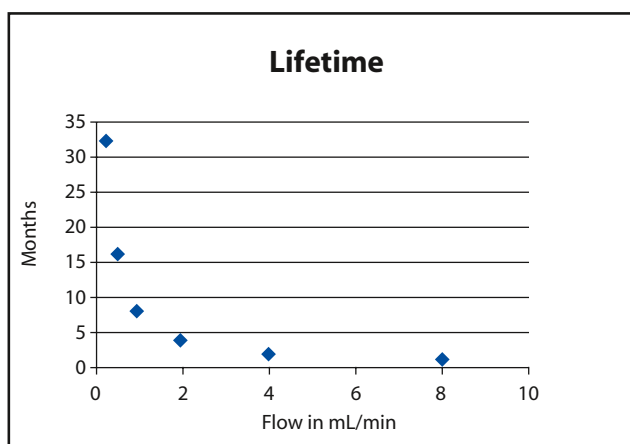

### Safety exhaust filter lifetime

The lifetime of the filter depends on:

- Composition of the mobile phase
- Flow rate
- Ambient temperature: a higher temperature leads to an increased vapor pressure of the volatile compound

The theoretical lifetime of an exhaust filter, depending on the flow rate, is shown here.

#### TECH TIP

Acetonitrile and water form an "azeotropic" mixture at 85% acetonitrile. Due to this physical property, the vapor of an acetonitrile/water mixture always corresponds to the 85/15 mixture. This explains the non-linearity of the lifetime for acetonitrile compared to methanol.

#### ADDITIONAL APPLICATIONS

Other applications for safety filters are possible, e.g. an exhaust filter for the GC split vent.

Contact Tech Support for more information.

Safety caps have built-in stopcocks to close off unused ports or to provide a convenient cut off of the liquid stream when tubing is removed.

- Stopcocks rated up to 72 psi gas at 25°C
- Prevent evaporation of volatile compounds
- Ideal for delivery of mobile phase or collection of waste LC solutions
- 2, 3, or 4 ports with 1/4-28 threads for 1/8" or 1/16" tubing
- Optional bottom-sealing barbed adapter to connect 1/8" (3.2 mm) ID tubing (shown at right – order separately from "Spare Parts")



## Safety cap kit

GL45 threads

### SPECS

- Connections: 1/4-28
- Materials:
  - Closure: PP (polypropylene)
  - Insert: PTFE
  - Stopcock: CTFE
  - O-ring: EPDM (alternative: FEP-coated)
  - Nuts: PPS
  - Ferrules: ETFE
  - Sleeves: PP
  - Bottle: Borosilicate 3.3, safety coated
  - Filter: Polyethylene

In addition to a GL45 safety cap with color-coded fittings, this kit offers a safety bottle (1 or 2 liter), PFA tubing, and a mobile phase filter.

GL45 safety caps seal to the bottle with an O-ring\* (EPDM) which is resistant to commonly used HPLC solvents. If improved chemical resistance is required (e.g. for halogenated hydrocarbons), order an FEP-coated O-ring separately.

Bottle size:	No. of ports	1 Liter	2 Liters
		Product No.	Product No.
Safety cap kit, GL45	2	JR-S-11111	JR-S-11211
	3	JR-S-11112	JR-S-11212
	4	JR-S-11113	JR-S-11213
		Qty/pkg	Product No.
O-ring, FEP-coated, GL45		1	JR-S-40008



### INCLUDES

- 1 safety cap, GL45 thread, with fittings
- 1 safety-coated glass bottle (option: 1L or 2L)
- 3 m, PFA tubing
- 1 No-Met mobile phase filter, 5 µm
- Additional nuts and ferrules for the bottom side connections

\* GL45 caps produced beginning in 2012 will seal with an O-ring behind a PTFE lip.

### Spare parts

### for safety cap kit

	Tubing OD	Qty/pkg	Product No.
1/4-28 flangeless nuts, PPS	1/16" OD	10	JR-55070-10
	1/8" OD	10	JR-55071-10
Inverted ferrules, ETFE	1/16" OD	10	JR-041-10
	1/8" OD	10	JR-051-10
Safety-coated bottle with solid cap			
1 Liter		1	JR-9000-0005
2 Liters		1	JR-9000-0010
Plug, 1/4-28, PEEK		1	JR-409
Barbed hose adapter for 1/8" ID tubing		1	JR-075013
Color-Sleeve adapters, 12 colors (2 each)		24	JR-5508X-24
No-Met biocompatible mobile phase filter, economy slip-on connection, 5 µm pore		1	JR-32174
PFA tubing, 1/8" x 1.59 mm x 3m		1	JR-T-4003-M3

### TECH TIP

PFA tubing provides less re-gassing of the mobile phase.

### RELATED PRODUCTS

Color-Sleeve adapters . . . . . 17



No-Met mobile phase filters . . . . . 52  
PFA tubing . . . . . 9  
Safety-coated bottles (for additional sizes) . 70  
Safety air inlet valves . . 59

## Safety caps

## INCLUDES

- 1 cap, GL45 thread, with stopcocks
- 1 O-ring, EPDM
- 1 set: 1/16" nut/ferrule/color-sleeve adapter
- 1, 2, or 3 sets: 1/8" nut/ferrule/color-sleeve adapter
- Additional nuts and ferrules for the bottom side connections



- Fit all bottles or canisters with GL45 thread
- Built-in stopcocks for easy interruption of solvent lines
- 2, 3, or 4 ports with 1/4-28 threads for 1/8" or 1/16" tubing

GL45 safety caps seal to the bottle with an O-ring\* (EPDM) which is resistant to commonly used HPLC solvents. If improved chemical resistance is required (e.g. for halogenated hydrocarbons), order an FEP-coated O-ring separately.

	No. of ports	Product No.
Safety cap, GL45	2	JR-S-11011
	3	JR-S-11012
	4	JR-S-11013
O-ring, FEP-coated, GL45		JR-S-40008

## GL45 threads

## SPECS

- Connections: 1/4-28
- Materials
  - Closure: PP (Polypropylene)
  - Insert: PTFE
  - Stopcock: CTFE
  - O-ring: EPDM (alternative: FEP-coated)
  - Nuts: PPS
  - Ferrules: ETFE
  - Sleeves: PP

\* GL45 caps produced beginning in 2012 will seal with an O-ring behind a PTFE lip.

## Safety caps

## INCLUDES

- 1 cap, DIN50 thread, with stopcocks
- 1 O-ring,
- 1 set: 1/16" nut/ ferrule/ color-sleeve adapter
- 1, 2, or 3 sets: 1/8" nut/ferrule/color-sleeve adapter
- Additional nuts and ferrules for the bottom side connections



- Fit all bottles or canisters with DIN50 thread
- Built-in stopcocks for easy interruption of solvent lines
- 2, 3, or 4 ports with 1/4-28 threads for 1/8" or 1/16" tubing

DIN50 safety caps seal to the bottle/canister with an O-ring behind a PTFE lip, so the O-ring is therefore not a wetted part.

	No. of ports	Product No.
Safety cap, DIN50	2	JR-S-15011
	3	JR-S-15012
	4	JR-S-15013

## DIN50 threads

## SPECS

- Connections: 1/4-28
- Materials
  - Closure: PP (Polypropylene)
  - Insert: PTFE
  - Stopcock: CTFE
  - O-ring: EPDM/NBR
  - Nuts: PPS
  - Ferrules: ETFE
  - Sleeves: PP

## NOTE

One port must be reserved for the exhaust filter or the air inlet valve.

## MORE INFORMATION

Chemical resistance 96-98

## RELATED PRODUCTS

Air inlet valves..... 59  
Exhaust filters ..... 60-61

## Spare parts

## for safety caps

	Tubing OD	Qty/pkg	Product No.
1/4-28 flangeless nuts, PPS	1/16"	10	JR-55070-10
	1/8"	10	JR-55071-10
	1/16"	10	JR-041-10
Inverted ferrules, ETFE	1/8"	10	JR-051-10
Barbed hose adapter for 1/8" ID tubing		1	JR-075013
Color-Sleeve adapters, 12 colors (2 each)		24	JR-5508X-24

VICI caps seal to the bottle with an O-ring (EPDM) which is resistant to commonly used HPLC solvents. If improved chemical resistance is required (e.g. for halogenated hydrocarbons), order an FEP-coated O-ring separately. Unused ports should be closed with plugs, available in PEEK or polypropylene. (See "Spare Parts," below)

- Prevent evaporation of volatile compounds
- 2, 3, or 4 ports with 1/4-28 threads for 1/8" or 1/16" tubing
- 6 or 7 ports available in some models
- Optional bottom-sealing barbed adapter to connect 1/8" (3.2 mm) ID tubing (shown at right – order separately from "Spare Parts")



## VICI cap kit

GL45 threads

### SPECS

- Connections: 1/4-28
- Materials
  - Closure: PP (polypropylene)
  - Insert: PTFE
  - O-ring: EPDM (alternative: FEP-coated)
  - Nuts: PPS
  - Ferrules: ETFE
  - Sleeves: PP
  - Bottle: Borosilicate 3.3, safety coated
  - Filter: Polyethylene

In addition to a GL45 VICI cap with fittings, this kit includes a safety bottle (1 or 2 liter), PFA tubing, and a mobile phase filter.

Bottle size:	No. of ports	1 Liter	2 Liters
		Product No.	Product No.
VICI cap kit, GL45	2	JR-S-11101	JR-S-11201
	3	JR-S-11102	JR-S-11202
	4	JR-S-11103	JR-S-11203
O-ring, FEP-coated, GL45		Qty/pkg	Product No.
		1	JR-S-40008



### INCLUDES

- 1 VICI cap, GL45 thread, with fittings
- 1 safety-coated glass bottle (option: 1L or 2L)
- 3 m, PFA tubing
- 1 No-Met mobile phase filter, 5 µm

## VICI caps

GL45 threads

### SPECS

- Connections: 1/4-28
- Materials
  - Closure: PP (Polypropylene)
  - Insert: PTFE
  - Stopcock: CTFE
  - O-ring: EPDM (alternative: FEP-coated)
  - Nuts: PPS
  - Ferrules: ETFE
  - Sleeves: PP

- Fit all bottles or canisters with GL45 thread
- 2, 3, 4, 6, and 7 ports with 1/4-28 threads for 1/8" or 1/16" tubing

GL45 caps seal to the bottle with an O-ring\* (EPDM) which is resistant to commonly used HPLC solvents. If improved chemical resistance is required (e.g. for halogenated hydrocarbons), order an FEP-coated O-ring separately.

	No. of ports	Product No.
VICI cap, GL45	2	JR-S-11001
	3	JR-S-11002
	4	JR-S-11003
	6	JR-S-11005
	7	JR-S-11006
O-ring, FEP-coated, GL45		JR-S-40008



### INCLUDES

- 1 cap, GL45 thread
- 1 O-ring, EPDM
- 1 set: 1/16" nut/ferrule/ color-sleeve adapter
- 1, 2, 3, 5, or 6 sets: 1/8" nut/ferrule/ color-sleeve adapter



PTFE insert for 6 port cap



PTFE insert for 7 port cap

\* GL45 caps produced beginning in 2012 will seal with an O-ring behind a PTFE lip.

## Spare parts

for VICI caps

	Tubing OD	Qty/pkg	Product No.
1/4-28 flangeless nuts, PPS	1/16"	10	JR-55070-10
	1/8"	10	JR-55071-10
	1/8"	10	JR-051-10
Inverted ferrules, ETFE	1/16"	10	JR-041-10
	1/8"	10	JR-051-10
Plug, 1/4-28	PEEK	1	JR-409
	Polypropylene	1	JR-201999
Barbed hose adapter for 1/8" ID tubing		1	JR-075013
Color-Sleeve adapters, 12 colors (2 each)		24	JR-5508X-24

### Additional spare parts for VICI cap kits

	Qty/pkg	Product No.
Safety-coated bottle with solid cap		
	1 Liter	JR-9000-0005
No-Met biocompatible mobile phase filter, economy slip-on connection, 5 µm pore		
	2 Liters	JR-9000-0010
PFA tubing, 1/8" x 1.59 mm x 3m	1	JR-T-4003-M3

## VICI caps

NEW

S40 threads

## INCLUDES

1 cap, S40 thread  
 1 set: 1/16" nut/ ferrule/  
 color-sleeve adapter  
 2 sets:  
 1/8" nut/ferrule/  
 color-sleeve adapter



- Fit all bottles with S40 thread
- 3 ports with 1/4-28 threads for 1/8" or 1/16" tubing

The latest addition to the Jour Safety product line, the S40 VICI cap seals to the bottle with an O-ring (EPDM) which is resistant to commonly used HPLC solvents.

	No. of ports	Product No.
VICI cap, S40	3	JR-S-14002

## SPECS

- Connections:  
1/4-28
- Materials  
 Closure: PP\*  
 (Polypropylene)  
 Insert: PTFE  
 O-ring: EPDM  
 Nuts: PPS\*  
 Ferrules: ETFE  
 Sleeves: PP\*  
 \*not a wetted part

## VICI caps

DIN50 threads

## INCLUDES

1 cap, DIN50 thread  
 1 O-ring,  
 1 set: 1/16" nut/ ferrule/  
 color-sleeve adapter  
 1, 2, or 3 sets:  
 1/8" nut/ferrule/ color-  
 sleeve adapter



- Fit all bottles or canisters with DIN50 thread
- 2, 3, or 4 ports with 1/4-28 threads for 1/8" or 1/16" tubing

DIN50 safety caps seal to the bottle/canister with an O-ring behind a PTFE lip, so the O-ring is therefore not a wetted part.

	No. of ports	Product No.
VICI cap, DIN50	2	JR-S-15001
	3	JR-S-15002
	4	JR-S-15003



## SPECS

- Connections:  
1/4-28
- Materials  
 Closure: PP  
 (Polypropylene)  
 Insert: PTFE  
 Stopcock: CTFE  
 O-ring: EPDM/  
 NBR  
 Nuts: PPS  
 Ferrules: ETFE  
 Sleeves: PP

## VICI caps

GL38 threads

## INCLUDES

1 cap, GL37/430 thread  
 1 O-ring  
 1 set: 1/16" nut/ ferrule/  
 color-sleeve adapter  
 1, 2, or 3 sets:  
 1/8" nut/ferrule/  
 color-sleeve adapter  
 1 adapter ring



- Fit all bottles or canisters with GL38/430 thread  
(See NOTES, below)
- 2, 3, or 4 ports with 1/4-28 threads for 1/8" or 1/16" tubing

Color-coded PP sleeves fit on all ports on 2 and 3 port caps. With 4 port caps, the use of PP sleeves is limited to two sleeves per cap. (Four colors are included.)

	No. of ports	Product No.
VICI cap, GL38	2	JR-S-13001
	3	JR-S-13002
	4	JR-S-13003
O-ring, FEP-coated, GL38		JR-S-40009

## SPECS

- Connections:  
1/4-28
- Materials  
 Closure: PP  
 (Polypropylene)  
 Insert: PTFE  
 O-ring: EPDM  
 (alternative:  
 FEP-coated)  
 Adapter  
 ring: PA/PP  
 Nuts: PPS  
 Ferrules: ETFE  
 Sleeves: PP

## TECH TIP

PFA tubing provides less re-gassing of the mobile phase.

## MORE INFORMATION

Chemical resistance 96-98

## RELATED PRODUCTS

Color-Sleeve adapters . 17  
 No-Met mobile  
 phase filters . . . . . 52  
 PFA tubing . . . . . 9  
 Safety-coated bottles  
 (for additional sizes) . 70  
 Safety air inlet valves . . 59

## NOTES

GL38 thread is available in different variations. Please check the inside diameter of the neck:

- For 21-24 mm ID (e.g. Wheaton® bottles), use the GL38 cap as is, without an adapter.
- For 24-28 mm ID (e.g. Nalgene® bottles), use the adapter ring (JR-S-40006) which is included with this VICI cap.



Waste caps seal the bottle/canister with an O-ring behind a PTFE lip, so the O-ring is therefore not a wetted part. Unused ports should be closed with plugs, available in PEEK or polypropylene.  
(See "Spare Parts" on page 66.)

- Prevent evaporation of volatile compounds
- 3 or 4 ports with 1/4-28 threads for 1/8" or 1/16" tubing
- Additional port with 1/2-20 thread for 1/4" OD tubing or M10 x 1 thread for a barbed adapter
- Thread-sealing barbed adapters available in three sizes, to connect 4 mm, 6 mm, or 8 mm ID tubing (shown at right – order separately)



## Waste caps **NEW**

GL45 threads

### SPECS

- Connections: 1/4-28, 1/2-20
- Materials
  - Closure: PP (Polypropylene)
  - Insert: PTFE
  - O-ring: NBR
  - Nuts: PPS/POM
  - Ferrules: ETFE/CTFE
  - Plug: PTFE
  - Barbed adapter: PA/PP

- Fits all bottles or canisters with GL45 thread
- 4 ports with with one M10 x 1 port for barbed/angled adapters

This waste cap seals the bottle/cannister with an O-ring behind a PTFE lip, so the O-ring is therefore not a wetted part.



### INCLUDES

- 1 cap
- 1 O-ring
- 1 set: 1/16" nut/ ferrule/ color-sleeve adapter
- 2 sets: 1/8" nut/ferrule/ color-sleeve adapter

	No. of ports	Product No.
Waste cap, GL45	4	JR-S-11023

## Waste caps

GLS80 threads

### SPECS

- Connections: 1/4-28, 1/2-20
- Materials
  - Closure: PP (Polypropylene)
  - Insert: PTFE
  - O-ring: NBR
  - Nuts: PPS/POM
  - Ferrules: ETFE/CTFE
  - Plug: PTFE
  - Barbed adapter: PA/PP

- Fits all bottles or canisters with GLS80 thread
- 4 ports with 1/4-28 thread for 1/8" or 1/16" tubing
- Fifth port with choice of threads

Waste cap, GLS80

Fifth port:

1/2-20 thread for 1/4" OD tubing, with 1/4" plug JR-S-18004

M10 x 1 thread for barbed adapter JR-S-18014

(Order barbed adapter from "Spare Parts" on page 66)

Product No.



### INCLUDES

- 1 cap
- 1 O-ring
- 1 set: 1/16" nut/ ferrule/ color-sleeve adapter
- 1 set: 1/4" nut/ferrule
- 3 sets: 1/8" nut/ferrule/ color-sleeve adapter
- 1 plug, 1/4" (1/2-20), where noted in product number chart

### NOTES

One port must be reserved for the exhaust filter or the air inlet valve.

### MORE INFORMATION

Spare parts for  
waste caps ..... 68  
Chemical resistance 96-98

### RELATED PRODUCTS

Waste caps  
B53 threads ..... 68  
B83 threads ..... 68



# Waste caps

## DIN60 threads

### INCLUDES

- 1 cap
- 1 O-ring
- 1 set: 1/16" nut/ ferrule/ color-sleeve adapter
- 1 set: 1/4" nut/ferrule
- 2 sets:
  - 1/8" nut/ferrule/ color-sleeve adapter
- 1 plug, 1/4" (1/2-20), where noted in product number chart



- Fits all bottles or canisters with DIN60 thread
- 3 ports with 1/4-28 thread for 1/8" or 1/16" tubing
- Fourth port with choice of threads

Product No.

Waste cap, DIN60

Fourth port:

1/2-20 thread for 1/4" OD tubing, with 1/4" plug JR-S-16003

M10 x 1 thread for barbed adapter JR-S-16013

(Order barbed adapter from "Spare Parts" on page 64)

### SPECS

- Connections: 1/4-28, 1/2-20
- Materials
  - Closure: PP (Polypropylene)
  - Insert: PTFE
  - O-ring: NBR
  - Nuts: PPS/POM
  - Ferrules: ETFE/CTFE
  - Plug: PTFE
  - Barbed adapter: PA/PP

# Waste caps

## S51 threads

### INCLUDES

- 1 cap
- 1 O-ring
- 1 set: 1/16" nut/ ferrule/ color-sleeve adapter
- 1 set: 1/4" nut/ferrule
- 2 sets:
  - 1/8" nut/ferrule/ color-sleeve adapter
- 1 plug, 1/4" (1/2-20), where noted in product number chart



- Fits all bottles or canisters with S51 thread
- 3 ports with 1/4-28 thread for 1/8" or 1/16" tubing
- Fourth port with choice of threads

Product No.

Waste cap, S51

Fourth port:

1/2-20 thread for 1/4" OD tubing, with 1/4" plug JR-S-15103

M10 x 1 thread for barbed adapter JR-S-15113

(Order barbed adapter from "Spare Parts" on page 66)

### SPECS

- Connections: 1/4-28, 1/2-20
- Materials
  - Closure: PP (Polypropylene)
  - Insert: PTFE
  - O-ring: NBR
  - Nuts: PPS/POM
  - Ferrules: ETFE/CTFE
  - Plug: PTFE
  - Barbed adapter: PA/PP

# Waste caps

## S61 threads

### INCLUDES

- 1 cap
- 1 O-ring
- 1 set: 1/16" nut/ ferrule/ color-sleeve adapter
- 1 set: 1/4" nut/ferrule
- 2 sets:
  - 1/8" nut/ferrule/ color-sleeve adapter
- 1 plug, 1/4" (1/2-20), where noted in product number chart



- Fits all bottles or canisters with S61 thread
- 3 ports with 1/4-28 thread for 1/8" or 1/16" tubing
- Fourth port with choice of threads

Product No.

Waste cap, S61

Fourth port:

1/2-20 thread for 1/4" OD tubing, with 1/4" plug JR-S-16103

M10 x 1 thread for barbed adapter JR-S-16113

(Order barbed adapter from "Spare Parts" on page 66)

### SPECS

- Connections: 1/4-28, 1/2-20
- Materials
  - Closure: PP (Polypropylene)
  - Insert: PTFE
  - O-ring: NBR
  - Nuts: PPS/POM
  - Ferrules: ETFE/CTFE
  - Plug: PTFE
  - Barbed adapter: PA/PP

### RELATED PRODUCTS

- Color-Sleeve adapters..... 17
- Exhaust filters .....60-61
- Air inlet valves..... 59

## Waste caps

B53 threads

## SPECS

- Connections:  
1/4-28, 1/2-20
- Materials  
Closure: PP  
(Polypropylene)  
Insert: PTFE  
O-ring: NBR  
Nuts: PPS/POM  
Ferrules: ETFE/  
CTFE  
Plug: PTFE  
Barbed  
adapter: PA/PP

- Fits all bottles or canisters with B53 thread
- 3 ports with 1/4-28 thread for 1/8" or 1/16" tubing
- Fourth port with choice of threads

Waste cap, B53

Fourth port:

1/2-20 thread for 1/4" OD tubing, with 1/4" plug

M10 x 1 thread for barbed adapter

(Order barbed adapter from "Spare Parts" below)

Product No.

JR-S-15303

JR-S-15313



## INCLUDES

- 1 cap
- 1 O-ring
- 1 set: 1/16" nut/ ferrule/  
color-sleeve adapter
- 1 set: 1/4" nut/ferrule
- 2 sets:  
1/8" nut/ferrule/  
color-sleeve adapter
- 1 plug, 1/4" (1/2-20),  
where noted in product  
number chart

## Waste caps

B83 threads

## SPECS

- Connections:  
1/4-28, 1/2-20
- Materials  
Closure: PP  
(Polypropylene)  
Insert: PTFE  
O-ring: NBR  
Nuts: PPS/POM  
Ferrules: ETFE/  
CTFE  
Plug: PTFE  
Barbed  
adapter: PA/PP

- Fits all bottles or canisters with B83 thread
- 4 ports with 1/4-28 thread for 1/8" or 1/16" tubing
- Fifth port with choice of threads

Waste cap, B83

Fifth port:

1/2-20 thread for 1/4" OD tubing, with 1/4" plug

M10 x 1 thread for barbed adapter

(Order barbed adapter from "Spare Parts" below)

Product No.

JR-S-18304

JR-S-18314



## INCLUDES

- 1 cap
- 1 O-ring
- 1 set: 1/16" nut/ ferrule/  
color-sleeve adapter
- 1 set: 1/4" nut/ferrule
- 3 sets:  
1/8" nut/ferrule/  
color-sleeve adapter
- 1 plug, 1/4" (1/2-20),  
where noted in product  
number chart

## Spare parts

for waste caps

	Material	Tubing OD	Qty/pkg	Product No.
1/4-28 flangeless nuts	PPS	1/16"	10	JR-55070-10
		1/8"	10	JR-55071-10
Inverted ferrules	ETFE	1/16"	10	JR-041-10
		1/8"	10	JR-051-10
1/2-20 nut	POM	1/4"	1	JR-CFL-4D
Collapsible ferrule	CTFE	1/4"	1	JR-CFL-CB4KF-S
Adapter				
1/2-20 to 1/4-28	POM	<b>NEW</b>	1	JR-S-40200
Barbed adapter, M10 x 1				
For 4 mm ID tubing	PA		1	JR-S-40104
For 6 mm ID tubing	PA		1	JR-S-40106
For 8 mm ID tubing	PA		1	JR-S-40108
Barbed hose adapter				
For 1/8" ID tubing	PA/PP		1	JR-075013
Plug, 1/4-28	PEEK		1	JR-409
	Polypropylene		1	JR-201999
Plug, 1/2-20	PTFE		1	JR-S-40007
Color-Sleeve adapters, 12 colors (2 each)			24	JR-5508X-24

## NOTES

One port must be reserved for the exhaust filter or the air inlet valve.

## MORE INFORMATION

Chemical resistance 96-98

## RELATED PRODUCTS

Waste caps

DIN60 threads ..... 67

GL45 threads ..... 66

GLS80 threads ..... 66

S51 threads ..... 67

S61 threads ..... 67

Color-Sleeve

adapters ..... 17

Exhaust filters ..... 60-61

Air inlet valves ..... 59

**Eco cap kit***GL45 threads***INCLUDES**

- 1 cap, GL45 thread, with three ports
- 1 safety-coated glass bottle (option: 1 L or 2 L)
- 1 male Luer plug, polypropylene

In addition to an Eco cap with polypropylene male Luer plug, this kit offers a safety bottle (1 or 2L).

The Eco cap consists of a three-hole ETFE insert with a polyethylene collar, compatible with all standard GL45 thread bottles. The Eco cap has no threading for fittings on the inside of the tubing ports – the tubing is just inserted.

**SPECS**

- Materials
  - Screw collar: PP (polypropylene)
  - Pouring ring: PP
  - Insert: ETFE
  - O-rings: EPDM
  - Luer plug: PP

*Product No.*

Eco cap kit

- |                |              |
|----------------|--------------|
| 1 Liter bottle | JR-9000-0013 |
| 2 Liter bottle | JR-9000-0014 |

**Eco cap***GL45 threads***INCLUDES**

- 1 cap, GL45 thread, with three ports
- 1 male Luer plug, polypropylene

- Competitively priced alternative to sparge and deliver mobile phases
- Also available as bottle kit

The Eco cap consists of a three-hole ETFE insert with a polyethylene collar, compatible with all standard GL45 thread bottles. The Eco cap has no threading for fittings on the inside of the tubing ports – the tubing is just inserted.

**SPECS**

- Materials
  - Screw collar: PP (polypropylene)
  - Pouring ring: PP
  - Insert: ETFE
  - O-rings: EPDM
  - Luer plug: PP

*Product No.*

- |         |              |
|---------|--------------|
| Eco cap | JR-9000-0007 |
|---------|--------------|

**Spare parts***for Eco cap**Product No.*

- |   |              |
|---|--------------|
| Male Luer plug, polypropylene<br>10/pkg | JR-070350-10 |
|---|--------------|

**NOTE**

These caps are not usable for building up a helium atmosphere within in the solvent bottle, but only for continuous helium sparging.

**TECH TIP**

Once the mobile phase is degassed, we recommend our impermeable tubing between bottle cap and pump to prevent regassing of helium degassed solvents.

**RELATED PRODUCTS**

Impermeable tubing . . 13  
Polymeric tubing . . . 9-13  
Waste caps . . . . . 66-68

## SAFETY PRODUCTS

## Safety-coated glass bottles

## SPECS

- Pressure rating  
0.8 bar (11.6 psi)
- Materials  
Bottle:  
Borosilicate  
3.3 glass  
Cap: PP  
(Polypropylene)  
Pouring ring: PP

- Special design for helium sparging of LC mobile phases
- Safety coating to minimize risk of injury
- Available in 0.25, 0.5, 1 and 2 Liter sizes

These glass bottles with GL45 screw threading are made of chemical resistant Borosilicate 3.3 and are coated on the outside with a layer of sturdy polymer, decreasing the chances of breakage if the bottle is dropped. Even if the bottle breaks, glass fragments will be kept within the polymer coating. Solvents will be kept within the coating for a short while. The coating also protects solvents against UV light up to 385 nm.

These safety-coated glass bottles cannot be autoclaved.

Bottles include a solid cap.

	Bottle size	Product No.
Safety-coated bottle	0.25 L	JR-9000-0011
	0.5 L	JR-9000-0012
	1 L	JR-9000-0005
	2 L	JR-9000-0010



## SolviFlex

## SPECS

- 1/18" OD tubing
- Materials  
Tubing: PFA  
Nuts: PPS  
Ferrules: ETFE  
Filter: PE  
(Polyethylene)

- Be "mobile" with mobile phase bottles
- Ideal for with VICI safety caps or VICI caps
- Customized version on request

The SolviFlex is made of flexible telephone cable-like PFA tubing which expands and contracts, providing flexibility to move mobile phase bottles up to 1.5 m radius on the bench. Approximately 3 m of tubing is coiled with straight stretches of tubing at each end.

Includes:

- ~3 m coiled 1/8" PFA tubing
- 2 sets: 1/4-28 nut/ferrules for 1/8" OD tubing
- 1 No-Met mobile phase filter with economy slip-on connection (JR-32174)

	Product No.
SolviFlex assembly	JR-6120-50030



**A typical use of SolviFlex**  
Shown installed on a VICI cap and bottle, along with an exhaust filter



## NOTES

- The bottles are not suitable to be stressed with pressures above 0.8 bar (11.6 psi), but only for continuous helium sparging.
- The safety coating on a bottle is not intended to replace safe working conditions and precautions in laboratories.

## RELATED PRODUCTS

No-Met biocompatible mobile phase filter . . .	52
Safety caps . . . . .	63
Safety cap kit . . . . .	62
Safety exhaust filter 60-61	
VICI caps . . . . .	64-65
VICI cap kit . . . . .	64

**Carbon dioxide filter****NEW****CO<sub>2</sub>**

In ion chromatography, a basic solution is often used as mobile phase. The buffer capacity of such solutions is small, so the pH of the solution can be affected by the carbon dioxide in the air.

Using a bottle in combination with a VICI cap and a carbon dioxide filter prevents any change in the pH of the mobile phase, thus increasing the stability of the analysis. In this application, the carbon dioxide filter also functions as an air inlet valve.

<i>Description</i>	<i>Product No.</i>
Carbon dioxide filter	JR-S-20100

**Humidity filter****NEW****H<sub>2</sub>O**

In applications requiring a non-aqueous organic solvent, the moisture in the ambient air can influence the results. The use of a humidity filter in combination with a VICI cap solves this problem.

<i>Description</i>	<i>Product No.</i>
Humidity filter	JR-S-20200

**RELATED PRODUCTS**

Safety exhaust filters. . . 60

# Frits

## REPRODUCIBILITY AND PRECISION

VICI Jour unidense frits, with uniform flow across the entire frit, offer unexcelled reproducibility and precision. Proprietary manufacturing processes ensure uniformity and precise, repeatable dimensional control.

PEEK-encased frits surround the frit with a ring of PEEK, reducing the actual porous surface to the ID of an HPLC column or filter.

## Frits index

### Unidense frits

Stainless.....	page 73
Titanium.....	73
Polyethylene.....	72

### PEEK-encased frits

Stainless.....	74
Titanium.....	75
Polyethylene.....	75

## Polyethylene frits – unidense type

### SPECS

- Material: UHMWPE
- Tolerances:  $\pm 0.08$  mm (.003")

- Biocompatible
- Metal free

Our completely metal-free polymer frits are made of UHMWPE (Ultra High Molecular Weight PolyEthylene), and offer a competitively priced alternative to titanium.

For single use only. Sold in packages of 5.

OD		Thickness		Porosity	Frit volume	Product No.
inches	mm	inches	mm	$\mu$ m	$\mu$ L	
.255"	6.50	..079"	2.00	10	26.5	JR-F-255-079-10PE-5



### OPTIONS AND MORE INFORMATION

See facing page.

### RELATED PRODUCTS

PEEK-encased polyethylene frits . . . 75





### Stainless steel frits – unidense type

For both standard frits and those produced to custom specifications, these frits offer unexcelled reproducibility and precision. Proprietary manufacturing processes ensure uniformity and precise, repeatable dimensional control.

Stainless steel frits are manufactured of SS316L, which has extra-low carbon content. The sintering process used results in a strong bonding of the sinter metal particles. Depending on the grade of powder used, the porosity of the frit can be precisely tuned.

Sold in packages of 5.

#### SPECS

- Material:  
Stainless steel  
Type 316L  
sintered
- Tolerances:  
See chart below

OD		Thickness		Porosity	Frit volume	Product No.
inches	mm	inches	mm	μm	μL	
1/16"	1.59	.030"	0.76	0.5	0.39	JR-.5FR1-5
				1	0.45	JR-1FR1-5
				2	0.53	JR-2FR1-5
				10	0.60	JR-10FR1-5
1/8"	3.18	.040"	1.02	0.5	2.10	JR-.5FR2-5
				1	2.43	JR-1FR2-5
				2	2.83	JR-2FR2-5
				10	3.24	JR-10FR2-5
1/4"	6.35	.040"	1.02	0.5	8.39	JR-.5FR4-5
				2	11.30	JR-2FR4-5
				10	12.91	JR-10FR4-5
3/8"	9.53	.040"	1.02	0.5	18.91	JR-.5FR6-5
				2	25.45	JR-2FR6-5
				10	29.09	JR-10FR6-5
1/2"	12.70	.040"	1.02	0.5	33.58	JR-.5FR8-5
				2	45.2	JR-2FR8-5
				10	51.66	JR-10FR8-5
1"	25.40	.060"	1.52	2	269.43	JR-2FR1K-5
				10	307.92	JR-10FR1K-5

#### OPTIONS

- Standard frits do not have chamfers. Chamfered frits are available on request.
- Other dimensions and porosities are available on request. Minimum order quantity: 100 pieces. Please contact your local distributor or VICI directly.

#### TOLERANCES

OD (mm)	Tolerance
<5	±0.05 mm (.002")
5-12	±0.08 mm (.003")
12.5-25.4	±0.20 mm (.008")

#### NOTES

- Stated porosities are only nominal and do not reflect the maximum pore size of a frit.
- See frit volume chart for help choosing the right frit and porosity for your application. . . . . 95
- Frit volumes are theoretical and are calculated by multiplying overall frit volume times the porosity proportion.

#### RELATED PRODUCTS

PEEK-encased frits	
Stainless . . . . .	74
Titanium . . . . .	75

### Titanium frits – unidense type

- Biocompatible
- Recommended for protein analysis
- High corrosion resistance
- Various pore sizes and diameters

Titanium is preferred over stainless steel for applications related to the analysis of sensitive substances in general and large biomolecules in particular. Evidently proteins have a tendency to adsorb on the stainless steel (iron) frit and consequently reduce recovery or may even decompose. Titanium also offers higher corrosion resistance and can be used with solvents containing halides.

Sold in packages of 5.

OD		Thickness		Porosity	Frit volume	Product No.
inches	mm	inches	mm	μm	μL	
1/16"	1.59	.030"	0.76	0.5	0.82	JR-.5FR1TI-5
				2	0.71	JR-2FR1TI-5
1/8"	3.18	.040"	1.02	2	2.83	JR-2FR2TI-5
1/4"	6.35	.040"	1.02	2	11.30	JR-2FR4TI-5
3/8"	9.53	.040"	1.02	2	25.45	JR-2FR6TI-5
1/2"	12.70	.040"	1.02	2	45.2	JR-2FR8TI-5
				10	51.66	JR-10FR8TI-5
1"	25.40	.060"	1.52	2	180.80	JR-2FR1KTI-5

#### SPECS

- Material:  
Titanium,  
sintered
- Tolerances:  
See chart at left



PEEK-encased frits surround the frit with a ring of PEEK, reducing the actual porous surface to the ID of an HPLC column or filter.

The compressed PEEK ring forms a seal between the bottom of the fitting and the column end, preventing the mobile phase and the sample from entering the previously accessible poorly swept areas.

### Stainless steel frits – PEEK-encased

#### SPECS

- Material:  
Frit Stainless  
SS316L  
Ring PEEK
- Tolerances:  
±0.05 mm  
(.002")

- For HPLC columns or in-line filters
- Eliminate poorly swept volume and improves peak symmetry
- Various porosities and dimensions available

Our PEEK-encased frits with 0.118" or 1/4" OD ring are intended for HPLC columns with IDs of 3.2, 3.9, and 4.6 mm. With the frit surrounded by a PEEK ring, the actual porous surface is reduced to the column ID.

The compressed PEEK ring forms a seal between the bottom of the fitting and the column end, preventing the mobile phase and the sample from entering the previously accessible poorly swept areas.

Sold in packages of 5.



Ring OD		Frit OD		Thickness		Porosity	Frit volume	Product No.
inches	mm	inches	mm	inches	mm	μm	μL	
.118"	3.00	.038"	0.97	.028"	0.70	0.5	0.13	JR-1110-05P-5
						2	0.18	JR-1110-2P-5
						5	0.20	JR-1110-5P-5
		.077"	1.96	.028"	0.70	0.5	.055	JR-1111-05P-5
						2	0.74	JR-1111-2P-5
						5	0.80	JR-1111-5P-5
						10	0.87	JR-1111-10P-5
.25"	6.35	.078"	2.00	.062"	1.59	0.5	1.30	JR-1104-05P-5
						2	1.74	JR-1104-2P-5
		.125"	3.20	.062"	1.59	0.5	3.32	JR-1102-05P-5
						2	4.47	JR-1102-2P-5
		.153"	3.90	.062"	1.59	0.5	4.94	JR-1101-05P-5
						2	4.47	JR-1101-2P-5
		.181"	4.60	.062"	1.59	0.5	6.87	JR-1100-05P-5
						2	9.24	JR-1100-2P-5
		.125"	3.20	.062"	1.59	0.5	1.75	JR-1103-05P-5
						2	2.25	JR-1103-2P-5

#### OPTIONS

- Other dimensions and porosities are available on request. Minimum order quantity: 100 pieces. Please contact your local distributor or VICI directly.

#### NOTES

- Stated porosities are only nominal and do not reflect the maximum pore size of a frit.
- See frit volume chart for help choosing the right frit and porosity for your application. ....95
- Frit volumes are theoretical and are calculated by multiplying overall frit volume times the porosity proportion.

#### RELATED PRODUCTS

Unidense type frits	
Polyethylene	72
Stainless	73
Titanium	73



### Titanium frits – PEEK-encased

- For HPLC columns or in-line filters
- Biocompatible – recommended for protein analysis
- Eliminate poorly swept volume and improves peak symmetry
- Higher corrosion resistance
- Various porosities and dimensions available

#### SPECS

- Material:  
Frit Titanium  
Ring PEEK
- Tolerances:  
±0.05 mm  
(.002")

In terms of biocompatibility, the frit is the critical spot for contributing metal ions to the mobile phase and for interaction with protein samples.

VICI Jour PEEK-encased frits with 1/5" or 1/4" OD ring are intended for HPLC columns with IDs of 3.2 and 4.6 mm. With the frit surrounded by a PEEK ring, the actual porous surface is reduced to the column ID.

The compressed PEEK ring forms a seal between the bottom of the fitting and the column end, preventing the mobile phase and the sample from entering the previously accessible poorly swept areas.

Sold in packages of 5.

Ring OD		Frit OD		Thickness		Porosity	Frit volume	Product No.
inches	mm	inches	mm	inches	mm	μm	μL	
1/5"	5.08	.102"	2.60	.062"	1.59	5	3.20	JR-1128-5P-5
.25"	6.35	.125"	3.20	.062"	1.59	0.5	3.32	JR-1127-05P-5
						2	4.47	JR-1127-2P-5
						0.5	6.87	JR-1125-05P-5
		.181"	4.60	.062"	1.59	2	2.94	JR-1125-2P-5
						5	10.04	JR-1125-5P-5
						10	10.56	JR-1125-10P-5



### Polyethylene frits – PEEK-encased

- Biocompatible
- Metal free

Our polymer frits are made of UHMWPE (Ultra High Molecular Weight PolyEthylene). They are a competitively priced alternative to titanium and are completely metal free.

With the frit surrounded by a PEEK ring, the actual porous surface is reduced to the column ID. The compressed PEEK ring forms a seal between the bottom of the fitting and the column end, preventing the mobile phase and the sample from entering the previously accessible poorly swept areas.

Sold in packages of 5.

Ring OD		Frit OD		Thickness		Porosity	Frit volume	Product No.
inches	mm	inches	mm	inches	mm	μm	μL	
1/5"	5.08	.098"	2.50	.059"	1.50	10	2.94	JR-1150-10P-5
.25"	6.35	.181"	4.60	.059"	1.50	10	11.31	JR-1151-10P-5

#### SPECS

- Material:  
Frit UHMWPE  
Ring PEEK
- Tolerances:  
±0.05 mm  
(.002")

#### OPTIONS AND MORE INFORMATION

See facing page.

#### TECH TIPS

- PEEK-encased titanium frits can also be used as spare parts frits for the PEEK in-line filter kit . . .54
- PEEK-encased polyethylene frits can be used as spare parts for:  
PEEK in-line filter kit . . . . .54  
PEEK columns . . . . .77



# Columns

## PEEK COLUMNS ■ GUARD COLUMNS

### Columns index

#### PEEK columns

2.1 mm ID .....	page 77
4.0 mm ID .....	77

#### In-line biocompatible

##### guard columns .....78

This PEEK guard column system extends the life of valuable analytical columns by trapping contaminants from the sample and mobile phase. Only PEEK and titanium contact the mobile phase and sample.

#### Column couplers and plugs

Column coupler, one-piece fingertight .....	79
Column end plugs .....	79
Column end plugs, logotype .....	79

##### Packing kits .....78

VICI Jour column packing kits are ready-to-use tools to pack our PEEK columns and PEEK guard columns. Includes needed slurry reservoir, packing adapter, and packing seals.





### 2.1 mm ID PEEK columns

- Biocompatible
- High pressure applications
- Packing devices available

We have designed this column for narrow-bore applications, and for those when a limit amount of sample is available for analysis.

Because 2.1 mm ID columns operate at 1/5th of the flow rate of their 4.6 mm counterparts, the same sample mass injected on a 2.1 mm column will produce five times the detector signal. This greater mass sensitivity provides the same analysis as on a 4.6 mm column with only 20% of the sample.

All our columns are delivered without frits. Purchase frits separately.


#### SPECS

- Material: PEEK, black and natural
- Tolerances:  $\pm 0.05$  mm (.002")
- Pressure rating maximum: 210 bar (3000 psi)  
During packing, using a VICI Jour packing kit: 350 bar (5000 psi)

Length	Product No.
30 mm	JR-65001NF
50 mm	JR-65002NF
100 mm	JR-65003NF
150 mm	JR-65005NF
250 mm	JR-65006NF

#### Spare parts

Description	Qty/pkg	Product No.
Column end fitting	1	JR-062
PEEK-encased frit		
Titanium, 5 $\mu$ m	5	JR-1128-5P-5
Polyethylene, 10 $\mu$ m	5	JR-1150-10P-5

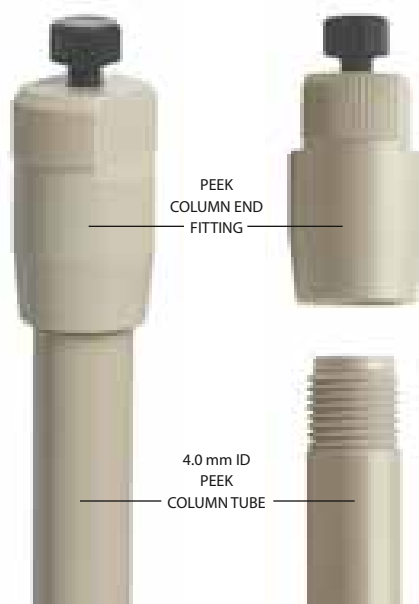

### 4.0 mm ID PEEK columns

- Biocompatible
- High pressure applications
- Packing devices available

Today's HPLC and IC applications are more demanding than ever. The limitations of traditional stainless steel pose significant problems for a growing number of important biotechnology and ion chromatography applications.

All our columns are delivered without frits. Purchase frits separately.

#### SPECS

- Material: PEEK, black and natural
- Tolerances:  $\pm 0.05$  mm (.002")
- Pressure rating maximum: 350 bar (5000 psi)  
During packing, using a VICI Jour packing kit: 560 bar (8000 psi)

Length	Product No.
50 mm	JR-68181NF
100 mm	JR-68182NF
125 mm	JR-68183NF
150 mm	JR-68184NF
250 mm	JR-68185NF
300 mm	JR-68186NF

#### Spare parts

Description	Qty/pkg	Product No.
Column end fitting	1	JR-66140
PEEK-encased frit		
Titanium, 5 $\mu$ m	5	JR-1125-2P-5
Polyethylene, 10 $\mu$ m	5	JR-1151-10P-5

#### RELATED PRODUCTS

Column packing kit.... 78  
PEEK-encased frits  
    Titanium ..... 75  
    Polyethylene frits .... 75

## In-line biocompatible guard column

## SPECS

- Material: PEEK, titanium
- Pressure rating: < 350 bar (< 5000 psi)


- Biocompatible
- Easy to fill
- Packing device for self-packing available

The VICI Jour in-line PEEK guard column system extends the life of valuable analytical columns by trapping contaminants from the sample and mobile phase, thereby preventing loss of retention.

Only PEEK and titanium contact the mobile phase and sample. This system is ideal for use with proteins or other biopolymers that may be contaminated or absorbed by steel surfaces. VICI Jour guard columns can be easily dry packed or slurry packed. Guard columns may be changed in minutes without the use of tools – all connections may be made by hand.

Guard column system includes:

- Guard column cartridge
- 2 Titanium frits
- Housing
- O-ring

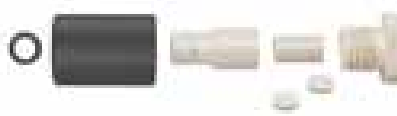

Guard column cartridge includes:

- Guard column cartridge
- 2 Titanium frits (2 µm or 5 µm)



Description		Product No.
Guard column system	2 µm	JR-68252TI2
	5 µm	JR-68252TI
Guard cartridge with 2 frits	2 µm	JR-68242TI2
	5 µm	JR-68242TI

## Spare parts

Description	Qty/pkg	Product No.
Replacement frit, titanium	2 µm	JR-F-197-062-2TINC
	5 µm	JR-F-197-062-5TI
Packing seal	1	JR-9

## Packing kits

## SPECS

- Material: PEEK
  - Pressure ratings during packing, using a VICI Jour packing kit:
- | Column ID | Packing pressure   |
|-----------|--------------------|
| 2.1 mm    | 350 bar (5000 psi) |
| 4.0 mm    | 560 bar (8000 psi) |
- See column pages for operating pressure ratings.

VICI Jour column packing kits are ready-to-use tools to pack our PEEK columns and PEEK guard columns. Every kit consists of the needed slurry reservoir, packing adapter, and packing seals.


Description	Product No.
<b>For 2.1 mm ID columns</b>	
Packing kit	JR-67-250-2A
Replacement PEEK packing adapter	JR-67140-1A
Replacement PEEK packing seal	JR-65011
<b>For 4.0 mm ID columns</b>	
Packing kit	JR-67-250-4
Replacement PEEK packing adapter	JR-67140-2B
Replacement PEEK packing seal	JR-65158
<b>For in-line biocompatible guard columns</b>	
Packing kit	JR-27305
Replacement CTFE packing seal	JR-9

## TECH TIP

Our packing kit is a "ready-to-use" tool to pack our PEEK guard cartridges. Every kit consists of the needed slurry reservoir, packing adapter, and packing seals.

## RELATED PRODUCTS

PEEK columns	
2.1 mm ID	77
4.0 mm ID	77

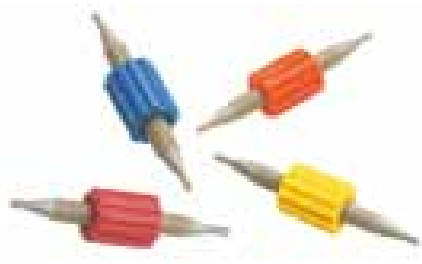

### One-piece fingertight column couplers

- Self-adjusting – zero dead volume for all fitting types
- Four bore sizes
- Biocompatible

Choose from a variety of coupler IDs, indicated by the color of the sleeve. A unique feature of this column coupler is that it adapts automatically to fit all pilot lengths – Valco, Waters, Upchurch, Rheodyne, etc. Since the tubing bottoms out in any fitting detail, added void volume is minimal.

Sold individually.

Color	ID	Product No.
Red	0.13 mm	JR-26501
Yellow	0.17 mm	JR-26502
Blue	0.25 mm	JR-26503
Orange	0.50 mm	JR-26504



**Coupler shown installed between two columns**  
(columns not included)

#### SPECS

- Material:
  - Tubing PEEK
  - Fittings PEEK
  - Sleeve PP (Polypropylene)
- Threads:
  - 10-32 male to 10-32 male
- Pressure rating:
  - <350 bar (5000 psi)


### Column end plugs – logotype

- Customizable with your logo and color (OEM version)
- Same features as column end plugs below

Ensure safe closure of valuable columns with your own customized column end plugs. Easy handling is guaranteed by the 5/16" wrench flat. Fit all columns with 10-32 thread.

Sold in packages of 10.

Color	Product No.
White	JR-4081-10
Black	JR-4082-10
Red	JR-4083-10
Green	JR-4084-10
Blue	JR-4085-10

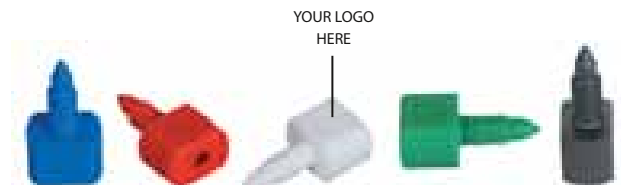

#### SPECS

- Material: Nylon
- Threads: 10-32
- Tolerances:
  - ±0.05 mm (.002")


### Column end plugs

- Seal analytical columns
- One size fits all
- Fingertight, with wrench flats as well

Fingertight column end plugs fit all columns with 10-32 threads. Easy handling is guaranteed by an additional 5/16" wrench flat, which allows the plug to be loosened by a wrench if it gets stuck. The long design and threading down to the ferrule ensure safe closing of valuable columns.

Sold in packages of 10.

Color	Product No.
White	JR-4071-10
Black	JR-4072-10
Red	JR-4073-10
Green	JR-4074-10
Blue	JR-4075-10

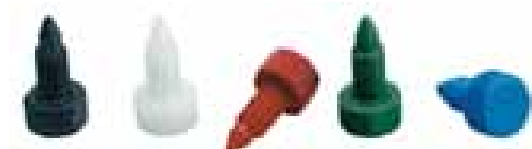

#### SPECS

- Material: Nylon
- Threads: 10-32
- Tolerances:
  - ±0.05 mm (.002")

#### RELATED PRODUCTS

Color-coded  
PEEK tubing ..... 6-8

# Sample loops

VALCO ■ CHEMINERT® ■ RHEODYNE®

Sample loops are accurately sized for each valve type. However, with small volume loops, the tolerance on the ID of the tubing ( $\pm 0.001$ " ) can have a significant effect on the volume. Therefore loop volume and loop appearance may differ from batch to batch.

## Sample loops index

### PEEK loops

PEEK loops for:

Cheminert.....	page 81
Rheodyne .....	81
Valco .....	81

### Stainless steel loops

Stainless loops for:

Cheminert.....	82
Rheodyne .....	83
Valco .....	82

### SPECS

- Material:
  - Loops PEEK
  - Nuts PEEK
  - Ferrules PEEK
  - Caps PE
  - (Polyethylene)
- Tolerances:
  - $\pm 0.025$  mm
  - (.001")

## PEEK loops

- Wide range of volumes available
- Biocompatible

VICI Jour PEEK sample loops are made from striped color-coded PEEK tubing and PEEK fittings. PEEK is a mechanically strong, chemically inert polymer ideal for HPLC applications. Since our PEEK sample loops contain no metals or other extractable compounds, no foreign materials will be introduced into the system. Softener-free polyethylene caps are used to avoid contamination of the tubing.

PEEK sample loops are particularly well suited for separations involving proteins, peptides, nucleic acids or other samples of biological origin as well as in IC systems where metal traces may interact with ionic samples or release transition metals that may deactivate ion exchange columns.

PEEK hex-head nuts and PEEK ferrules provide reliable connections, and avoid twisting of the tube.

### DISCLAIMER

The maximum holding pressure for any type of connections involving tubing and a ferrule varies considerably with the tubing material the ferrule material, the clearance between tubing OD and ferrule ID, and the shape of the fitting detail.

### MORE INFORMATION

ValvTool..... 88



**PEEK loops***for Cheminert valves*

Fit all Cheminert C1, C1CF, C2, C2V, and C3 type valves.

Sample size	Pressure * bar	psi	Product No.	Sample size	Pressure * bar	psi	Product No.
5 µl	400	5800	JR-SLC-85-5	250 µl	240	3500	JR-SLC-85-250
10 µl	400	5800	JR-SLC-85-10	300 µl	240	3500	JR-SLC-85-300
20 µl	400	5800	JR-SLC-85-20	500 µl	240	3500	JR-SLC-85-500
50 µl	400	5800	JR-SLC-85-50	1 ml	240	3500	JR-SLC-85-1000
100 µl	400	5800	JR-SLC-85-100				
200 µl	240	3500	JR-SLC-85-200				

**PEEK loops***for Rheodyne valves*

Fit all Rheodyne metal and metal-free valves.

Sample size	Pressure * bar	psi	Product No.	Sample size	Pressure * bar	psi	Product No.
1 µl	435	6300	JR-SLR-2190	50 µl	350	4500	JR-SLR-2202
2 µl	420	6100	JR-SLR-2191	100 µl	350	4500	JR-SLR-2204
3 µl	420	6100	JR-SLR-2192	200 µl	240	3500	JR-SLR-2206
4 µl	420	6100	JR-SLR-2193	250 µl	240	3500	JR-SLR-2207
5 µl	400	5800	JR-SLR-2194	500 µl	240	3500	JR-SLR-2210
10 µl	400	5800	JR-SLR-2196	1 ml	240	3500	JR-SLR-2222
20 µl	386	5600	JR-SLR-2195				
30 µl	386	5600	JR-SLR-2200				
40 µl	386	5600	JR-SLR-2201				

\* Recommended maximum permanent working pressure

**TECH TIPS**

- See general information about PEEK loops on facing page.
- The tighter ID tolerance of the tubing minimizes the tolerance of the loop volume.
- The ValvTool tightens nuts with tube or loop installed, providing easy access to many hard-to-reach areas. .... 88

**OPTIONS**

Other sample sizes are available on request. Please contact your distributor or VICI Jour directly.

**MORE INFORMATION**

Chemical resistance 96-98  
About PEEK ..... 99

**PEEK loops***for Valco valves*

Fit all 1/16" Valco W and UW type valves.

Sample size	Pressure * bar	psi	Product No.
5 µl	400	5800	JR-SLV-75-5
10 µl	400	5800	JR-SLV-75-10
20 µl	400	5800	JR-SLV-75-20
50 µl	400	5800	JR-SLV-75-50
100 µl	400	5800	JR-SLV-75-100
250 µl	240	3500	JR-SLV-75-250
500 µl	240	3500	JR-SLV-75-500
1 ml	240	3500	JR-SLV-75-1000

**Spert parts***for PEEK Cheminert, Rheodyne, and Valco loops*

Sold in packages of 5.

Description	Product No.
PEEK nuts for use with double ferrules, 1/16" hex head long, 10-32	JR-5510-5
PEEK single ferrules, 1/16" (Use 2 per nut, as shown at right.)	JR-5005-5



## SAMPLE LOOPS

### SPECS

- Material:
  - Loops SS316
  - Nuts SS303
  - Ferrules SS316
  - Caps PE (Polyethylene)
- Tolerances:
  - ±0.05 mm (.002")
- Pressure rating:
  - ≤ 1375 bar (≤ 20,000 psi)

- Wide range of volumes available
- Ultra high pressure resistance

VICI Jour stainless steel sample loops are made of high quality SS316. They are cut using an electrolytic cutting process, cleaned with ultrasound and high pressure steam to remove contaminants, dried, and capped with softener-free caps to avoid recontamination of the tubing.

All loops are supplied with stainless hex-head nuts and stainless ferrules for reliable connections.

### Stainless steel loops

*for Cheminert valves*

Fit all Cheminert C1, C1CF, C2, C2V, and C3 type valves.

*Sample size Product No.*

5 µl	JR-SLC-86-5
10 µl	JR-SLC-86-10
20 µl	JR-SLC-86-20
50 µl	JR-SLC-86-50
100 µl	JR-SLC-86-100
250 µl	JR-SLC-86-250
500 µl	JR-SLC-86-500
1 ml	JR-SLC-86-1000



### Stainless steel loops

*for Valco valves*

Fit all 1/16" Valco W and UW type valves.

*Sample size Product No.*

5 µl	JR-SLV-76-5
10 µl	JR-SLV-76-10
20 µl	JR-SLV-76-20
50 µl	JR-SLV-76-50
100 µl	JR-SLV-76-100
250 µl	JR-SLV-76-250
500 µl	JR-SLV-76-500
1 ml	JR-SLV-76-1000



### Spart parts *for stainless Cheminert and Valco loops*

Sold in packages of 5.

Description	Product No.
1/16" stainless nuts, 10-32	JR-ZN1-5
1/16" stainless ferrules	JR-ZF1S6-5

**TECH TIPS AND  
OPTIONS**  
See facing page.



**Stainless steel loops***for Rheodyne valves*

- Fit all Rheodyne valves
- Special loops for valve Model 7725

*Sample size Product No.*

5 µl	JR-SLR-111
10 µl	JR-SLR-112
20 µl	JR-SLR-113
50 µl	JR-SLR-114

100 µl	JR-SLR-115
200 µl	JR-SLR-116
500 µl	JR-SLR-117
1 ml	JR-SLR-118

**Loops for Rheodyne valve 7725, 1/16"***Sample size Product No.*

5 µl	JR-SLR-111-77
10 µl	JR-SLR-112-77
20 µl	JR-SLR-113-77
50 µl	JR-SLR-114-77

100 µl	JR-SLR-115-77
200 µl	JR-SLR-116-77
500 µl	JR-SLR-117-77
1 ml	JR-SLR-118-77

**Spert parts***for stainless Rheodyne loops*

Sold in packages of 5.

<i>Description</i>	<i>Product No.</i>
1/16" Rheodyne type stainless nuts	
For valve 7725, long head	JR-57-5
X-long	JR-58-5
Stainless ferrules	
1/16"	JR-59-5
2-step Rheodyne type, 1/16"	JR-60-5

**TECH TIPS**

- Not recommended for use with halides.
- We offer nuts in three lengths to facilitate wrench access on valves. For loops that will be replaced frequently, we recommend our longest nut – particularly when space is restricted. See long nuts . . . . .31
- The ValvTool tightens nuts with tube or loop installed, providing easy access to many hard-to-reach areas. . . . .88

**OPTIONS**

Other sample sizes are available on request. Please contact your distributor or VICI Jour directly.

**MORE INFORMATION**

About stainless steel. . 100

# Accessories

## TOOLS ■ STARTER KITS ■ BACK-PRESSURE REGULATOR

Our low pressure fittings are ideally suited for flow injection analysis, low pressure liquid chromatography, and stream sampling devices. They may be safely used at pressures up to 500 psi and temperatures to 50°C.

### Accessories index

#### Accessories for tubing

Easy-Flange kit.....	page 85
Easy-Flange combi kit.....	85
PEEK tubing elbows.....	86
Tubing clip - the LC tubing organizer .....	86
LC tubing organizer.....	86
Tubing cutters for polymeric tubing	
Clean-cut.....	86
Guillotine.....	86
Tubing cutters for stainless tubing	
Cutting wheel type .....	87
Plier type .....	87

#### Accessories for fittings

Color-It fingertight adapters.....	17
Work great with all our PEEK and stainless hex head nuts for 1/16" tubing.	
Color-Sleeve adapters.....	17
Designed to be used with our black PPS nuts and one-piece fittings.	

#### Back-pressure regulators .....88

Our back-pressure regulator improves baseline stability by preventing bubble formation in the flow cell.

#### Tools


Ferrule removal kit.....	89
File .....	87
Hex key set .....	89
Tweezer .....	89
Wrench.....	89
ValvTool .....	88
Tubing cutters	
For polymeric.....	86
For stainless .....	87

#### Starter kits



Easy flange combi kit.....	85
Flange-rolling tool, PTFE tubing, fittings, and accessories	
PEEK starter kit.....	90
PEEK fittings and tubing	
Stainless starter kit.....	90
Stainless steel fittings and tubing	
ChromBox kit .....	91
PEEK and stainless fittings and tubing	

#### RELATED PRODUCTS

Color-Sleeve adapters turn PPS nuts into fingertight fittings and color-code different solvent lines. .... 17

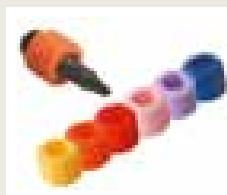

Color-It fingertight adapters work great with all our PEEK and stainless hex head nuts for 1/16" tubing..... 17

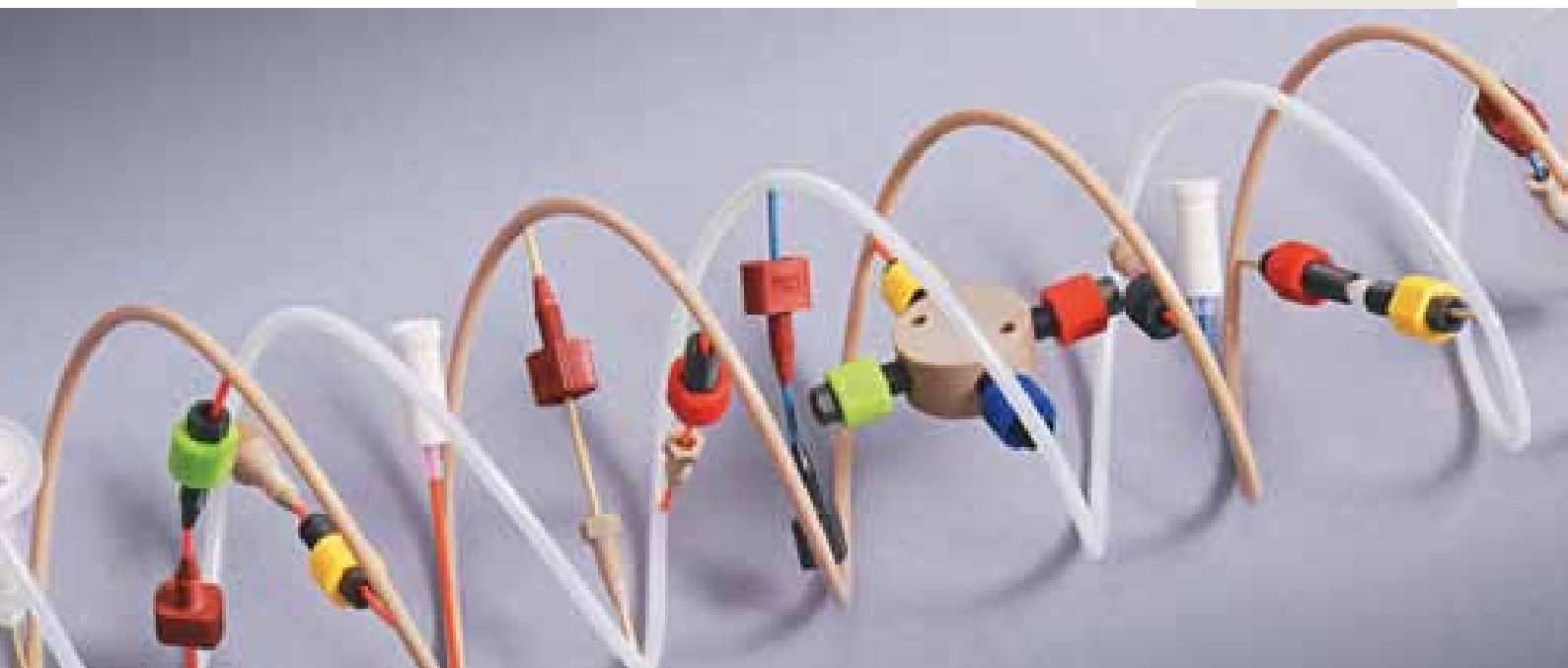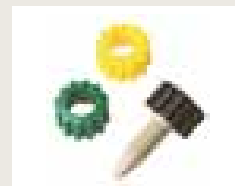



### Easy-Flange kit

- Handy tool for perfect flanges
- No heat or electricity required
- No contamination

The Easy-Flange flange-rolling tool for flanged fittings is suitable for PTFE and PTFE-like tubes with OD between 1/16" and 1/8". The tool forms the flange by applying a mechanical force and does not require electrical power.

The quality of the flange is considerably improved over other available methods, since the tubing is not stressed by heat and because the specially-designed negative conical profile of the flange-forming part yields an ideal shape for optimized sealing properties.

Description		Product No.
Easy-Flange kit		JR-201540
Kit includes:	Qty	Product No.
Flanging clamp	1 pair	JR-201531
Flanging discs with:		
0.5 mm SS pin for PEEK tubing	1	JR-202235
0.8 mm polymer pin	1	JR-201541
0.8 mm titanium pin	1	JR-201554
1.3 mm polymer pin	1	JR-201536
1.3 mm titanium pin	1	JR-201537
Clean-Cut tubing cutter	1	JR-797
Tubing, PTFE, 1/16" x 0.75 mm ID	2 m	JR-T-4036
Plastic box	1	JR-201540L
Flanger/distance gauge	1	JR-201532



### Easy-Flange combi kit

In addition to the contents of the Easy-Flange kit, the combi kit includes 1/8" PTFE tubing, nuts/rings for tube end fittings, and color-coded adapters.

Description		Product No.
Easy-Flange combi Kit		JR-201539
Includes:	Qty	Product No.
Easy-Flange kit above	1	JR-201540
PPS nuts, 1/4-28, flanged		
White, for 1/16" tubing	10	JR-201094WH-10
Black, for 1/8" tubing	10	JR-201580BK-10
PP washers		
For 1/16" tubing	10	JR-201926-10
For 1/8" tubing	10	JR-201928-10
PTFE tubing		
1/16" x 0.75 mm ID	3m	JR-T-4036-M3
1/8" x 1/16" ID	3m	JR-T-680 0-M3

#### OPTIONS

Additional sizes of flanging discs are available.

Flanging disc with

1.0 mm

polymer pin JR-201662

1.5 mm


polymer pin JR-201663

#### MORE INFORMATION

Flanged fittings. . . . 41, 43

Polymeric

tubing. . . . . 9-13

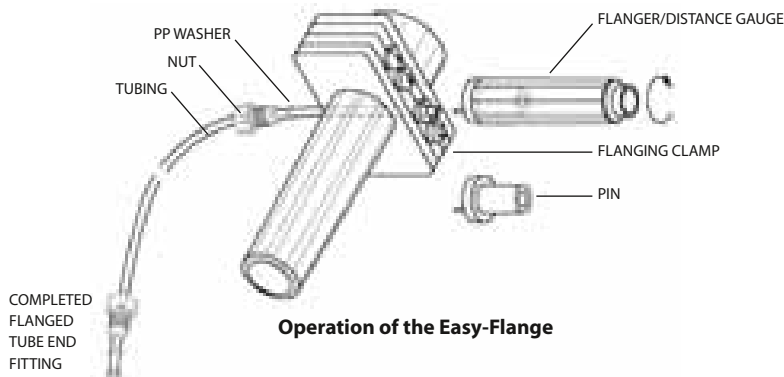


Operation of the Easy-Flange

#### PEEK tubing elbows

Tubing elbows (90° and 180°) are ideal for routing 1/16" PEEK tubing through an LC system. These elbows are proportioned to bend PEEK tubing at the optimum radius for maximum chemical resistance and burst pressure. Installation is simple – just snap the tubing into the elbow.

Sold in packages of 5.

Description	Product No.
90° elbow	JR-357090
180° elbow	JR-357180

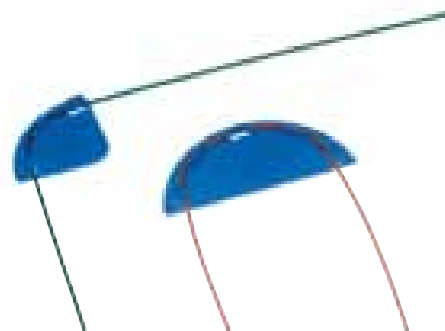

#### Tubing clip – the LC tubing organizer

- Fixes 1/16" and 1/8" tubing to beakers, flasks, or bottles
- Stainless steel spring for long life

This convenient device holds both 1/16" and 1/8" polymer tubes precisely where they should be in your beakers, flasks and bottles etc. up to 4 mm wall thickness. The stainless steel spring will keep your tubing clip working for a long time, and helps keep the solvent line in the solvent.

Description	Product No.
Tubing clip	JR-9001



#### Clean-Cut tubing cutter

*for polymeric tubing*

- Cuts PEEK, PTFE, ETFE, and other polymeric tubing
- No distortion to OD or ID
- Ideal for the most demanding LC applications
- Safety lock secures the blade

Right angles and clean cuts are essential in LC to avoid leaks and dead volume. This handy tool makes it easy to get burr-free, perpendicular ends on polymeric tubing without distorting the outside diameter and/or closing the inside diameter of the tubing.

It is designed to cut PTFE, ETFE and polymers in general and PEEK capillary tubing in particular. The compact design and unique safety locking mechanism allow every chromatographer to keep one in their pocket.

Description	Product No.
Clean-Cut tubing cutter	JR-797
Replacement blade	JR-798



#### Guillotine tubing cutter

*for polymeric tubing*

- Cuts PEEK, PTFE, ETFE, and other polymeric tubing
- Ideal for less demanding LC applications

Description	Product No.
Guillotine tubing cutter	JR-794
Replacement blade	JR-795





### Stainless steel tubing cutter – cutting wheel type

- For 1/16" and 1/8" SS tubing
- Leaves smooth, uniform cuts
- Minimal deburring and/or reaming required
- Easy replacement of the cutting wheel

This is the ideal tool for cutting 1/16" and 1/8" stainless steel tubing, leaving smooth, uniform cuts requiring a minimum of deburring and/or reaming. The easily replaced cutting wheel scores tubing, leaving cuts which allow tube-to-tube and tube-to-frit connections.

Description	Product No.
Tubing cutter for metal tubing	JR-792
Replacement cutting wheel	JR-793



### Stainless steel tubing cutter – plier type

- For easy cutting of 1/16" tubing off a coil
- Reach "hard-to-get" places in an HPLC system
- Delivered in a plastic sheath

These pliers can be used for quick cuts for non-critical applications or at "hard to get" places. However, for scientific equipment, we strongly recommend our pre-cut and cleaned stainless tubing.

Description	Product No.
Tubing pliers for metal tubing	JR-796



### File

This deburring file is handy tool for cleaning up metal tubing ends.

Description	Product No.
Deburring file	JR-802

#### TECH TIP

Between the sample injector, column, and detector, use our pre-cut stainless steel tubing, which has been electrolytically cut and steam cleaned. . . . . 14

#### MORE INFORMATION

Pre-cut stainless steel tubing. . . . . 14

## Backpressure regulator

### SPECS

- For 1/16" OD tubing
- Materials
  - Wetted parts: PEEK/CTFE
  - Housing: PEEK/CTFE
  - Membrane: Polyimide
- 10-32 threads
- Dead volume: 0.25 µl (theoretical)

- Improves baseline stability by preventing bubble formation in the flow cell
- Adjust backpressure without disconnecting anything
- Completely biocompatible fluid path, with low dead volume
- Easy mount with 17 mm panel hole, 9 mm maximum thickness

Even properly degassed mobile phases may contain some dissolved gases that can release bubbles in the detector flow cell resulting in baseline noise and drift. The VICI Jour backpressure regulator provides a quick and convenient way to improve the stability of detector baseline by providing a constant back-pressure on the flow cell, which stops the release of dissolved gases from the mobile phase.

Because of the unique design, the backpressure regulator has a very small internal volume. The peak shape is not affected when the regulator is installed between two detectors in a multi-detector system, or between a detector and a fraction collector or fraction collector system. Typically, 90% or more of a given component's efficiency will be retained upon passing through the regulator.

A spring-loaded diaphragm provides the force that generates the back-pressure. Spring tension may be adjusted using the backpressure adjustment screw on the top of the regulator to vary the backpressure without disconnecting any lines.

Backpressure regulators are available in two ranges: 1-20 bar (15-300 psi) and 20-103 bar (300-1500 psi), and two housing types: PEEK and CTFE.

Range		Housing	Product No.
1-20 bar	15-300 psi	PEEK	JR-BPR1
		CTFE	JR-BPR1KF
20-103 bar	300-1500 psi	PEEK	JR-BPR2
		CTFE	JR-BPR2KF



## ValvTool

- For making connections in hard-to-reach areas
- Ideal for nuts with 1/4" hex and 5/16" heads

The ValvTool is a time-saving device which provides easy access to many hard-to-reach areas. The unique design with its slotted wrench allows tightening of nuts where a loop or a capillary might otherwise make it difficult. The ValvTool is ideal for Valco and Rheodyne® valves and fittings, as well as most of the HPLC and PEEK fittings on the market with 1/4" and 5/16" heads.

Description	Product No.
ValvTool wrench 1/4" and 5/16"	JR-800



**ValvTool shown loosening a nut on a valve**  
(Valve and loop not included)



**Ferrule removal kit**  
For 360 µm and 1/32" ferrules



**Ferrule removal kit**  
For 1/16" and 1/8" ferrules

**Ferrule removal kits**

- Remove polymeric ferrules stuck in fitting details
- One version for 1/32" and 360 micron ferrules
- Another version for 1/16" and 1/8" ferrules

Ferrule removal kits are easy-to-handle tools for removing ferrules stuck in a fitting detail. Two versions cover four different ID's of ferrules.

Description	For use with	Product No.
Ferrule removal kits	360 µm and 1/32" ferrules	JR-FRK1
	1/16" and 1/8" ferrules	JR-FRK2

**Tightening tools for PEEK fittings**

These handy tools make it fast and easy to tighten hex-head fittings. The green tool is for any 1/32" fitting with a 3/16" hex head nut, and the blue version fits the 1/4" hex common in fittings for 1/16" tubing.

Color	For use with	Prod No
Green	1/32" fittings (6-40 threads)	JR-CNFT
Blue	1/16" fittings	JR-ZNFT



**Hex key set**

The hex key set has a wrench to fit any socket head screw on any VICI valve or actuator. Includes the following sizes: .050", 1/16", 5/64", 3/32", 7/64", 1/8", 9/64", and 5/32".

Description	Product No.
Hex ket set	JR-HKS3



**Tweezer**

Ideal for picking up small fittings, etc.

Description	Product No.
Tweezer	JR-806

**Wrenches**

Open end wrenches for nuts with 1/4", 5/16", 3/8", or 7/16" heads.

Description	Product No.
Open end wrench	
	1/4" and 5/16" JR-804
	3/8" and 7/16" JR-805



**RELATED PRODUCTS**  
Deburring file ..... 87  
Tools for tubing .... 86-87

## PEEK starter kit

- All parts biocompatible
- Ideal kit for all HPLC users

The PEEK starter kit is ideal for anyone working in the field of biochromatography. It consists of fittings, ferrules, tubing, and mobile phase filters, as well as a union, our Clean-Cut tubing cutter, and other useful accessories like our tubing elbows – all boxed together for a very competitive price.

For a number of interesting LC applications involving proteins, peptides, nucleic acids, or other samples of biological origin, metal systems may interact with samples or release transition metals that will deactivate columns. Replacing stainless steel parts (such as tubing fittings, ferrules, filters, etc.) with PEEK creates a biocompatible environment for samples and mobile phase.

Description		Product No.	
PEEK starter kit		JR-35P	
Kit includes:		Qty	Product No.
PEEK fittings	One-piece 10-32 fittings	10	JR-55020-5
	Fingertight 10-32 nuts	5	JR-5580-5
	Hex-head long 10-32 nuts	5	JR-5510-5
	1/16" double ferrules	20	JR-5004-10
	Union body, 0.25 mm bore	1	JR-061
PEEK tubing, 1/16" OD			
	0.25 mm ID x 3 m, blue stripe	1	JR-T-6001-M3
	0.50 mm ID x 3 m, orange stripe	1	JR-T-6002-M3
Filters	PEEK in-line filter kit	1	JR-68256
	Last Drop mobile phase filter, 5 µm, PTFE	1	JR-9000-0521
Accessories	Tubing elbow, 90°	2	JR-357090
	Tubing elbow, 180°	2	JR-357180
Tools	Clean-Cut polymeric tubing cutter	1	JR-797
	Tweezer	1	JR-806



## Stainless steel starter kit

UHPLC

The kit consists of stainless fittings, ferrules, tubing, and solvent filters, as well as a union, our cutting-wheel type metal tubing cutter, and other useful items like our ValvTool. All boxed together for a very competitive price.

Description		Product No.	
Stainless steel starter kit		JR-35S	
Kit includes:		Qty	Product No.
Stainless fittings	1/16" nuts, 10-32	10	JR-ZN1-5
	1/16" nuts, 10-32, 19 mm long	10	JR-LZN1-5
	1/16" ferrules	20	JR-ZF1S6-5
	1/16" union with nuts/ferrules, 0.25 mm bore	1	JR-ZU1CS6
Stainless tubing			
	1/16" x 0.25 mm x 3 m	1	JR-T-625-10-M3
	1/16" x 0.50 mm ID x 3 m	1	JR-T-625-20-M3
	1/16" x 0.18 mm ID x 100 mm	2	JR-T-97030
	1/16" x 0.25 mm ID x 200 mm	2	JR-T-97060
Filters	Mobile phase, 2 µm, 1/8" pipe OD	1	JR-3675-2
	Stainless in-line, 0.25 mm ID	1	JR-68230-5
	Last Drop, 10 µm	1	JR-9000-0532
Tools	ValvTool 1/4" x 5/16" wrench	1	JR-800
	Open end wrench, 1/4" x 5/16"	1	JR-804
	Tubing cutter	1	JR-792
	Tweezer	1	JR-806



## RELATED PRODUCTS

### IN THE PEEK KIT

PEEK fittings	
One-piece 10-32 fittings	20
Nuts for double ferrule	
Fingertight	23
Hex-head, long	23
Double ferrules	23
Union body	26
PEEK tubing	
Color-coded striped	6
Filters	
PEEK in-line filter kit	54
Last Drop mobile phase filter	51
Accessories	
Tubing elbows	86
Tools	
Clean-Cut polymeric tubing cutter	86
Tweezer	89

### IN THE STAINLESS KIT

Stainless fittings	
1/16" nuts	31
1/16" ferrules	32
1/16" union	34
Stainless tubing	15
Filters	
Mobile phase	53
Stainless in-line	57
Last Drop	51
Tools	
ValvTool	88
Open end wrench	89
Tubing cutter	87
Tweezer	89



## ChromBox kit

### Includes

- PEEK, PTFE, and stainless tubing
- PEEK and stainless nuts, ferrules, and unions
- Tubing cutters and wrenches
- More accessories, all in a convenient carrying case

The ChromBox is an all-purpose analytical spare parts kit for all our customers using HPLC. It contains the most often used spare parts made both of high quality Type 316L stainless steel as well as PEEK and PTFE.

Description			Product No.	
ChromBox kit			JR-34N	
ChromBox kit includes:				
PEEK fittings	One-piece 10-32 fittings	5	JR-55020-5	
	Fingertight 10-32 nuts	5	JR-5580-5	
	Long hex-head 10-32 nuts	5	JR-5510-5	
	1/16" double ferrules	10	JR-5004-10	
	Union, 0.25 mm bore, body only	1	JR-061	
SS fittings	1/16" nuts, 10-32, 13 mm long	5	JR-MZN1-5	
	1/16" nuts, 10-32, 19 mm long	5	JR-LZN1-5	
	1/16" ferrules	10	JR-ZF1S6-5	
	1/16" union with nuts and ferrules, 0.50 mm bore	1	JR-ZU1MS6	
Tubing				
PEEK	1/16" x 0.25 mm ID x 3 m, blue stripe	1	JR-T-6001-M3	
	1/16" x 0.50 mm ID x 3 m, orange stripe	1	JR-T-6002-M3	
PTFE	1/8" x 1/16" ID x 3 m	1	JR-T-6800-M3	
Stainless	1/16" x 0.25 mm x 3 m	1	JR-T-625-10-M3	
	1/16" x 0.50 mm ID x 3 m	1	JR-T-625-20-M3	
Filters	PEEK in-line filter kit	1	JR-68256	
	Stainless Last Drop, 2 µm	1	JR-9000-0530	
Accessories	Tubing elbow, 90°	1	JR-357090	
	Tubing elbow, 180°	1	JR-357180	
	Tubing clip	2	JR-9001	
Tools	Clean-Cut polymeric tubing cutter	1	JR-797	
	Plier style SS tubing cutter	1	JR-796	
	ValvTool 1/4" x 5/16" wrench	1	JR-800	
	Open end wrench, 1/4" x 5/16"	1	JR-804	
	Deburring file	1	JR-802	
	Tweezer	1	JR-806	

## RELATED PRODUCTS

### IN THE CHROMBOX KIT

PEEK fittings	
One-piece 10-32 fittings	20
Nuts for double ferrule	
Fingertight	23
Hex-head, long	23
Double ferrules	23
Union body	26
Stainless fittings	
1/16" nuts	31
1/16" ferrules	32
1/16" union	34
Tubing	
Color-coded striped PEEK	6
PTFE	12
Stainless	15
Filters	
PEEK in-line filter kit	54
Stainless in-line	57
Accessories	
Tubing elbows	86
Tubing clip	86
Tools	
Clean-Cut polymeric tubing cutter	86
Plier style stainless tubing cutter	87
ValvTool	88
Open end wrench	89
Deburring file	87
Tweezer	89



# General reference

CHEMICAL RESISTANCE ■ MATERIALS ■ GLOSSARY

## Reference index

**General conversion factors** .....page 92

### Tubing charts

Maximum working pressures

Stainless steel tubing .....93

PEEK tubing .....94

Polymeric tubing .....94

Tubing tolerances.....95

Tubing volume/length conversion.....93

**Volume – frits**.....95

**Standard threads** .....95

### Chemical resistance charts

PEEK at high temperature .....96

PEEK and other polymers.....96

### Properties of materials

Metals.....100

Polymers .....99

### Miscellaneous

Glossary .....102

Safety .....101

Trademarks .....101

### Conversion tables

Length .....107

Pressure .....107

Temperature.....107

### Index

By product numbers .....108

Alphabetical by product name .....110

**Order information** .....112

### General conversion factors

inches to mm	inches x 25.4 = mm
feet to meters	feet x 0.3048 = meters
meters to feet	meter x 3.28 = feet

bar to psi	bar x 14.5 = psi
psi to bar	psi x 0.069 = bar
psi to kPa	psi x 6.9 = kPa

### RELATED PRODUCTS

Tubing

ETFE..... 9

FEP

Natural ..... 11

Solid-color..... 11

PEEK ..... 5-8

Dual-layer ..... 7

Natural ..... 5

Solid-color..... 8

Striped ..... 6

PFA..... 9

PTFE ..... 12

Natural ..... 12

Solid-color..... 12

Stainless steel

Premium grade .... 15

Welded..... 15



## Tubing volume/length conversion chart

1/16" OD tubing sold by VICI typically has an ID tolerance of  $\pm 0.025$  mm. This is much tighter than normal commercial grades of tubing; however, it is enough to result in noticeable error if exact volumes are not measured, according to the equation

$$V = \Pi (\frac{1}{2} D)^2 L.$$

Tubing ID		Volume/length	
mm	inches	μl/cm	μl/in
0.050	.002"	0.02	0.05
0.064	.0025"	0.03	0.08
0.075	.003"	0.04	0.12
0.10	.004"	0.08	0.21
0.13	.005"	0.13	0.32
0.17	.0067"	0.23	0.58
0.18	.007"	0.25	0.63
0.25	.010"	0.51	1.29
0.38	.015"	1.14	2.90
0.50	.020"	2.03	5.15
0.64	.025"	3.17	8.04
0.75	.030"	4.56	11.58

Tubing ID		Volume/length	
mm	inches	μl/cm	μl/in
1.00	.040"	8.11	20.59
1.40	.055"	15.39	38.93
1.52	.060"	18.24	46.33
1.59	.062"	19.86	49.47
1.65	.065"	21.38	54.38
1.70	.067"	22.70	57.78
1.78	.070"	24.83	63.06
2.00	.079"	31.42	80.32
2.10	.083"	34.64	88.66
2.16	.085"	36.61	92.99
2.40	.094"	45.24	113.72
3.18	.125"	79.4	201
3.76	.148"	111	282
4.35	.171"	149	376

## Maximum recommended working pressures for stainless steel tubing

PREMIUM GRADE STAINLESS TUBING							
OD	ID (mm)	bar	psi	OD	ID (mm)	bar	psi
1/32"	0.13	1735	25,160	1/8"	0.68	1622	23,520
	0.18	1601	23,225		0.75	1572	22,800
	0.25	1395	20,230		1.00	1406	20,400
	0.50	734	10,645		1.50	1075	15,600
1/16"	0.13	1902	27,580		1.70	960	13,920
	0.18	1835	26,610		2.16	662	9,600
	0.25	1735	25,160	1/4"	4.65	554	8,040
	0.50	1395	20,230				
	0.75	1067	15,480				
	1.00	734	10,645				

WELDED STAINLESS TUBING			
OD	ID (mm)	bar	psi
1/16"	0.13	1871	27,140
	0.18	1794	26,025
	0.25	1717	24,910
	0.50	1384	20,080
	0.75	1076	15,610
	1.00	743	10,780
1/8"	2.10	680	9,860
1/4"	4.65	538	7,800

### WORKING PRESSURE FOR STAINLESS STEEL TUBING

The working pressure for stainless tubing is calculated from the burst pressure by dividing by 2.5. The theoretical burst pressure(s) were calculated using Lamé's formula and a safety factor of 2. Additionally, these results assume the tubing is not subject to any other stresses (e.g. temperature different from ambient, bends) that would influence the calculations. The end user must decide the safe maximum working pressure based on this information and the specific application.

## GENERAL REFERENCE

## Maximum recommended working pressures for PEEK tubing

PEEK TUBING											
OD	ID (µm)	bar	psi	OD	ID (mm)	bar	psi	OD	ID (mm)	bar	psi
360 µm	50	358	5200	1/32"	0.050	455	6600	1/16"	0.050	480	7000
	75	345	5000		0.064	455	6600		0.064	480	7000
	100	300	4400		0.075	455	6600		0.075	470	6800
	125	275	4000		0.100	455	6600		0.100	460	6700
	150	240	3500		0.13	414	6000		0.13	460	6700
					0.15	400	5800		0.18	440	6400
					0.18	386	5600		0.25	420	6100
					0.20	369	5350		0.50	345	5000
					0.25	338	4900		0.75	262	3800
					0.38	269	3900		1.00	179	2600
					0.50	207	3000		1.40	55	800
								1/8"	0.75	379	5500
									1.59	248	3600
									2.00	179	2600
								1/4"	3.17	227	3300

## Maximum recommended working pressures for polymeric tubing

FEP TUBING											
OD	ID (mm)	bar	psi	OD	ID (mm)	bar	psi	OD	ID (mm)	bar	psi
1/32"	0.10	82	1200	1/16"	0.10	124	1800	1/8"	1.59	69	1000
	0.13	78	1100		0.13	124	1800		2.10	41	600
	0.15	78	1100		0.15	117	1700	3/16"	3.18	41	600
	0.18	72	1050		0.18	117	1700	1/4"	4.35	45	650
	0.25	68	980		0.20	117	1700	1.0 mm	0.50	62	900
	0.38	51	750		0.25	117	1700	2.0 mm	1.00	62	900
	0.50	35	500		0.50	96	1400	3.0 mm	1.00	82	1200
					0.75	76	1100		2.00	41	600
					1.00	52	750				

ETFE TUBING			
OD	ID (mm)	bar	psi
1/32"	0.25	152	2200
1/16"	0.17	200	2900
	0.25	186	2700
	0.50	152	2200
	0.75	117	1700
	1.00	83	1200
1/8"	1.59	110	1600
3.0 mm	1.00	138	2000
6.0 mm	4.00	69	1000

PFA TUBING			
OD	ID (µm)	bar	psi
360 µm	50	103	1500
	75	96	1400
	100	90	1300
	150	75	1100
OD	ID (mm)	bar	psi
1/16"	0.50	100	1450
	0.75	76	1100
	1.00	55	800
1/8"	1.59	72	1050
	1.70	65	940
3/16"	3.20	41	600
1/4"	4.35	45	650

PTFE TUBING			
OD	ID (mm)	bar	psi
1/16"	0.18	62	900
	0.25	55	800
	0.50	48	700
	0.75	37	550
	1.00	25	350
1/8"	1.59	35	500
	2.40	18	250
1/4"	4.75	17	250
2.0 mm	1.70	10	150
4.0 mm	3.00	17	250
4.76 mm	3.76	14	200

## WORKING PRESSURE FOR POLYMERIC TUBING

The quoted maximum working pressure(s) were calculated from the burst pressure and the factor 2.5.

## Tolerances – PEEK tubing

Due to tighter tolerances, VICI Jour PEEK tubing is the best choice for chromatography applications in which exact flow rates or tube volume are of importance. The tighter tolerance ensures a more consistent internal volume of transfer lines.

PEEK TUBING		
Tubing OD	OD (±) (mm)	ID (±) (mm)
360 µm	0.025	0.013
1/32"	0.025	0.013
1/16" *	0.025	0.013
1/16" **	0.025	0.025
1/8"	0.050	0.050
* For ID up to 0.25 mm ** For ID over 0.25 mm		

## Tolerances – polymeric tubing

POLYMERIC TUBING		
Tubing OD	OD (±) (mm)	ID (±) (mm)
360 µm	0.025	0.013
1/32"	0.025	0.025
1/16"	0.050	0.025
1/8"	0.10	0.10
3/16"	0.10	0.10
1/4"	0.10	0.10

## Frit volume




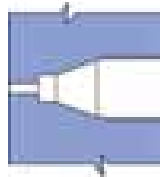

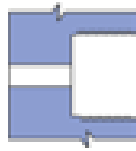
### Stainless, titanium, and polymer frits

Internal frit volumes are theoretical and are calculated by multiplying the total frit volume by the proportion of the pores in the frit (porosity proportion). Stated frit porosities are only nominal and do not reflect the maximum pore size of a frit. The chart below will help you in selecting the right frit for your application.

Porosity in µm	Max pore size in µm	Particles retained in µm	Porosity in %
0.2	1	≥ 1.0	22
0.5	2	≥ 1.5	26
1	3	≥ 3.0	30
2	5	≥ 3.5 - 5.0	35
5	15	≥ 15.0	38
10	30	≥ 30	40

## Standard threads

Following are illustrations of standard threads for VICI Jour fittings. Place your fitting on the drawings below to identify the thread size.

2-56	6-40	10-32	5/16-24	1/4-28	1/2-20
for 360 µm OD tubing	for 1/32" OD tubing	for 1/16" OD tubing	for 1/8" OD tubing	for 1/16" and 1/8" OD tubing	for 1/4" OD tubing
					

## GENERAL REFERENCE

## Chemical resistance of PEEK at high temperature

Chemical	°C	°F	Rating
Acetic acid	200	420	□
Ammonia gas	200	420	■
Carbon monoxide gas	200	420	■
Ethylene glycol	200	420	□
Ethylene glycol (50%)	140	284	■
Hydrogen sulphide gas	200	420	■
Liquid ammonia	200	420	□
Methane	200	420	■
Methylethylketone	200	420	■
Nitrobenzene	200	420	■
Phosphoric acid (50%)	200	420	■
Sodium hydroxide solution	200	420	■
Sulphuric acid (50%)	200	420	□
Sulphur dioxide gas	200	420	■

## Chemical characteristics of PEEK®

PEEK exhibits excellent resistance to a wide range of organic and inorganic chemicals. The compatibility of PEEK with many chemicals at 20°C (68°F) has been investigated and the results for unreinforced grades are given in the table below.

PEEK is compatible with almost any of the solvents used in HPLC. The only solvents which will attack PEEK are concentrated nitric acid and sulphuric acids. However, PEEK tubing can safely withstand 20-30% nitric acid when passivating a system.

Methylene chloride, DMSO, and THF may cause swelling in PEEK. The highest temperature we recommend for PEEK is 100°C. The tubing will maintain the stated pressure rating up to this temperature.

## Chemical resistance of PEEK and other polymers

## A

Resistance at 20°C	PEEK	PA	PE	PP	PPS	FEP	ETFE	PFA	PTFE	NBR	EPDM
Acetaldehyde	■	■	–	–	–	■	–	■	–	■	■
Acetic acid (20 %)	■	■	■	■	■	■	■	■	■	□	■
Acetic acid (80 %)	■	■	■	■	■	■	–	■	■	–	–
Acetic acid (glacial)	■	■	■	■	■	■	■	■	■	■	■
Acetone	■	■	□	■	■	■	■	■	■	■	■
Acetonitrile	■	■	–	–	■	■	■	■	■	□	■
Acrylic acid	■	–	–	–	■	■	–	■	–	–	–
Ammonia, anhydrous	■	□	–	–	■	■	–	■	–	–	–
Ammonia (10 %)	■	□	□	■	■	■	■	■	■	–	–
Ammonia (Liquid)	□	■	–	–	■	■	–	■	–	–	–
Ammonium hydroxide	■	–	■	■	■	■	■	■	■	–	–
Aqua regia	■	■	–	–	■	■	–	■	–	■	■
Aromatic hydrocarbons	■	–	□	■	■	■	–	■	–	–	–

## B

Resistance at 20°C	PEEK	PA	PE	PP	PPS	FEP	ETFE	PFA	PTFE	NBR	EPDM
Benzene	■	■	□	■	■	■	■	■	■	–	–
Benzoic acid	■	–	–	–	■	■	–	■	–	–	–
Benzaldehyde	■	–	–	–	■	■	–	■	–	–	–
Bromine/dibromoethane	■	–	–	–	■	■	–	■	–	–	–
Bromine (dry)	■	–	–	–	■	■	–	■	–	–	–
Bromine (wet)	■	□	–	–	■	■	–	■	–	■	■
Boric acid	■	–	–	–	■	■	–	■	–	–	–
Butanol	■	–	■	■	■	■	■	■	■	–	–

## LEGEND

- Suitable
- Marginal; dependent on application
- Not recommended
- No data available

## C

Resistance at 20°C	PEEK	PA	PE	PP	PPS	FEP	ETFE	PFA	PTFE	NBR	EPDM
Calcium hydroxide	■	–	–	–	■	■	–	■	–	–	–
Carbon tetrachloride	■	–	–	–	■	■	–	■	–	–	–
Chlorine (gas)	■	■	–	–	■	■	–	■	–	■	■
Chlorine (liquid)	■	■	–	–	■	■	–	■	–	–	–
Chloroacetic acid	■	–	□	□	■	■	■	■	■	–	–
Chlorobenzene	■	■	–	–	■	■	–	■	–	–	–
Chloroform	■	■	□	□	■	■	■	■	■	–	–
Cyclohexane	■	■	□	■	■	■	■	■	■	–	–
Cyclohexanone	■	■	■	■	■	■	■	■	■	■	■

## D

Resistance at 20°C	PEEK	PA	PE	PP	PPS	FEP	ETFE	PFA	PTFE	NBR	EPDM
Diethylamine	■	–	■	■	–	■	■	■	■	–	–
Diethylether	■	■	–	–	■	■	–	■	–	–	–
Diethylformamide	■	■	■	■	■	■	■	■	■	■	■
Dioxane	■	■	–	–	■	■	■	■	–	–	–

## E

Resistance at 20°C	PEEK	PA	PE	PP	PPS	FEP	ETFE	PFA	PTFE	NBR	EPDM
Ethanol	■	■	□	■	■	■	■	■	■	■	■
Ether	■	■	□	■	■	■	■	■	■	–	–
Ethyl acetate	■	■	□	■	■	■	■	■	■	■	■
Ethylene chloride	–	■	□	□	■	■	■	■	■	–	–
Ethylene glycol	■	■	■	■	■	■	■	■	■	■	■

## H

Resistance at 20°C	PEEK	PA	PE	PP	PPS	FEP	ETFE	PFA	PTFE	NBR	EPDM
Heptane	■	□	□	□	■	■	■	■	■	■	■
Hexane	■	■	□	□	■	■	■	■	■	–	–
Hydrobromic acid (100%)	■	–	□	□	■	■	–	■	■	–	–
Hydrobromic acid (20%)	■	–	□	■	–	■	■	■	■	–	–
Hydrochloric acid (100%)	■	–	–	□	■	■	■	■	■	–	–
Hydrochloric acid (20%)	■	–	■	□	■	■	■	■	■	–	–
Hydrofluoric acid (100%)	■	–	–	□	■	■	■	■	■	–	–
Hydrofluoric acid (20%)	–	–	■	■	■	■	■	■	■	–	–
Hydrogen peroxide (100%)	■	–	□	□	■	■	–	■	■	–	–
Hydrogen peroxide (50%)	■	–	□	□	–	■	■	■	■	–	–
Hydrogen peroxide (10%)	■	–	■	■	■	■	■	■	■	–	–

## I

Resistance at 20°C	PEEK	PA	PE	PP	PPS	FEP	ETFE	PFA	PTFE	NBR	EPDM
Iso-octane	■	–	–	■	■	■	–	■	■	–	–
Isopropanol	■	□	■	■	■	■	–	■	■	□	■
Isopropyl ether	–	■	■	□	–	■	–	■	■	–	–

## K

Resistance at 20°C	PEEK	PA	PE	PP	PPS	FEP	ETFE	PFA	PTFE	NBR	EPDM
Ketones	■	–	□	□	■	■	–	■	■	–	–

## GENERAL REFERENCE

## M

Resistance at 20°C	PEEK	PA	PE	PP	PPS	FEP	ETFE	PFA	PTFE	NBR	EPDM
Methanol	■	■	■	■	■	■	■	■	■	■	■
Methyl dichloride	–	□	–	■	–	■	–	■	–	–	–
Methyl ethylketone	■	■	□	□	■	■	■	■	■	–	–
Methylene chloride	□	□	□	□	■	■	■	■	■	■	■

## N

Resistance at 20°C	PEEK	PA	PE	PP	PPS	FEP	ETFE	PFA	PTFE	NBR	EPDM
Nitric acid (100%)	■	–	□	■	■	■	■	■	■	–	–
Nitric acid (20%)	■	–	□	■	■	■	■	■	■	–	–

## P

Resistance at 20°C	PEEK	PA	PE	PP	PPS	FEP	ETFE	PFA	PTFE	NBR	EPDM
Pentane	■	–	■	■	–	■	–	■	■	–	–
Perchloric acid	■	–	□	□	–	□	–	□	■	–	–
Phenol (dilute)	■	–	–	–	■	■	–	■	–	–	–
Phenol (concentrated)	■	–	–	–	■	■	–	■	–	–	–
Phosphoric acid (100%)	■	–	□	■	■	■	■	■	■	–	–
Phosphoric acid (40%)	■	–	■	■	■	■	■	■	■	–	–
Potassium hydroxide (dilute)	■	–	–	–	■	■	–	■	–	–	–
Potassium hydroxide (70%)	■	–	–	–	■	■	–	■	–	–	–
Propanol	■	□	–	–	–	■	–	■	–	■	■
Pyridine	■	–	□	■	■	■	■	■	■	–	–

## S

Resistance at 20°C	PEEK	PA	PE	PP	PPS	FEP	ETFE	PFA	PTFE	NBR	EPDM
Sodium hydroxide (80%)	■	□	□	■	■	■	–	■	■	–	–
Sodium hydroxide (20%)	■	□	■	■	■	■	■	■	■	–	–
Sulfuric acid (100%)	■	■	□	□	■	■	■	■	■	–	–
Sulfuric acid (75%)	■	□	■	■	■	■	■	■	■	–	–
Sulfuric acid (40%)	■	■	■	■	■	■	■	■	■	–	–

## T

Resistance at 20°C	PEEK	PA	PE	PP	PPS	FEP	ETFE	PFA	PTFE	NBR	EPDM
Tetrahydrofuran	□	–	□	□	■	■	■	■	■	–	–
Toluene	■	–	□	□	■	■	■	■	■	–	–
Trichloroacetic acid	■	–	■	■	■	■	■	■	■	–	–
Trichloroethane	■	–	–	□	–	■	–	■	■	–	–
Trichloroethylene	■	–	–	–	■	■	–	■	–	–	–
Triethylamine	–	–	–	■	–	■	■	■	■	–	–

## W

Resistance at 20°C	PEEK	PA	PE	PP	PPS	FEP	ETFE	PFA	PTFE	NBR	EPDM
Water	■	■	■	■	■	■	■	■	■	■	■
Water (distilled)	■	–	■	■	■	■	■	■	■	–	–

## X

Resistance at 20°C	PEEK	PA	PE	PP	PPS	FEP	ETFE	PFA	PTFE	NBR	EPDM
Xylene	■	–	–	–	■	■	–	■	–	–	–

## LEGEND

- Suitable
- Marginal; dependent on application
- Not recommended
- No data available



## Properties of polymers

### CTFE

Chlorotrifluoroethylene, is the generic name for the material produced as Kel-F®. It is very resistant to all chemicals except THF and some halogenated solvents, and is resistant to all inorganic corrosive liquids, including oxidizing acids. CTFE can be used at temperatures up to 100°C. Swells in ketones.

### EPDM

Ethylene Propylene Diene Monomer. Ethylene-propylene synthetic rubbers offer excellent heat resistance as well as resistance to aging from oxidation, ozone, and weather due to their stable, saturated polymer backbone structure. As non-polar elastomers, they have good electrical resistivity as well as resistance to polar solvents such as water, acids, alkanes, phosphates, esters, and many ketones and alcohols.

### ETFE

Ethyltrifluoroethylene is the generic name for the material such as Tefzel®. A fluoropolymer used for sealing surfaces, it is resistant to most chemical attack; however, some chlorinated chemicals will cause a physical swelling of ETFE tubing.

### FEP

Fluorinated ethylene propylene is another member of the fluorocarbon family with similar chemical properties. It is generally more rigid than PTFE, with somewhat increased tensile strength. It is typically more transparent than PTFE, slightly less porous, and less permeable to oxygen. FEP is not as subject to compressive creep at room temperature as PTFE, and because of its slightly higher coefficient of friction is easier to retain in a compression fitting.

### FFKM

FFKM (perfluorinated rubber) is a perfluoroelastomer, a class of materials such as Kalrez®. It is chemically resistant to polar solvents, organic solvents, inorganic and organic acids and bases, fuels, oils, lubricants, inorganic salts, aldehydes, metal halogen compounds, chlorine, sodium hydroxide, aromatics, alcohols, steam, and strong oxidizing agents. The normal temperature service range is -40°C to 315°C, with up to 343°C in intermittent service.

### PEEK™

Considered relatively inert and biocompatible, polyetheretherketone tubing can withstand temperatures up to 100°C. Under the right circumstances, .005" - .020" ID tubing can be used up to 5000 psi for a limited time, and 0.030" to 3000 psi. Larger IDs are typically good to 500 psi. These limits will be substantially reduced at elevated temperatures and in contact with some solvents or acids.

Its mechanical properties allow PEEK to be used instead of stainless in many situations and in some environments where stainless would be too reactive. However, PEEK can be somewhat absorptive of solvents and analytes, notably methylene chloride, DMSO, THF, and high concentrations of sulfuric and nitric acid. This tubing is highly prone to "kinking," or sealing off, if held in a sharp bend over time.

### PFA

Perfluoroalkoxy is a fluorocarbon with chemical and mechanical properties similar to FEP. More rigid than either PTFE or FEP. Commonly used for injection molded parts.

### PPS

Polyphenylene sulphide is the generic name for the material produced as Fortron®, Ryton®, and Techtron®. It is very resistant to all solvents, acids, and bases.

### PTFE

Polytetrafluoroethylene is the generic name for the class of materials such as Teflon®. It offers superior chemical resistance but is limited in pressure and temperature capabilities. Because it's so easy to handle, it is often used in low pressure situations where stainless steel might cause adsorption. PTFE tubing is relatively porous, and compounds of low molecular weight can diffuse through the tubing wall.

### Polyacetal

Polyacetal or polyoxymethylene (POM) is the generic name for the material produced as Delrin®. It is an engineering polymer with high mechanical strength and rigidity, excellent dimensional stability and excellent resistance to moisture and solvents. It has a wide useful temperature range in air of -50°C to +90°C. It is not suitable for use with acids, alkalis, or oxidizing agents, and has limited resistance to dichloromethane and methyl ethyl ketone.

*Continued, next page.*

## GENERAL REFERENCE

### Polyamide

Polyamide is a thermoplastic polymer also known as nylon. Nylon has very good physical properties, but moisture can have a significant effect. It exhibits very good heat resistance and excellent chemical and wear resistance.

### Polyethylene

Polyethylene is a semi-crystalline material with excellent chemical resistance, plus good fatigue and wear resistance. Polyethylenes provide good resistance to organic solvents and staining, and have low moisture absorption rates.

### Polypropylene

Widely used polymer for non-wetted parts. Attacked by strong oxidizers as well as aromatic and chlorinated hydrocarbons.

### PVDF

PVDF, or polyvinylidene fluoride, has excellent resistance to most mineral and organic acids, aliphatic and aromatic hydrocarbons, and halogenated solvents. Poor resistance to acetone, MEK, THF, and potassium and sodium hydroxide. Often supplied as Kynar®.

## Properties of metals

### Stainless steel, Type 316

This is the standard tubing material for chromatography, suitable for a wide variety of applications. It is cold drawn seamless, not welded, with close tolerances held on both ID and OD. We neither recommend nor offer Type 304 stainless steel for analytical applications.

Austenitic stainless steels may be used for most chromatographic applications. Type 316 is most commonly used for HPLC because of its superior chloride ion resistance.

### Stainless steel, Type 303

Recommended for GC use and general purpose connections, combining excellent machining characteristics with good resistance to corrosion and high temperature oxidation. Susceptible to attack by chlorides, iodides, and bromides.

### Titanium

Although it is more difficult to machine than common alloys containing aluminum and vanadium, VICI uses Grade 2 pure titanium in order to avoid possible contamination of the sample stream with these metals. Good for organic and inorganic salts except aluminum and calcium chlorides, and all alkalis except boiling concentrated potassium hydroxide.

Good with dilute, low temperature formic, lactic, sulfuric, hydrochloric, and phosphoric acids, but rapidly attacked by hydrofluoric acid. Good with dilute nitric acid at low temperatures; corrodes at high concentrations and temperatures. Can ignite with fuming nitric acid. Attacked by oxalic acid, concentrated phosphoric acid, hot trichloroacetic acid, and zinc chloride.

## Safety

1. Never tighten or loosen a fitting or valve connection while it is pressurized. Provisions should be made within the system to release pressure via suitable valve components.
2. Do not exceed pressure or temperature specifications. Note that in many cases, the system pressure is limited by the tubing used, not the fittings.
3. The use of toxic or hazardous fluids requires extra caution during operation or maintenance. The user is responsible for insuring safe operation and for understanding the nature of the fluids and chemistry involved.
4. The use of thread lubricants or sealants is required only on tapered pipe threads. These sealants and lubricants may have different temperature limits or chemical compatibility than the valves or fittings.

## Trademarks

Trademarks of materials and products found in VICI Jour products or referred to in this website. Please contact us regarding errors, omissions, or discrepancies.

Cheminert®	Valco Instruments Co. Inc. and VICI AG International
Delrin®	E.I. duPont de Nemours
Duran®	Schott AG
Fortron®	Fortron Industries LLC
Kalrez®	DuPont Dow Elastomers
Kel-F®	3M Company
Kynar®	Elf Atochem North America Inc.
Nanovolume®	Valco Instruments Co. Inc.
PEEK™	Victrex plc
PEEK HT™	Victrex plc
Rheodyne®	Rheodyne LLC
Ryton®	Chevron Phillips Chemical Company LLC
SCAT®	SCAT-Europe GmbH
Techtron®	Quadrant Engineering Plastic Products
Teflon®	E.I. duPont de Nemours
Valco®	Valco Instruments Co. Inc.
VICI®	Valco Instruments Co. Inc. and VICI AG International
VICI Jour®	Valco Instruments Co. Inc. and VICI AG International

### CAUTION

The improper selection or use of components or systems described herein can cause personal injury or property damage.

The system designer and user are solely responsible for the selection of products suitable for the specific requirements of the application, as well as proper installation, operation, and maintenance of these products.

Compatibility with hazardous fluid streams, environmental conditions, and mechanical requirements are the responsibility of the user.

## A

**Adapter**

A type of fitting which provides a method of joining two components of differing thread types or systems.

**Analytical column**

A long narrow tube packed or coated with one of many available chemically diverse compounds that can separate the components in a sample according to their boiling point, polarity, molecular size, or combination thereof. A column of some kind is used with most chromatographic techniques.

## B

**Biocompatibility**


Defines the materials used in a system (i.e. fittings, tubing, and valves) that do not change the bioactivity of the biological substances that come into contact with the surface of these materials. Note that in chromatographic systems, the tubing and column contribute over 99% of the surface area and the valves and fittings are insignificant.

**Bore**

The diameter of the minimum orifice through the fitting. See capillary bore, through-type bore, and large bore.

**Bulkhead fitting**

A type of fitting in which the fitting body is inserted through an instrument panel or mounting bracket, to which it is affixed with a mounting nut. The Valco fitting body is uniquely undercut so that it "bites" into the panel when the mounting nut is tightened, eliminating the need for a lock washer.

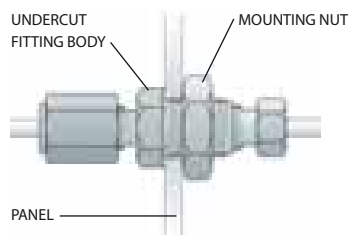**Butt connection**

A type of connection in which the two tube ends are directly and squarely in contact, usually effected with a through-type union. Typically used with fused silica connections, or small bore metal tubing.

## C

**Cap**

A cap is used to dead-end a piece of tubing with a nut and ferrule attached.

**Capillary bore**

The smallest available standard orifice in a given fitting design (usually 0.25 mm). Typically denoted by suffix "C" in the product number.

**Collapsible ferrule**

A one-piece ferrule engineered to collapse as it is tightened. The collapse takes place in two very narrow areas, and results in a very effective seal with virtually no distortion of the tubing ID and no dead volume. Patent No. 6,575,501

**Compression fitting**

A style of fitting in which a threaded nut compresses a tapered ferrule onto tubing as the nut is tightened. Valco metal ferrules cut a ring into the tubing wall while polymer types rely on surface compression to form a seal.

**Connecting volume**

The volume between two or more connections. This may be cleanly swept, thus not contributing to peak distortion, or may be "dead volume" such as that found in fittings with larger bores than the connecting tubing.


**Cross**

A type of distribution fitting which connects four pieces of tubing, arranging them in the pattern of a cross.

## D

**Dead volume**

Any volume which a component introduces to a system that is not cleanly swept and relies on diffusion to clear the space. See connecting volume.

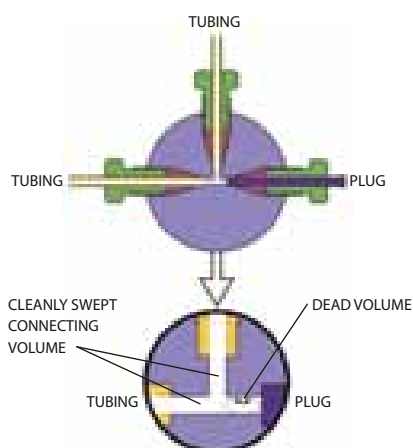

**Detail**

See *fitting detail*.

**Distribution fitting**

A generic term for tees, crosses, and manifolds, used to provide multiple access points to “distribute” a gas or liquid through a system. CAUTION! Using a distribution fitting in reverse to coalesce multiple streams may create dead volume. Special manifolds are available for this application.

**F****FIA**

Flow Injection Analysis.

A simple and versatile analytical technique for automating wet chemical analyses based on the manipulation of a sample zone formed from the injection of the sample into a continuous stream of fluid used as a carrier.

**Female nut**

A nut with female threads; also called external nut.

**Ferrule**

One of the components of a compression fitting; the conical piece of metal or plastic that compresses onto the tube as it is forced into a tapered seat. Valco metal ferrules are unique in that they attach to and seal at the tube by cutting a shallow ring into it, instead of by actually swaging it. This is preferable since it introduces no flow restriction.

**Filter**

A type of union or reducing union which traps the particulates in a stream. The filtering element is typically a mesh screen or sintered frit.


**Fingertight fitting**

A fingertight fitting is a knurled or winged nut which is hand tightened to compress the ferrule. Note: Polymeric ferrules rely on surface compression to form a seal and hold the tube by friction.

**Fitting detail**

A fitting detail includes the thread to accept the nut, a tapered ferrule seat, and the pilot for the tube. The nut, ferrule, and tube go into the fitting detail to make the seal.

Valco fitting details have the perfect shape, with a smooth transition from the ferrule seat into the cylindrical and threaded part. Many other fitting details are poor in design where the ferrule is forced in a too shallow ferrule seat. The result is a “nicked” ferrule which can no longer apply its sealing force to the front of the ferrule.

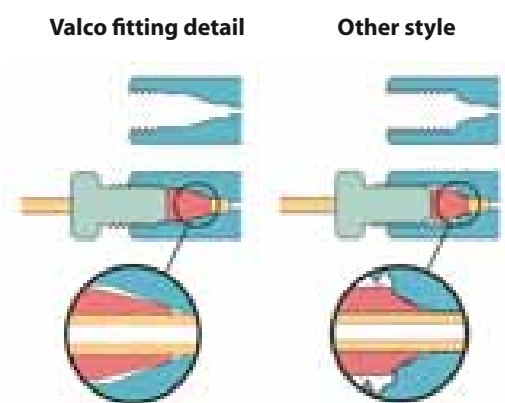**Flanged fitting**

A type of fitting used with fluoropolymer tubing (PTFE, FEP) in which a flange is made at the tube end. Connections are made at the flange either by compressing the flange into a flat detail (typically 1/4-28 threaded) or by butting two flanges together. A special flanging tool forms the flanges.

**Flangeless fitting**

Similar in application to the flanged fitting, but the flange is not required. A ferrule system is used which grips/compresses the tube. This fitting type can be used with virtually any polymeric tubing since the tube end does not have to be formed, but simply square cut. Typically used in 1/4-28 threaded fittings, it is usually interchangeable with flanged fittings.

**Frit**

A filter element typically made of stainless, Hastelloy, Titanium, or polymers, usually 0.75 mm or 1 mm thick. Frits may provide better filtration than screens, but because they are thicker there is greater mixing potential, and they typically result in increased pressure drop.

**G****GC**

Gas Chromatography.

An analytical method incorporating an injection system, analytical column, controlled temperature zone, and detector. An inert carrier gas moves the sample through the column, which separates the sample components into discrete bands which are measured as they pass through the detector.

**Guard column**

A column used in series between the injector and analytical column to prevent certain types of components from entering the analytical column.

**H****HPLC**

High Performance Liquid Chromatography.

An analytical system consisting of an injector, pump, analytical column, and detector. Using a liquid mobile phase, the sample is pumped through the column, where it is separated into discrete sample component bands which are detected and measured as the bands elute from the column.

**I****ID**

Internal diameter.

**Inert**

Technically, unreactive with other substances; however, in the instrumentation field, "inert" is a relative term. Often polymers are termed inert but are soluble in some fluids and can react with some compounds.

**Internal fitting**

A type of compression fitting in which the fitting body has female threads; an internal nut has male threads.

**L****LC**

Liquid Chromatography.

Any of a variety of low to medium pressure techniques which use a liquid mobile phase as the carrier to move sample. Similar to HPLC.

**Large bore**

A bore that is larger than the standard for a given fitting; a fitting ordered with a large bore will have a larger flow orifice than the standard or capillary bore fitting of the same design. Denoted by suffix "L" in the product number.

**Luer adapter**

An adapter that connects a tapered luer fitting (square nib) of a syringe to a tube or tube fitting.

**M****Male nut**

A nut with male threads.

**Make up**

The point at which a ferrule, nut, and tube are assembled in the fashion which will effect a leak-free seal. In most compression fittings, that is accomplished by compressing the tube with the small end of the ferrule. With Valco metal ferrules, the ferrule usually makes up on the tube by cutting a shallow ring in it.

**Manifold**

A type of distribution fitting in which a single source is directed to multiple outlets, or vice versa. Caution! Using a common distribution fitting in reverse to merge multiple streams may create dead volume. Special manifolds are available for this application.

**Microbore column**

A liquid chromatography column of narrow bore (typically 2 mm or less) for improved resolution.

**N****NPT**

National Pipe Thread; a standardized tapered pipe fitting. See pipe thread.

**Nanobore**

Nanobore generally refers to components with bore sizes less than 250  $\mu\text{m}$  (0.010").

**Nanovolume®**

The trademark registered to Valco Instruments Co. Inc. and VICI AG International, applied to our nanobore components with bore sizes less than 250  $\mu\text{m}$  (0.010").

**Nut**

The tensioning component of a compression fitting. As the threaded nut is tightened into the fitting detail, it pushes the ferrule forward into the tapered ferrule seat, causing it to make up on the tube.

**O****OD**


Outside diameter.

**One-piece fingertight**

A one-piece fingertight fitting is a knurled or winged nut with integrated ferrule. It is convenient to use since the ferrule is retrieved from the port with the nut. Note: Polymer type ferrules rely on surface compression to form a seal and hold the tube by friction.

**P****Pilot**

The tubing which extends beyond the ferrule in a made-up fitting, or the integral portion of a ZRF internal reducing ferrule which extends beyond the ferrule. See also pilot depth, below.

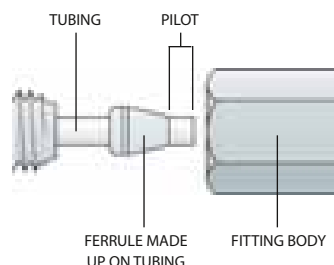**Pilot depth**

The length of the tubing diameter cavity beyond the tapered ferrule seat within a fitting detail. Valco fitting pilot depths are tightly controlled to facilitate the interchangeability of components without the risk of leaks or dead volume.

**Pipe thread**

The external or internal threads of a fitting designed to effect a metal-to-metal seal on the conical thread faces. This type of fitting does not "bottom out" in the detail. Typically used with PTFE tape or other compound to lubricate the threads; however, since the diffusion rate of air components through the PTFE tape is considerable, pipe fittings should not be used in systems where leakage rates are critical.

**Plug**

A plug is used to block a fitting detail in a union, tee, cross, or valve.

**R****Reducing ferrule**

A ferrule which allows a smaller tube to be used in a fitting detail designed for a larger tube. Caution should be taken if standard reducing ferrules without integral pilots are used, since dead volume may be created in the fitting pilot depth.

**Reducing union**

A fitting which joins two tubes of different ODs. The bore of the fitting should typically match the ID of the smaller tube.

**S****SFE**

Supercritical Fluid Extraction.

An extraction technique using a fluid in its supercritical state as the extraction medium. Some liquids and mixtures maintained above a critical temperature and pressure exhibit properties of both the liquid and gas phases of the element. These are defined as supercritical. CO<sub>2</sub> is a common supercritical fluid. Extreme caution must be used with supercritical CO<sub>2</sub>, since uncontrolled expansion (leaks) can be very hazardous due to the substantial stored energy.

**SFC**

Supercritical Fluid Chromatography.

An analytical technique using a supercritical fluid (see SFE) as the mobile phase/carrier.

**Standard bore**

A bore which was chosen as the standard for a particular fitting, typically based on the most common tubing ID used with that fitting.



## T

**Tee**

A type of distribution fitting which connects three pieces of tubing, arranging them in the pattern of a "T".

**Through-type bore**

A bore which is slightly larger than the OD of the tubing which is used with the given fitting. A union with a through-type bore allows the tube ends to butt directly together, or for one tube to run completely through the fitting. Denoted by suffix "T" in the product number. In order to assure correct pilot lengths, we recommend that ferrules be made up on the tubing in a standard union.

## U

**Union**

A fitting for connecting two pieces of tubing of the same OD.

*See also reducing union.*

**Female union**

A type of compression fitting in which the fitting body has female threads. Also referred to as internal union.

**Male union**

A type of compression fitting in which the fitting body has male threads. Also referred to as external union.

**Unswep volume**

The volume of any portion of a fitting which is in the flowpath but which is a different diameter than the primary flow orifice through the tubing/fitting assembly, or any area not directly swept by the fluid flow. This can also be known as "dead volume" if it is very poorly swept.

## V

**VICI**

VICI® is an acronym for Valco Instruments Company Inc.

VICI now defines multiple companies which manufacture various products for the analytical industry.

## W

**Wetted surfaces**

The surfaces which are contacted by the sample stream.

## Y

**Y ("Wye")**

A type of distribution fitting which connects three pieces of tubing, arranging them in the pattern of a "Y". Occasionally referred to as a "wye".


## Z

**ZDV**

Zero dead volume. See definition below.

**Zero dead volume**

Describes a connection which does not add volume to the system beyond what an extension of tubing would in its place.

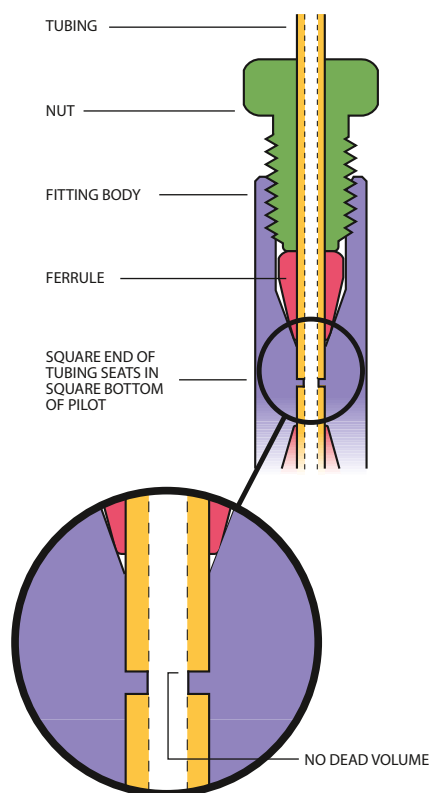**Zero volume**

While often used interchangeably with zero dead volume, it ideally describes a fitting design in which there is no internal volume, such as a through-type union designed to butt-fit two pieces of tubing.



## Temperature conversions

°C	°F	°C	°F	°C	°F	°C	°F	°C	°F
-40	-40	65	149	170	338	275	527	500	932
-35	-31	70	158	175	347	280	536	525	977
-30	-22	75	167	180	356	285	545	550	1022
-25	-13	80	176	185	365	290	554	575	1067
-20	-4	85	185	190	374	295	563	600	1112
-15	5	90	194	195	383	300	572	625	1157
-10	14	95	203	200	392	305	581	650	1202
-5	23	100	212	205	401	310	590	675	1247
0	32	105	221	210	410	315	599	700	1292
5	41	110	230	215	419	320	608	725	1337
10	50	115	239	220	428	325	617	750	1382
15	59	120	248	225	437	330	626	775	1427
20	68	125	257	230	446	335	635	800	1472
25	77	130	266	235	455	340	644	825	1517
30	86	135	275	240	464	345	653	850	1562
35	95	140	284	245	473	350	662	875	1607
40	104	145	293	250	482	375	707	900	1652
45	113	150	302	255	491	400	752	925	1697
50	122	155	311	260	500	425	797	950	1742
55	131	160	320	265	509	450	842	975	1787
60	140	165	329	270	518	475	887	1000	1832

## Length

µm	inches
100	.004"
150	.006"

mm	inches
0.050	.002"
0.064	.0025"
0.075	.003"
0.100	.004"
0.13	.005"
0.15	.006"
0.18	.007"
0.20	.008"
0.25	.010"
0.38	.015"
0.40	.016"
0.50	.020"
0.75	.030"
1.00	.039"
1.50	.059"
2.00	.079"
3.00	.118"
4.00	.157"
4.60	.181"
6.00	.236"
6.35	.250"
6.40	.252"
7.00	.276"
10.00	.394"
12.70	.500"

inches	mm
1/32"	0.79
1/16"	1.59
3/32"	2.38
1/8"	3.18
1/4"	6.35
5/16"	7.94
3/8"	9.50
1/2"	12.70
1"	25.40

feet	mm
1 ft	304.8

## Pressure conversions

bar	psi	kPa	MPa	bar	psi	kPa	MPa	bar	psi	kPa	MPa
0.1	1.5	10	0.01	6	87.0	600	0.6	100	1,450	10,000	10
0.2	2.9	20	0.02	7	101.5	700	0.7	200	2,900	20,000	20
0.3	4.4	30	0.03	8	116.0	800	0.8	300	4,350	30,000	30
0.4	5.8	40	0.04	9	130.5	900	0.9	400	5,800	40,000	40
0.5	7.3	50	0.05	10	145	1,000	1	500	7,250	50,000	50
0.6	8.7	60	0.06	20	290	2,000	2	600	8,700	60,000	60
0.7	10.2	70	0.07	30	435	3,000	3	700	10,150	70,000	70
0.8	11.6	80	0.08	40	580	4,000	4	800	11,600	80,000	80
0.9	13.1	90	0.09	50	725	5,000	5	900	13,050	90,000	90
1	14.5	100	0.1	60	870	6,000	6	1,000	14,500	100,000	100
2	29.0	200	0.2	70	1,015	7,000	7	1,100	15,950	110,000	110
3	43.5	300	0.3	80	1,160	8,000	8	1,200	17,400	120,000	120
4	58.0	400	0.4	90	1,305	9,000	9	1,300	18,850	130,000	130
5	72.5	500	0.5								

## GENERAL REFERENCE

JR-5FR1-5.....73	JR-082.....35	JR-201928-10.....41	JR-4676-10TF.....51	JR-5511-5.....23	JR-68264.....55	JR-BPR1KF.....88
JR-5FR1TI-5.....73	JR-083-5.....35	JR-201999.....64	JR-4676-2.....51	JR-5520-10.....19	JR-68264P.....55	JR-BPR2KF.....88
JR-5FR2-5.....73	JR-086.....35	JR-201999-10.....48	JR-4676-2.5TF.....51	JR-5525-10.....19	JR-68330-05.....57	JR-C360C.....39
JR-5FR4-5.....73	JR-087-5.....35	JR-202118-10.....43	JR-4676-5TF.....51	JR-5570-5.....22	JR-68330-2.....57	JR-C360CFS.....39
JR-5FR6-5.....73	JR-096.....35	JR-202119-10.....43	JR-4677-10.....51	JR-5580-5.....23	JR-68330-5.....57	JR-C360CPK.....38
JR-031-10.....41	JR-097-5.....35	JR-202120-10.....43	JR-4677-10TF.....51	JR-5595.....29	JR-68331-05.....57	JR-C360ET.....38
JR-032-10.....41	JR-1030.....28	JR-202121-10.....43	JR-4677-2.....51	JR-56-5.....31	JR-68331-2.....57	JR-C360NFFS.....39
JR-040 RE-10.....42	JR-1032.....28	JR-202122-10.....43	JR-4677-2.5TF.....51	JR-57-5.....31	JR-68331-5.....57	JR-C360NFPK.....38
JR-040BK-10.....42	JR-1033.....28	JR-202123-10.....43	JR-4677-5TF.....51	JR-58000-5.....19	JR-69230-05.....56	JR-C360NFS6.....39
JR-040BL-10.....42	JR-1040.....28	JR-202124-10.....43	JR-5003-5.....22	JR-58-5.....31	JR-69230-2.....56	JR-C360PPK.....38
JR-040GR-10.....42	JR-1042.....28	JR-202125-10.....43	JR-5004-5.....23	JR-59-5.....32	JR-69230-5.....56	JR-C360QTPK2.....38
JR-040WH-10.....42	JR-1043.....28	JR-202142-10.....43	JR-5010-5.....22	JR-6029.....46	JR-69330-05.....56	JR-C360QTPK4.....38
JR-040WL-10.....42	JR-1061.....26	JR-202143-10.....43	JR-55010-5.....17	JR-6029.....46	JR-69330-2.....56	JR-C360QTPK6.....38
JR-040YL-10.....42	JR-1064.....44	JR-202145-10.....43	JR-55011-5.....17	JR-6030.....46	JR-69330-5.....56	JR-C360QUPK2.....38
JR-041-10.....41	JR-1065.....44	JR-202235.....85	JR-55012-5.....17	JR-6032.....46	JR-792.....87	JR-C360QUPK4.....38
JR-042-10.....41	JR-1066.....26	JR-2030.....28	JR-55013-5.....17	JR-6033.....46	JR-793.....87	JR-C360QUPK6.....38
JR-049-10.....41	JR-1067.....26	JR-2032.....28	JR-55014-5.....17	JR-6034.....46	JR-794.....86	JR-C360QXPk2.....38
JR-050 RE-10.....42	JR-1068.....44	JR-2033.....28	JR-55015-5.....17	JR-6039.....46	JR-795.....86	JR-C360QXPk4.....38
JR-050BK-10.....42	JR-10FR1-5.....73	JR-2040.....28	JR-55016-12.....17	JR-6040.....46	JR-796.....87	JR-C360QXPk6.....38
JR-050BL-10.....42	JR-10FR1K-5.....73	JR-2042.....28	JR-55020-5.....20	JR-6041.....46	JR-797.....86	JR-C360RU.5FS2.....39
JR-050GR-10.....42	JR-10FR2-5.....73	JR-2043.....28	JR-55021-5.....20	JR-6042.....46	JR-798.....86	JR-C360RU.5FS4.....39
JR-050WH-10.....42	JR-10FR4-5.....73	JR-26501.....79	JR-55021-5.....26	JR-6043.....46	JR-800.....88	JR-C360RU.5FS6.....39
JR-050WL-10.....42	JR-10FR6-5.....73	JR-26502.....79	JR-55022-5.....20	JR-6044.....46	JR-802.....87	JR-C360RU.5S62.....39
JR-050YL-10.....42	JR-10FR8TI-5.....73	JR-26503.....79	JR-55023-5.....20	JR-60-5.....32	JR-804.....89	JR-C360RU.5S64.....39
JR-051-10.....41	JR-1100-05P-5.....74	JR-26504.....79	JR-55024-5.....20	JR-6120-50030.....70	JR-805.....89	JR-C360RU.5S66.....39
JR-060-5.....43	JR-1100-2P-5.....74	JR-27305.....78	JR-5502-5.....20	JR-6140.....13	JR-80-5.....35	JR-C360RU1FS6.....39
JR-061.....26	JR-1101-05P-5.....74	JR-2FR1-5.....73	JR-55025-5.....20	JR-6141-10.....13	JR-806.....89	JR-C360RU1PK6.....38
JR-0611.....29	JR-1101-2P-5.....74	JR-2FR1K-5.....73	JR-5503-5.....23	JR-65001NF.....77	JR-9.....78	JR-C360RU1S66.....39
JR-0611 TI05-3.....54	JR-1102-05P-5.....74	JR-2FR1KI-5.....73	JR-5504-5.....20	JR-65002NF.....77	JR-9000-0005.....62	JR-C360UFS2.....39
JR-0611-SS05-3.....54	JR-1102-2P-5.....74	JR-2FR1TI-5.....73	JR-5505.....29	JR-65003NF.....77	JR-9000-0007.....69	JR-C360UFS4.....39
JR-0611-SS2-3.....54	JR-1103-05P-5.....74	JR-2FR2-5.....73	JR-55050-10.....25	JR-65005NF.....77	JR-9000-0010.....62	JR-C360UFS6.....39
JR-0611-TI2-3.....54	JR-1103-2P-5.....74	JR-2FR2TI-5.....73	JR-55050-10.....41	JR-65006NF.....77	JR-9000-0011.....70	JR-C360UPK2.....38
JR-0612.....29	JR-1104-05P-5.....74	JR-2FR4-5.....73	JR-55051-10.....41	JR-65011.....78	JR-9000-0012.....70	JR-C360UPK4.....38
JR-0617.....29	JR-1104-2P-5.....74	JR-2FR4TI-5.....73	JR-55060-10.....41	JR-65158.....78	JR-9000-0013.....69	JR-C360UPK6.....38
JR-062.....29	JR-1110-05P-5.....74	JR-2FR6-5.....73	JR-55061-10.....41	JR-660100.....49	JR-9000-0014.....69	JR-C360US62.....39
JR-0622.....29	JR-1110-2P-5.....74	JR-2FR6TI-5.....73	JR-55070-10.....41	JR-660110.....49	JR-9000-0520.....51	JR-C360US64.....39
JR-0623.....29	JR-1110-5P-5.....74	JR-2FR8TI-5.....73	JR-55071-10.....41	JR-660200.....49	JR-9000-0520B.....52	JR-C360US66.....39
JR-0625.....29	JR-1111-05P-5.....74	JR-32171.....52	JR-55072-10.....41	JR-660210.....49	JR-9000-0520F.....51	JR-C5M1PK.....46
JR-0626.....29	JR-1111-10P-5.....74	JR-32172.....52	JR-55073-10.....41	JR-660300.....49	JR-9000-0521.....51	JR-C5M2PK.....46
JR-064.....44	JR-1111-2P-5.....74	JR-32174.....52	JR-5507-5.....20	JR-660310.....49	JR-9000-0521F.....51	JR-C9M1PK.....46
JR-0640.....45	JR-1111-5P-5.....74	JR-32178.....52	JR-55080-10.....17	JR-660320.....49	JR-9000-0522.....51	JR-C9M2PK.....46
JR-0641-SS05.....54	JR-1125-05P-5.....75	JR-32179.....52	JR-55080-5.....17	JR-66140.....77	JR-9000-0522F.....51	JR-CBN.....26
JR-0641-SS2.....54	JR-1125-10P-5.....75	JR-34N.....91	JR-55081-10.....17	JR-67140-1A.....78	JR-9000-0522H.....51	JR-CBULPK.....44
JR-0641-TI05.....54	JR-1125-2P-5.....75	JR-357090.....86	JR-55081-5.....17	JR-67140-2B.....78	JR-9000-0522HF.....51	JR-CBUMPk.....44
JR-0641-TI2.....54	JR-1125-5P-5.....75	JR-357180.....86	JR-55082-10.....17	JR-67-250-2A.....78	JR-9000-0525GC.....52	JR-CCPK-1.....24
JR-065.....44	JR-1127-05P-5.....75	JR-35P.....90	JR-55082-5.....17	JR-67-250-4.....78	JR-9000-0526G.....52	JR-C-EN.5FPKB.....37
JR-0650.....45	JR-1127-2P-5.....75	JR-35S.....90	JR-55083-10.....17	JR-681821NF.....77	JR-9000-0528G.....52	JR-C-EN1FPKB.....37
JR-066.....26	JR-1128-5P-5.....75	JR-367008-10.....53	JR-55083-5.....17	JR-68181NF.....77	JR-9000-0530.....51	JR-CFE-S10-5.....50
JR-067.....26	JR-1150-10P-5.....75	JR-367008-2.....53	JR-55084-10.....17	JR-68183NF.....77	JR-9000-0530F.....51	JR-CFE-S2-5.....50
JR-068.....44	JR-1151-10P-5.....75	JR-367008-20.....53	JR-55084-5.....17	JR-68184NF.....77	JR-9000-0531.....51	JR-CFE-S75-5.....50
JR-0680.....45	JR-1FR1-5.....73	JR-367016-10.....53	JR-55085-5.....20	JR-68185NF.....77	JR-9000-0531F.....51	JR-CFL-2W-5.....45
JR-070101-10.....48	JR-1FR2-5.....73	JR-367016-2.....53	JR-55085-10.....17	JR-68186NF.....77	JR-9000-0532.....51	JR-CFL-4D.....45
JR-070111-10.....48	JR-201094BK-10.....43	JR-367016-20.....53	JR-55085-5.....17	JR-68230-05.....57	JR-9000-0532F.....51	JR-CFLAPK.....29
JR-070112-10.....48	JR-201094WH-10.....43	JR-3675-2.....53	JR-55086-10.....17	JR-68230-2.....57	JR-9000-0602.....53	JR-CFL-CB1PK.....41
JR-070113-10.....48	JR-201144-10.....43	JR-3678-2.....53	JR-55086-5.....17	JR-68230-5.....57	JR-9000-0603.....53	JR-CFL-CB1PK.....41
JR-070123-10.....48	JR-20116-10.....42	JR-3678-25.....53	JR-55087-10.....17	JR-68231-05.....57	JR-9000-0604.....53	JR-CFL-CB1PK.....49
JR-070124-10.....48	JR-20124-10.....42	JR-4071-10.....79	JR-55087-5.....17	JR-68231-2.....57	JR-9000-0604H.....53	JR-CFL-CB2KF.....41
JR-070211-10.....48	JR-201531.....85	JR-4072-10.....79	JR-55088-10.....17	JR-68231-5.....57	JR-9000-0640.....53	JR-CFL-CB2KF.....49
JR-070212-10.....48	JR-201532.....85	JR-4073-10.....79	JR-55088-5.....17	JR-68242TI.....78	JR-9000-0641.....53	JR-CFL-CB2KF-S.....45
JR-070213-10.....48	JR-201536.....85	JR-4074-10.....79	JR-55089-10.....17	JR-68242TI2.....78	JR-9000-0642.....53	JR-CFL-CB4KF-S.....45
JR-070224-10.....48	JR-201537.....85	JR-4075-10.....79	JR-55089-5.....17	JR-68247.....54	JR-9000-0670.....28	JR-CM1XPK.....47
JR-070225-10.....48	JR-201539.....85	JR-408.....24	JR-5508X-12.....17	JR-68252TI.....78	JR-9001.....86	JR-CM2XPK.....47
JR-070350-10.....69	JR-201540.....41	JR-4081-10.....79	JR-5508X-24.....17	JR-68252TI2.....78	JR-B-C5M1PK.....46	JR-CMLAPK.....29
JR-071123-10.....48	JR-201540.....85	JR-4082-10.....79	JR-55090-10.....17	JR-68253.....54	JR-B-C5M2PK.....46	JR-C-NEU1FPK.....37
JR-071124-10.....48	JR-201540L.....85	JR-4083-10.....79	JR-55090-5.....17	JR-68256.....54	JR-B-C9M1PK.....46	JR-C-NEU.5FPK.....37
JR-071323-10.....48	JR-201541.....85	JR-4084-10.....79	JR-55091-10.....17	JR-68257.....54	JR-B-C9M2PK.....46	JR-C-NEU.5XFPK.....37
JR-071324-10.....48	JR-201554.....85	JR-4085-10.....79	JR-55091-5.....17	JR-68257PE.....54	JR-B-CM1XPK.....47	JR-CNFT.....19
JR-072350-10.....48	JR-201580BK-10.....43	JR-409.....24	JR-55100-5.....19	JR-68258.....55	JR-B-CM2XPK.....47	JR-C-NL.15L-5.....37
JR-072450-10.....48	JR-201580WH-10.....43	JR-410.....24	JR-5510-5.....23	JR-68262.....55	JR-BPR1.....88	JR-C-NL.15S-5.....37
JR-075013.....48	JR-201926-10.....41	JR-4676-10.....51	JR-55110-5.....19	JR-68263.....55	JR-BPR1.....88	JR-C-NL.20L-5.....37

JR-C-NL20S-5.....	37	JR-S-18314.....	68	JR-SLV-75-1000....	81	JR-T-5720-**-**	6	JR-T-6801-**-**	12	JR-T-97100.....	14	JR-ZU1XC.....	34
JR-C-NL25L-5.....	37	JR-S-20002.....	60	JR-SLV-75-20.....	81	JR-T-5764-**-**	6	JR-T-6802-**-**	11	JR-T-97105.....	14	JR-ZUFLPK.....	29
JR-C-NL25S-5.....	37	JR-S-20003.....	59	JR-SLV-75-250.....	81	JR-T-5802-**-**	6	JR-T-6802C-**-**	10	JR-T-97110.....	14	JR-ZX1.....	33
JR-C-NL30L-5.....	37	JR-S-20004.....	59	JR-SLV-75-5.....	81	JR-T-5803-**-**	6	JR-T-6803-**-**	11	JR-TD-5804-**-**	7	JR-ZX1C.....	33
JR-C-NL30S-5.....	37	JR-S-20005.....	60	JR-SLV-75-50.....	81	JR-T-5804-**-**	6	JR-T-6803C-**-**	10	JR-TD-6007-**-**	7	JR-ZX1M.....	33
JR-C-NL35L-5.....	37	JR-S-20006.....	59	JR-SLV-75-500.....	81	JR-T-5993-**-**	5	JR-T-6804-**-**	11	JR-TD-6008-**-**	7		
JR-C-NL35S-5.....	37	JR-S-20007.....	59	JR-SLV-76-10.....	82	JR-T-5993-**-**	8	JR-T-6805-**-**	12	JR-TD-6009-**-**	7		
JR-C-NNFFPK.....	37	JR-S-20008.....	59	JR-SLV-76-100.....	82	JR-T-5994-**-**	8	JR-T-6806-**-**	11	JR-TD-6010-**-**	7		
JR-CU4LPK.....	45	JR-S-20009.....	59	JR-SLV-76-1000....	82	JR-T-5995-**-**	8	JR-T-6806C-**-**	10	JR-TD-6011-**-**	7		
JR-F-197-062 2TI.....	78	JR-S-20100.....	71	JR-SLV-76-20.....	82	JR-T-5998-**-**	5	JR-T-6807-**-**	12	JR-TSS.505-**-**	15		
JR-F-197-062-2TINC.....	78	JR-S-20200.....	71	JR-SLV-76-200.....	82	JR-T-5998-**-**	8	JR-T-6808-**-**	12	JR-TSS.507-**-**	15		
JR-F-255-079-10PE-5.....	72	JR-S-40004.....	60	JR-SLV-76-5.....	82	JR-T-5998C-**-**	6	JR-T-6809.....	11	JR-TSS.510-**-**	15		
JR-FPKBSS.5-10.....	25	JR-S-40006.....	65	JR-SLV-76-50.....	82	JR-T-5999-**-**	6	JR-T-6810.....	12	JR-TSS.520-**-**	15		
JR-FPKBSS1-10.....	25	JR-S-40007.....	68	JR-SLV-76-500.....	82	JR-T-6000-**-**	6	JR-T-6811-**-**	12	JR-TSS105-**-**	15		
JR-FRK1.....	89	JR-S-40008.....	62	JR-T-078-**-**	9	JR-T-6001-**-**	6	JR-T-6812-**-**	11	JR-TSS110-**-**	15		
JR-FRK2.....	89	JR-S-40009.....	65	JR-T-080-**-**	9	JR-T-6002-**-**	6	JR-T-6812C-**-**	10	JR-TSS115-**-**	15		
JR-HKS3.....	89	JR-S-40104.....	68	JR-T-082-**-**	9	JR-T-6003-**-**	6	JR-T-6813-**-**	11	JR-TSS120-**-**	15		
JR-LZN1-5.....	31	JR-S-40106.....	68	JR-T-083-**-**	9	JR-T-60031-**-**	6	JR-T-6835-**-**	11	JR-TSS130-**-**	15		
JR-LZNF1KF-5.....	21	JR-S-40108.....	68	JR-T-084-**-**	9	JR-T-60032-**-**	6	JR-T-96005.....	14	JR-TSS140-**-**	15		
JR-LZNF1PK-5.....	21	JR-SLC-85-10.....	81	JR-T-085-**-**	9	JR-T-6004-**-**	5	JR-T-96010.....	14	JR-TSS230-**-**	15		
JR-LZNF1PKG-5.....	21	JR-SLC-85-100.....	81	JR-T-086-**-**	9	JR-T-60041-**-**	5	JR-T-96015.....	14	JR-TSS240-**-**	15		
JR-LZP1PK.....	24	JR-SLC-85-1000....	81	JR-T-087-**-**	9	JR-T-60042-**-**	5	JR-T-96020.....	14	JR-TSS260-**-**	15		
JR-MZN1-5.....	31	JR-SLC-85-20.....	81	JR-T-088-**-**	9	JR-T-6006.....	5	JR-T-96025.....	14	JR-TSS267-**-**	15		
JR-MZNF1KF-5.....	21	JR-SLC-85-200.....	81	JR-T-36100-**-**	6	JR-T-6007-**-**	8	JR-T-96030.....	14	JR-TSS285-**-**	15		
JR-MZNF1PK-5.....	21	JR-SLC-85-250.....	81	JR-T-36125-**-**	6	JR-T-6008-**-**	8	JR-T-96035.....	14	JR-ZBU1CFPK.....	26		
JR-MZNF1PKG-5.....	21	JR-SLC-85-300.....	81	JR-T-36150-**-**	6	JR-T-6009-**-**	8	JR-T-96040.....	14	JR-ZBU1FPK.....	26		
JR-MZP1PK.....	24	JR-SLC-85-5.....	81	JR-T-3650-**-**	6	JR-T-6010-**-**	8	JR-T-96105.....	14	JR-ZBU1MFPK.....	26		
JR-PTAS.....	49	JR-SLC-85-50.....	81	JR-T-3675-**-**	6	JR-T-6011-**-**	8	JR-T-96110.....	14	JR-ZC.5.....	32		
JR-S-11001.....	64	JR-SLC-85-500.....	81	JR-T-37100-**-**	8	JR-T-6020-**-**	8	JR-T-96115.....	14	JR-ZC.5FPK.....	24		
JR-S-11002.....	64	JR-SLC-86-10.....	82	JR-T-37125-**-**	8	JR-T-6030-**-**	8	JR-T-96120.....	14	JR-ZC1.....	32		
JR-S-11003.....	64	JR-SLC-86-100.....	82	JR-T-37150-**-**	8	JR-T-6040-**-**	8	JR-T-96125.....	14	JR-ZC1PK.....	24		
JR-S-11005.....	64	JR-SLC-86-1000....	82	JR-T-3750-**-**	8	JR-T-6060-**-**	5	JR-T-96130.....	14	JR-ZC2.....	32		
JR-S-11006.....	64	JR-SLC-86-20.....	82	JR-T-3775-**-**	8	JR-T-6060-**-**	8	JR-T-96135.....	14	JR-ZF.5-5.....	32		
JR-S-11011.....	63	JR-SLC-86-250.....	82	JR-T-4001-**-**	9	JR-T-6105-**-**	11	JR-T-96140.....	14	JR-ZF.5PK-5.....	22		
JR-S-11012.....	63	JR-SLC-86-5.....	82	JR-T-4002-**-**	9	JR-T-6130-**-**	13	JR-T-96145.....	14	JR-ZF1-5.....	32		
JR-S-11013.....	63	JR-SLC-86-50.....	82	JR-T-4003-**-**	9	JR-T-6210-**-**	11	JR-T-96205.....	14	JR-ZF1S6-5.....	33		
JR-S-11023.....	66	JR-SLC-86-500.....	82	JR-T-4004.....	9	JR-T-625-04-**-**	15	JR-T-96210.....	14	JR-ZF2-5.....	32		
JR-S-11101.....	64	JR-SLR-111.....	83	JR-T-4005.....	9	JR-T-625-05-**-**	15	JR-T-96215.....	14	JR-ZF2PK-5.....	22		
JR-S-11102.....	64	JR-SLR-111-77.....	83	JR-T-4007-**-**	9	JR-T-625-10-**-**	15	JR-T-96220.....	14	JR-ZGF.5PK.....	37		
JR-S-11103.....	64	JR-SLR-112.....	83	JR-T-4009.....	9	JR-T-625-20-**-**	15	JR-T-96225.....	14	JR-ZGF1PK.....	37		
JR-S-11111.....	62	JR-SLR-112-77.....	83	JR-T-4011-**-**	12	JR-T-625-30-**-**	15	JR-T-96230.....	14	JR-ZN.5-5.....	31		
JR-S-11112.....	62	JR-SLR-113.....	83	JR-T-4011C-**-**	12	JR-T-625-40-**-**	15	JR-T-96235.....	14	JR-ZN.5FPK-5.....	22		
JR-S-11113.....	62	JR-SLR-113-77.....	83	JR-T-4036-**-**	12	JR-T-626-00-**-**	15	JR-T-96240.....	14	JR-ZN1-5.....	31		
JR-S-11201.....	64	JR-SLR-114.....	83	JR-T-4037-**-**	12	JR-T-628-00.....	15	JR-T-96305.....	14	JR-ZN2-5.....	31		
JR-S-11202.....	64	JR-SLR-114-77.....	83	JR-T-4039-**-**	12	JR-T-6310-**-**	11	JR-T-96310.....	14	JR-ZN2PK-5.....	22		
JR-S-11203.....	64	JR-SLR-115.....	83	JR-T-4041.....	12	JR-T-6320-**-**	11	JR-T-96315.....	14	JR-ZNF1FKFK-5.....	21		
JR-S-11211.....	62	JR-SLR-115-77.....	83	JR-T-4136C-**-**	12	JR-T-6604C-**-**	10	JR-T-96320.....	14	JR-ZNF1FPK-5.....	21		
JR-S-11212.....	62	JR-SLR-116.....	83	JR-T-4183-**-**	12	JR-T-6605-**-**	11	JR-T-96325.....	14	JR-ZNF1FPKG-5.....	21		
JR-S-11213.....	62	JR-SLR-116-77.....	83	JR-T-4183C-**-**	12	JR-T-6605C-**-**	10	JR-T-96330.....	14	JR-ZNF1KF-5.....	21		
JR-S-13001.....	65	JR-SLR-117.....	83	JR-T-436100-**-**	9	JR-T-6606-**-**	11	JR-T-96335.....	14	JR-ZNF1KF-F-5.....	21		
JR-S-13002.....	65	JR-SLR-117-77.....	83	JR-T-43650-**-**	9	JR-T-6606C-**-**	10	JR-T-96340.....	14	JR-ZNF1PK-5.....	21		
JR-S-13003.....	65	JR-SLR-118.....	83	JR-T-43675-**-**	9	JR-T-6607-**-**	11	JR-T-97005.....	14	JR-ZNF1PK-F-5.....	21		
JR-S-14002.....	65	JR-SLR-118-77.....	83	JR-T-5602-**-**	8	JR-T-6607C-**-**	10	JR-T-97010.....	14	JR-ZNF1PKG-5.....	21		
JR-S-15001.....	65	JR-SLR-2190.....	81	JR-T-5603-**-**	8	JR-T-6610-**-**	11	JR-T-97015.....	14	JR-ZNF1PKG-F-5.....	21		
JR-S-15002.....	65	JR-SLR-2191.....	81	JR-T-5604-**-**	8	JR-T-6610C-**-**	10	JR-T-97020.....	14	JR-ZNF2S6-5.....	34		
JR-S-15003.....	65	JR-SLR-2192.....	81	JR-T-5605-**-**	8	JR-T-6615-**-**	11	JR-T-97025.....	14	JR-ZNFT.....	19		
JR-S-15011.....	63	JR-SLR-2193.....	81	JR-T-5606-**-**	8	JR-T-6620-**-**	11	JR-T-97030.....	14	JR-ZP.5.....	32		
JR-S-15012.....	63	JR-SLR-2194.....	81	JR-T-5607-**-**	8	JR-T-6620C-**-**	10	JR-T-97035.....	14	JR-ZP.5FPK.....	24		
JR-S-15013.....	63	JR-SLR-2195.....	81	JR-T-5608-**-**	6	JR-T-6702-**-**	11	JR-T-97040.....	14	JR-ZP1.....	32		
JR-S-15103.....	67	JR-SLR-2196.....	81	JR-T-5608-**-**	8	JR-T-6704-**-**	11	JR-T-97045.....	14	JR-ZP1FPK.....	24		
JR-S-15113.....	67	JR-SLR-2200.....	81	JR-T-5610-**-**	8	JR-T-6704C-**-**	10	JR-T-97050.....	14	JR-ZP2.....	32		
JR-S-15303.....	68	JR-SLR-2201.....	81	JR-T-5620-**-**	8	JR-T-6705-**-**	11	JR-T-97055.....	14	JR-ZRU1.5FPK.....	27		
JR-S-15313.....	68	JR-SLR-2202.....	81	JR-T-5664-**-**	8	JR-T-6705C-**-**	10	JR-T-97060.....	14	JR-ZRU21.....	34		
JR-S-16003.....	67	JR-SLR-2204.....	81	JR-T-5702-**-**	6	JR-T-6706-**-**	11	JR-T-97065.....	14	JR-ZRU21FPK.....	27		
JR-S-16013.....	67	JR-SLR-2206.....	81	JR-T-5703-**-**	6	JR-T-6706C-**-**	10	JR-T-97070.....	14	JR-ZT1.....	33		
JR-S-16103.....	67	JR-SLR-2207.....	81	JR-T-5704-**-**	6	JR-T-6707-**-**	11	JR-T-97075.....	14	JR-ZT1C.....	33		
JR-S-16113.....	67	JR-SLR-2210.....	81	JR-T-5705-**-**	6	JR-T-6707C-**-**	10	JR-T-97080.....	14	JR-ZT1M.....	33		
JR-S-18004.....	66	JR-SLR-2222.....	81	JR-T-5707-**-**	6	JR-T-6708-**-**	11	JR-T-97085.....	14	JR-ZU1CS6.....	34		
JR-S-18014.....	66	JR-SLV-75-10.....	81	JR-T-5710-**-**	6	JR-T-6794-**-**	11	JR-T-97090.....	14	JR-ZU1MS6.....	34		
JR-S-18304.....	68	JR-SLV-75-100.....	81	JR-T-5715-**-**	8	JR-T-6800-**-**	12	JR-T-97095.....	14	JR-ZU1S6.....	34		

## GENERAL REFERENCE

**A**

Adapters	
High pressure PEEK .....	29
Low pressure polypropylene .....	48
Air inlet valve .....	59
Air outlet filter .....	59

**B**

Backpressure regulator .....	88
Barbed hose adapter	
for 1/8" ID tubing .....	63
Bottles, safety-coated glass .....	70
Bulkhead unions, PEEK	
High pressure .....	26
Low pressure .....	44
Bushing for high pressure	
1/4-28 connection .....	25

**C**

Caps	
Nanovolume® 360 µm ...	38-39
PEEK .....	24
Stainless .....	32
Carbon dioxide filter .....	71
ChromBox kit .....	91
Clean-Cut tubing cutter .....	86
Color-It fingertight adapters .....	17
Color-Sleeve	
fingertight adapters .....	17
Column couplers .....	79
Column end plugs	
Color-coded .....	79
Logotype .....	79
Columns, PEEK	
2.1 mm .....	77
4.0 mm .....	77
Crosses	
Nanovolume® .....	38
PEEK	
High pressure .....	28
Low pressure .....	46
Stainless .....	33
Cutting wheel SS tubing cutter ..	87

**E**

ETFE tubing .....	9
Easy-Flange kit .....	85
Easy-Flange combi kit .....	85
Eco cap .....	69
Eco cap kit .....	69
Economy Last Drop	
mobile phase filters .....	51
Elbows for PEEK tubing .....	86
Exhaust filter .....	60

**F**

FEP tubing	
Natural color .....	11
Solid color-coded .....	10
Ferrule removal kits .....	89
Ferrules	
High pressure	
Collapsible .....	37
Nanovolume® .....	37
PEEK	
Collapsible .....	37
Double-ended .....	23
Single .....	22
UHPLC .....	NEW 22

**F, continued**

Ferrules, continued	
High pressure, cont'd	
Stainless	
Metric .....	35
Rheodyne .....	32
Valco .....	32
Low pressure	
Collapsible .....	41
Inverted .....	41
File .....	87
Filters	
UHPLC stainless	
In-line filter .....	57
Precolumn filter .....	56
High pressure PEEK	
In-line filter kit .....	54
Precolumn filter kit ...	55
Sure-Guard .....	54
Low pressure	
Economy Last Drop .....	51
In-line .....	50
Last Drop .....	51
Last Drop glass .....	52
Mobile phase filters ...	51-52
No-Met biocompatible ...	52
Safety	
Air outlet filter .....	59
Carbon dioxide filter .....	71
Humidity filter .....	71
Filter/spargers	
Last Drop .....	53
Stainless steel .....	53
Fingertight adapters	
Color-It .....	17
Color-Sleeve .....	17
Fingertight one-piece column	
couplers .....	79
Flanged nuts	
Fingertight color-coded .....	43
Hex-head nuts .....	41
Flangeless nuts	
Fingertight color-coded .....	42
Fingertight PEEK .....	42
Hex-head nuts .....	41
Frits	
for PEEK filter kits .....	55
for stainless in-line filter .....	57
PEEK-encased	
Stainless .....	74
Titanium .....	75
Polyethylene .....	75
Unidense	
Stainless .....	73
Titanium .....	73
Polyethylene .....	72

**G**

Guard columns,	
in-line biocompatible .....	78
Guillotine tubing cutter .....	86

**H**

Hex key set .....	89
Humidity filter .....	71

**I**

Impermeable	
polymeric tubing .....	13
In-line biocompatible	
guard column .....	78
In-line filter cartridges .....	50
In-line filter, stainless steel .....	57
In-line filter kit	
PEEK .....	54
Sure-guard .....	54

**K**

Kits	
Chrombox kit .....	91
Easy-Flange kit .....	85
Easy-Flange combi kit .....	85
PEEK starter kit .....	90
Stainless steel starter kit .....	90

**L**

Last Drop biocompatible glass	
mobile phase filters .....	52
Last Drop filter/sparger .....	53
Last Drop mobile phase filters ...	51
Liners for 1/32" Nanovolume®	
fittings .....	37
Loops	
PEEK .....	80-81
Stainless steel .....	82-83

**M**

Manifolds, low pressure PEEK ...	46
Micro valves for GC and LC .....	49
Mixing tees	
Biocompatible high pressure	28
Low pressure .....	47
Mobile phase filters	
Biocompatible glass .....	52
Economy Last Drop .....	51
Last Drop .....	51
No-Met biocompatible .....	52

**N**

Nanovolume® fittings	
for 1/32" tubing .....	37
for 360 µm tubing .....	38-39
No-Met biocompatible	
mobile phase filters .....	52
Nuts	
High pressure	
PPS .....	25
PEEK	
for double ferrules .....	23
for single ferrules .....	22
Stainless	
Rheodyne compatible	31
Valco .....	31
Low pressure	
Flanged	
Fingertight	
color-coded .....	43
Hex-head .....	41
Flangeless	
Fingertight	
color-coded .....	42
Fingertight PEEK .....	42
Hex-head .....	41
Nanovolume®	
1/32" .....	37
360 micron .....	38-39

**O**

One-piece fittings	
Column couplers .....	79
Nanovolume® nut-ferrule .....	38-39
No-twist PEEK .....	21
PPS hex-head .....	19
PEEK fingertight	
Color-coded .....	20
Natural .....	20
PEEK hex-head	
5000 psi .....	19
UHPLC .....	NEW 19

**P**

PEEK columns .....	77
PEEK starter kit .....	89
PEEK tubing	
Dual layer color-coded .....	7
Natural .....	5
Solid color-coded .....	8
Striped color-coded .....	6
PEEK tubing elbows .....	86
PFA tubing .....	9
PTFE tubing	
Natural color .....	12
Solid color-coded .....	12
Packing kits	
for PEEK columns .....	78
for guard columns .....	78
Peristaltic tubing adapter set .....	49
Plier type SS tubing cutter .....	87
Plugs	
High pressure	
PEEK .....	24
Stainless .....	32
Low pressure	
Polypropylene .....	48
Precolumn filter kit, PEEK .....	55
Precolumn filter, stainless .....	56

**R**

Reducing unions	
Nanovolume® .....	37-39
PEEK, high pressure .....	27
PEEK, low pressure	
Imperial .....	45
Imperial to metric NEW ..	45
Stainless	
Imperial .....	34
Imperial to metric NEW ..	35

**S**

Safety	
Air inlet valves .....	59
Air outlet filters .....	59
Exhaust filters .....	60-61
Caps .....	63
Cap kit .....	62
Carbon dioxide filter .....	71
Humidity filter .....	71
Safety-coated glass bottles .....	70
Sample loops	
PEEK .....	80-81
Stainless steel .....	82-83

**S, continued**

SolviFlex coiled tubing.....	70
Stainless steel in-line filter .....	57
Stainless steel mobile phase filter/ helium sparger.....	53
Stainless steel starter kit .....	90
Stainless steel tubing	
Pre-cut.....	14
Pre-cut flexible .....	14
Premium grade .....	15
Standard grade .....	15
Stainless steel tubing cutters	
Cutting wheel type.....	87
Plier type .....	87
Starter kits	
Chrombox kit .....	91
PEEK starter kit .....	90
Stainless steel starter kit.....	90
Sure-Guard in-line filter kit .....	54

**T**

Tees	
Biocompatible mixing .....	28
Nanovolume®.....	38
PEEK	
High pressure.....	28
Low pressure .....	46
Mixing .....	47
Stainless .....	33

**T, continued**

Thermo-formed	
polymeric tubing.....	13
Tubing	
ETFE.....	9
FEP	
Natural color.....	11
Solid color-coded .....	10
Impermeable	
polymeric tubing .....	13
PEEK	
Dual layer color-coded ....	7
Natural .....	5
Solid color-coded .....	8
Striped color-coded.....	6
PFA.....	9
PTFE	
Natural color.....	12
Solid color-coded .....	12
Stainless	
Pre-cut .....	14
Pre-cut flexible.....	14
Premium grade .....	15
Standard grade .....	15
Thermo-formed	
polymeric tubing .....	13
Tubing clip .....	86

**T, continued**

Tubing cutters	
Clean-Cut (for polymers).....	86
Cutting wheel type	
(for stainless) .....	87
Guillotine type (for polymers)	86
Plier type (for stainless.....	87
Tubing elbows .....	86
Tweezer.....	89

**U**

UHPLC	
Filters	
Stainless in-line filter .....	57
Stainless precolumn	
filter..... <b>NEW</b> ..	56
PEEK	
Ferrules..... <b>NEW</b> ..	22
HT™ one-piece hex-head	
fittings..... <b>NEW</b> ..	19
Stainless steel	
Caps.....	32
Crosses.....	33
Ferrules.....	32
Nuts.....	31
Plugs .....	32
Reducing unions.....	34, 35
Tees .....	33
Unions .....	34

**U, continued**

Unidense frits .....	72-73
Unions	
Body only, for	
low-pressure flanged.....	43
Bulkhead ... <i>see Bulkhead unions</i>	
Nanovolume®.....	37-39
PEEK	
High pressure.....	26
Low pressure .....	44
Reducing ... <i>see Reducing unions</i>	
Stainless .....	34

**V**


VICI caps .....	64-65
VICI cap kit .....	64
ValvTool.....	88

**W**

Washers	
for flanged connections .....	41
Waste caps .....	66-68
Wrenches	
Open end.....	89
ValvTool .....	88



- Valco valves
- Cheminert valves
- Actuators
- Valco fittings
- Cheminert fittings
- Tubing
- Flow, pressure, and on/off devices
- Dynacal® permeation devices
- Instrumentation
- GC columns
- GC purifiers

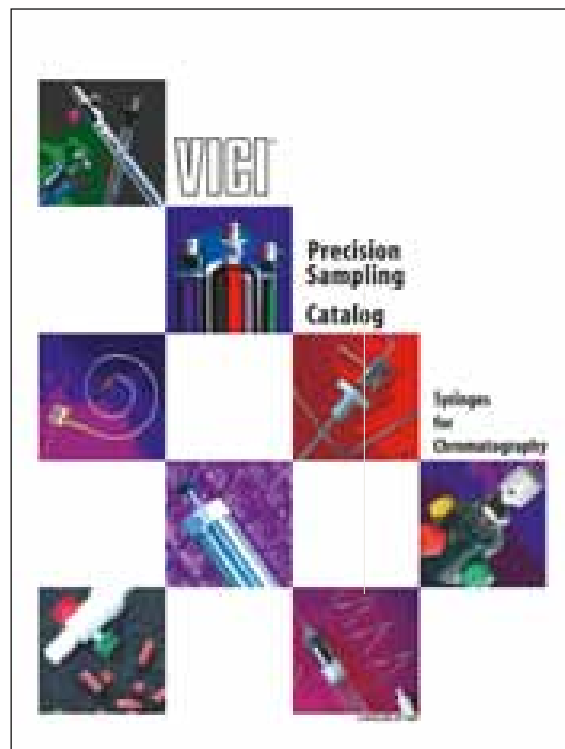

- Gas tight syringes
- HPLC injection syringes
- Magnum syringes up to 500 mL
- Syringes for sampling from gas cylinders
- Mininert™ – Valves for vials
- Customized probes and shaped tubing



## TO ORDER

Our sales personnel will work with you by email, phone, or fax to determine the best products for your specifications.

Our Schenkon, Switzerland, Brockville, Ontario, and Baton Rouge, Louisiana offices serve different regions of the world.



### Europe, Middle East, Asia, Africa, Australia, and other countries

#### VICI AG International

##### Sales

Phone: Int + 41· 41· 925-6200

Fax: Int + 41· 41· 925-6201

Email: [order@vici.ch](mailto:order@vici.ch)

##### Technical/Service

Phone: Int + 41· 41· 925-6200

Fax: Int + 41· 41· 925-6201

Email: [support@vici.ch](mailto:support@vici.ch)

Parkstrasse 2  
CH-6214 Schenkon / Switzerland  
[info@vici.ch](mailto:info@vici.ch)

EXW shipping point, our factory,  
Schenkon, Switzerland.

We accept:



### North, South, and Central America (except Canada)

#### VICI Precision Sampling

##### Sales

Toll-free: 800-828-1653

Phone: 225-927-1128

Fax: 225-923-1331

Email: [precsamp@vici.com](mailto:precsamp@vici.com)

##### Technical/Service

Phone: 225-927-1128

Fax: 225-923-1331

Email: [precsamp@vici.com](mailto:precsamp@vici.com)

P.O. Box 15886  
Baton Rouge, Louisiana / USA  
[precsamp@vici.ch](mailto:precsamp@vici.ch)

EXW shipping point, our factory,  
Schenkon, Switzerland.

We accept:



### Canada

#### VICI Valco Instruments Canada

##### Sales

Toll-free: 866-297-2626

Phone: 613-342-2600

Fax: 613-342-0111

Email: [canada@vici.com](mailto:canada@vici.com)

##### Technical/Service

Phone: 613-342-2600

Fax: 613-342-0111

Email: [canada@vici.com](mailto:canada@vici.com)

26 Water Street East  
Brockville, Ontario / Canada K6V 1A1  
[canada@vici.com](mailto:canada@vici.com)

F.O.B. shipping point, our factory,  
Brockville, Canada.

We accept:



## Conditions of Sale

There is no minimum order.  
Terms are net 30 days.  
Prices subject to change without notice.

## Notice

VICI reserves the right to correct errors and change prices, designs, or specifications without notice or liability. The information in this catalog is correct to the best of our knowledge but is not guaranteed to be so. VICI assumes no responsibility with respect thereto.

## Repairs

For repair return authorization or information about factory refurbishments, contact our Service Department.

## Returns

No returns will be accepted more than 90 days after shipment for any reason. Before 90 days, no returns will be accepted without prior authorization. If it is necessary to return material to us, please contact our Sales Department for a return authorization number and forwarding instructions. Inspect shipments upon receipt and report

## Technical Service

Our capable staff is waiting to assist you with your technical questions and application needs.

shortages and incorrect or damaged material to us immediately.

*Important!* Damaged shipments must remain with the original packaging for freight company inspection.

Returned material will be subject to a restocking charge of 20% for catalog items. Special orders cannot be returned unless defective.

**CANADA**

VICI Valco Instruments Canada Corp  
26 Water Street East  
CA-Brockville, ON K6V 1A1  
Phone: +1 613 342 2600  
Toll free: +1 866 297 2626  
Fax: +1 613 342 0111  
Email: canada@vici.com  
Web: www.vici.com

**CHINA**

J & J Industries  
Room 917/918, Prime Tower  
22 Chaoyangmen Wai Da Jie  
CN-Beijing 100020  
Phone: +86 10 65888666  
Fax: +86 10 65886577  
Email: sales@jjindustries.com.cn  
Web: www.jjindustries.com.cn

**DENMARK**

Mikrolab Aarhus A/S  
Axel Kiers Vej 34  
DK-8270 Højbjerg  
Phone: +45 86 29 61 11  
Fax: +45 86 29 61 22  
Email: ml@mikrolab.dk  
Web: www.mikrolab.dk

**FINLAND**

Krotek O.Y.  
Salmintie 4B  
FI-02140 Espoo  
Phone: +358 9 412 88 08  
Fax: +358 9 412 88 09  
Email: krotek@krotek.fi  
Web: www.krotek.fi

**FRANCE**

Interchim SA  
211 Bis, Avenue J.F. Kennedy  
BP 1140  
FR-03103 Montluçon Cedex  
Phone: +33 4 700 388 55  
Fax: +33 4 700 382 60  
Email: interchrom@interchim.com  
Web: www.interchim.com

**GERMANY**

Machery-Nagel GmbH + Co. KG  
Neumann-Neander-Strasse 6-8  
DE-52355 Düren  
Phone: +49 2421 969 0  
Fax: +49 2421 969 122  
Email: sales-de@mn-net.com  
Web: www.mn-net.com

**IRELAND**

Carl Stuart Ltd.  
Tallaght Business Park  
Whitestown  
IE-Dublin 24  
Phone: +353 1 452 34 32  
Fax: +353 1 452 39 67  
Email: info@carlstuart.com  
Web: www.carlstuart.com

**ITALY**

C.P.S. Analitica S.r.l.  
Via Neera 8/A  
IT-20141 Milano  
Phone: +39 02 895 4201  
Fax: +39 02 895 42022  
Email: cps@cps.it  
Web: www.cps.it

**ITALY**

Labservice Analytica S.r.l.  
Via Emilia 51/C  
IT-40011 Anzola Emilia/Bologna  
Phone: +39 051 732351  
Fax: +39 051 732759  
Email: info@labservice.it  
Web: www.labservice.it

**JAPAN**

GL Sciences Inc.  
22-1 Nishishinjuku 6-chome  
Shinjuku-Ku  
JP-163-1130 Tokyo  
Phone: +81 3 5323 6620  
Fax: +81 3 5323 6621  
Email: world@gls.co.jp  
Web: www.gls.co.jp

**NETHERLANDS**

Da Vinci Europe Laboratory Solutions B.V.  
P.O. Box 12.103  
NL-3004 GC Rotterdam  
Phone: +31 10 258 18 70  
Fax: +31 10 258 18 79  
Email: solutions@davincieurope.com  
Web: www.davincieurope.com

**NORWAY**

Teknolab AS  
Trollasveien 34  
NO-1414 Trollasen  
Phone: +47 66 81 34 70  
Fax: +47 66 81 34 71  
Email: mail@teknolab.no  
Web: www.teknolab.no

**RUSSIA**

Fizlabpribor Ltd.  
Podolskiy Kursantov str. 3/2  
RU-117535 Moscow  
Phone: +7 495 740 5406  
Fax: +7 495 316 9600  
Email: info@fizlabpribor.ru  
Web: www.fizlabpribor.ru

**SINGAPORE**

Nareen Technologies PTE Ltd.  
196, Pandan Loop, # 07-18  
Pantech Industrial Complex  
SG-Singapore 128384  
Phone: +65 6777 4001  
Fax: +65 6779 4401  
Email: enquiry@nareen.com.sg  
Web: www.nareen.com

**SOUTH KOREA**

Young-Wha Scientific Co., Ltd.  
5th Floor, Gu-Jung B/D  
577-7 Shinsa-Dong, Kangnam-Ku  
KR-135-891 Seoul, South Korea  
Phone: +82 2 2140 5400  
Fax: +82 2 2140 5405  
Email: ywsc@youngwha.com  
Web: www.youngwha.com

**SPAIN**

Scharlab S.L.  
Gato Pérez, 33, Polg. Industrial Mas  
D'en Cisa  
ES-08181 Sentmenat  
Phone: +34 93 745 64 00  
Fax: +34 93 349 80 23  
Email: purchase02@scharlab.com  
Web: www.scharlab.com

**SWEDEN**

Scantec Lab AB  
Box 238  
SE-43324 Partille  
Phone: +46 31 336 90 00  
Fax: +46 31 336 90 36  
Email: scantec.info@scanteclab.se  
Web: www.scanteclab.se

**SWITZERLAND**

BGB Analytik AG  
Rohrmattstrasse 4  
CH-4461 Böckten  
Phone: +41 61 991 00 46  
Fax: +41 61 991 00 25  
Email: sales@bgb-analytik.ch  
Web: www.bgb-analytik.ch

**THAILAND**

Amani Corporation Ltd.  
201/220 M6 Pravate, Pravate  
TH-10250 Bangkok  
Phone: +66 (2) 322 44 22  
Fax: +66 (2) 322 67 05  
Email: amani@ksc.th.com

**UNITED KINGDOM**

Thames Restek UK Ltd.  
Units 8/16, Ministry Wharf  
Wycombe Road, Saunderton  
GB-Buckinghamshire HP14 4HW  
Phone: +44 1494 563377  
Fax: +44 1494 564990  
Email: sales@thamesrestek.co.uk  
Web: www.thamesrestek.co.uk

**USA**

VICI Precision Sampling  
8275 W. El Cajon Drive  
US-70815, Baton Rouge / LA  
Phone: +1 225 927 1128  
Toll free: +1 800 828 1653  
Fax: +1 225 923 1331  
Email: precsamp@vici.com  
Web: www.vici.com

**NOTE**

Our products could also be available from your local distributor. Please contact VICI AG International for other countries.

[www.vici-jour.com](http://www.vici-jour.com)

**VICI Jour®**

VICI AG International

Parkstrasse 2  
6214 Schenkon  
Switzerland

Phone +41 41 925 62 00

Fax +41 41 925 62 01

Email [info@vici.ch](mailto:info@vici.ch)



**printed carbon neutral**

Energy efficient and CO<sub>2</sub> compensated print  
SC2012030701 - grempel.ch

© 2012 by VICI AG International

Printed in Switzerland