


Baskets & Basket Shafts

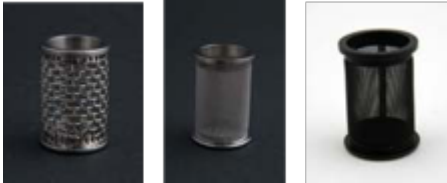

QLA offers both clip on style and O-ring style Baskets for most dissolution baths. Both meet the associated Original Equipment Manufacturer's dimensional specifications. Clip style baskets comply with USP (711) specifications. All QLA baskets are made from electropolished, sintered 316 Stainless Steel. The sintering process essentially creates a welded joint at each overlapping wire junction. Sintered baskets are stronger and last longer. They retain original dimensions and geometry and resist wobble. Most conventional baskets are not sintered and weaken over time. Critical basket geometry is verified on an optical comparator. All QLA 316 SS baskets are serialized for identification purposes.